



Sigma 8KS

Sigma 8KBS

from serial no. 143558



Refrigerated Centrifuge



Operating Manual

Please retain for later use!

In case of inquiries, please state the following numbers:

Order number:

Serial number:

© Copyright by
Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode am Harz
Germany

Tel.: +49 (0) 5522 / 5007-0
Fax: +49 (0) 5522 / 5007-12
Web: www.sigma-zentrifugen.de
E-mail: info@sigma-zentrifugen.de

1	General information.....	9
1.1	Importance of the operating manual.....	9
1.2	Intended use.....	9
1.3	Warranty and liability.....	9
1.4	Copyright.....	10
1.5	Standards and regulations.....	10
1.6	Scope of supply.....	10
2	Layout and mode of operation.....	11
2.1	Layout of the centrifuge.....	11
2.1.1	Functional and operating elements.....	11
2.1.2	Name plate.....	14
2.1.3	Mains plug.....	14
2.2	Mode of operation.....	15
2.2.1	Centrifugation principle.....	15
2.2.2	Area of application.....	15
2.2.2.1	Speed, radius, and relative centrifugal force.....	16
2.2.2.2	Density.....	16
3	Safety.....	17
3.1	Marking of the unit.....	17
3.2	Explanation of the symbols and notes.....	18
3.3	Responsibility of the operator.....	19
3.4	Operating personnel.....	19
3.5	Informal safety instructions.....	19
3.6	Safety instructions.....	20
3.6.1	Electrical safety.....	20
3.6.2	Mechanical safety.....	20
3.6.3	Fire prevention.....	21
3.6.4	Chemical and biological safety.....	21
3.6.5	Safety instructions for centrifugation.....	22
3.6.6	Resistance of plastics.....	22
3.6.7	Rotors and accessories.....	23
3.6.7.1	Lifting and carrying rotors.....	23
3.6.7.2	Service life.....	23
3.7	Safety devices.....	24
3.7.1	Lid lock device.....	24
3.7.2	Standstill monitoring system.....	24
3.7.3	System check.....	24
3.7.4	Earth conductor check.....	24
3.7.5	Imbalance monitoring system.....	24
3.7.6	Temperature monitoring system.....	24
3.7.7	Rotor monitoring system.....	24
3.8	Measures in the event of hazards and accidents.....	25
3.9	Remaining hazards.....	25

Table of contents

4	Storage and transport	26
4.1	Dimensions and weight.....	26
4.2	Storage conditions	26
4.3	Notes on transport	26
4.4	Packaging.....	27
5	Set-up and connection	31
5.1	Installation site	31
5.2	Set-up and alignment of the centrifuge	31
5.2.1	Set-up.....	32
5.2.2	Alignment	32
5.3	Installation of the cover holder	34
5.4	Power supply	35
5.4.1	Type of connectiont	35
5.4.2	Customer-provided fuses.....	35
5.4.3	Power isolating device	35
5.5	Condensate drain	35
5.6	Special equipment: water cooling system	37
6	Using the centrifuge	38
6.1	Initial start-up	38
6.2	Switching the centrifuge on	38
6.2.1	Opening and closing the lid.....	38
6.2.2	Installation of rotors and accessories	39
6.2.2.1	Installation of the rotor.....	39
6.2.2.2	Installation of angle rotors with a hermetically sealed lid	39
6.2.2.3	Installation of accessories	40
6.2.2.4	Adapters	41
6.2.2.5	Tubes.....	41
6.2.2.6	Blood bag systems.....	42
6.3	Control system "Spincontrol S"	43
6.3.1	User interface	43
6.3.2	Manual mode.....	44
6.3.2.1	Starting a centrifugation run	44
6.3.2.2	Interrupting a centrifugation run	44
6.3.2.3	Interrupting a deceleration process	44
6.3.2.4	Selection, display, and modification of data.....	44
6.3.2.5	Standard menu	45
6.3.2.6	Process library menu	51
6.3.2.7	Parameters menu	52
6.3.2.8	Setup menu	56
6.3.2.9	Curve menu	59
6.3.2.10	Option: Barcode menu	60
6.3.2.11	Help menu	62
6.3.2.12	Changing the contrast.....	63

Table of contents

6.3.3	Program mode.....	63
6.3.3.1	Saving a program.....	64
6.3.3.2	Loading a program.....	64
6.3.3.3	Executing a program.....	65
6.3.3.4	Deleting a program	65
6.3.3.5	Automatic program rotation.....	66
6.3.4	Options for data input and output.....	67
6.4	Switching the centrifuge off.....	67
7	Malfunctions and error correction.....	68
7.1	General malfunctions.....	68
7.1.1	Emergency lid release	69
7.2	Table of error codes.....	70
7.3	Service contact	71
8	Maintenance and service.....	72
8.1	Maintenance	72
8.1.1	Centrifuge.....	72
8.1.2	Condenser (only refrigerated centrifuges).....	73
8.1.3	Accessories	73
8.1.3.1	Plastic accessories	74
8.1.4	Rotors, buckets and carriers.....	74
8.1.5	Load bearing bolts.....	74
8.1.6	Glass breakage	75
8.2	Sterilisation and disinfection of the rotor chamber and accessories	75
8.2.1	Autoclaving.....	76
8.3	Service	77
8.4	Return of defective parts.....	78
9	Disposal.....	80
9.1	Disposal of the centrifuge	80
9.2	Disposal of the packaging.....	80
10	Technical data.....	81
10.1	Ambient conditions	82
10.2	Technical documentation	82

Table of contents

11	Appendix	83
11.1	Range of accessories	83
11.1.1	Maximum speed for tubes.....	85
11.1.2	Rotor radii.....	85
11.2	Speed-gravitational-field-diagram	86
11.3	Acceleration and deceleration curves	87
11.4	Table of the service life of rotors and accessories.....	89
11.5	Resistance data.....	90
11.6	Serial Control Interface	95
11.7	EC declaration of conformity	119
11.8	Layout plan	121
12	Index	123

1 General information

1.1 Importance of the operating manual

A fundamental requirement for the safe and trouble-free operation of the centrifuge is to be familiar with the fundamental safety instructions and all possible hazards.

The operating manual includes important information concerning the safe operation of the centrifuge.

This operating manual and, in particular, the notes on safety and hazards must be observed by all persons operating the centrifuge.

In addition, the local rules and regulations for the prevention of accidents must be complied with.

1.2 Intended use

Centrifuges are power-driven machines that separate liquids from solid matter, liquid mixtures, or solid mixtures by centrifugal force. They are solely intended for this purpose. Any other use beyond this area of application is regarded as improper use. Sigma Laborzentrifugen GmbH cannot be held liable for any damage resulting from such improper use.

The intended use also includes

- observation of all the notes and instructions included in the operating manual and
- compliance with the care, cleaning, and maintenance instructions.

1.3 Warranty and liability

The warranty and liability are subject to our "General Conditions" that were distributed to the operator upon the conclusion of the contract.

Warranty and liability claims are excluded if they are due to:

- improper use.
- non-compliance with the safety instructions and hazard warnings in the operating manual.
- improper installation, start-up, operation, or maintenance of the centrifuge.

1 General information

1.4 Copyright

The copyright concerning the operating manual remains with Sigma Laborzentrifugen GmbH.

The operating manual is solely intended for the operator and their personnel. It includes instructions and information that must not be

- duplicated,
- distributed, or
- communicated in any other way.

Non-compliance may be prosecuted under criminal law.

1.5 Standards and regulations

EC declaration of conformity
(see chapter 11.7 - "EC declaration of conformity")

1.6 Scope of supply

The centrifuge comprises:

- | | |
|---|-------------------------------------|
| • 1 square spanner, size 8 (door) | Part no. 930 114 |
| • 1 open spanner, size 8/10
(adjustable feet front side) | Part no. 930 015 |
| • 1 open spanner, size 24
(adjustable feet rear side) | Part no. 930 024 |
| • 1 wrench, angulate, size 17/19 (rotor) | Part no. 26 448 |
| • 1 tube wrench (emergency release) | Part no. 930 110 |
| • 1 wrench, hex socket, size 4 (rotor) | Part no. 930 050 |
| • 1 holder for the rotor cover
with hexagon socket head screw (M6x16)
(see chapter 5.3 - "Installation of the cover
holder") | Part no. 28 598
Part no. 964 216 |
| • 1 tube of grease for load bearing bolts | Part no. 70 284 |
| • 1 hose connector for condensate drain (installed) | Part no. 80 415 |

Documentation

Operating manual incl. EC declaration of conformity
(see chapter 11.7 - "EC declaration of conformity")

Accessories

according to your order, our order confirmation, and your delivery note.

2 Layout and mode of operation

2.1 Layout of the centrifuge

2.1.1 Functional and operating elements

- 1 Lid
- 2 User interface (see chapter 6.3.1 - "User interface")
- 3 Mains power switch



Fig. 1: Total view of the centrifuge

- 4 Locks of the front door
- 5 Name plate (see chapter 2.1.2 - "Name plate")



Fig. 2: Right side of the centrifuge

2 Layout and mode of operation

- 6 RS-232 interfach
- 7 Option: barcode scanner



Fig. 3: Left side of the centrifuge

- 8 Mains cable
- 9 Adjustable foot
- 10 Option: serial interface (see chapter 6.3.4 - "Options for data input and output")
- 11 Castor



Fig. 4: Rear view of an air cooled centrifuge

2 Layout and mode of operation

- 12 Cooling water connection (see chapter 5.6 - "Special equipment: water cooling system")



Fig. 5: Rear view of a water cooled centrifuge¹

¹ Not for Sigma 8KB

2 Layout and mode of operation

2.1.2 Name plate

- 1 Manufacturer and registered office
- 2 Type
- 3 Serial number
- 4 Max. speed
- 5 Max. kinetic energy
- 6 Max. density
- 7 Nominal voltage
- 8 Input fuse
- 9 Symbol for special disposal (see chapter 9 - "Disposal")
- 10 CE mark in compliance with the directive 2006/42/EC
- 11 Part number
- 12 Year of manufacture
- 13 Power consumption

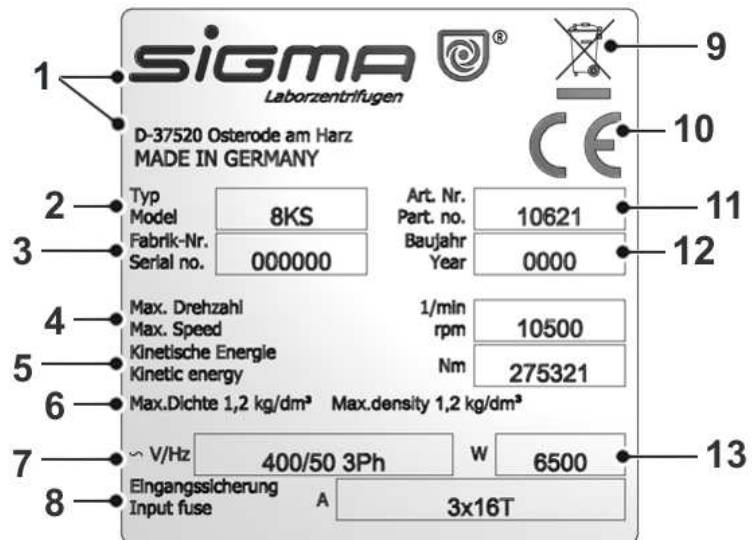


Fig. 6: Example of a name plate

2.1.3 Mains plug

The centrifuge is equipped with a 5-pin Cekon plug (16 A).



Fig. 7: 5-pin Cekon plug, 16A



NOTE

If the centrifuge is operated within a different power supply network, it may be necessary to replace the mains power plug.
 The rotary field is of no importance when replacing the plug.



DANGER

Work on the power supply system must only be performed by certified electricians!

2.2 Mode of operation

2.2.1 Centrifugation principle

Centrifugation is a process for the separation of heterogeneous mixtures of substances (suspensions, emulsions, or gas mixtures) into their components. The mixture of substances, which rotates on a circular path, is subject to centripetal acceleration that is several times greater than the gravitational acceleration.

Centrifuges use the mass inertia inside the rotor chamber for separating the substances. Due to their higher inertia, particles or media with a higher density travel outwards. In doing so, they displace the components with a lower density, which in turn travel towards the centre.

The centripetal acceleration of an object inside a centrifuge, as the effect of centripetal force, depends on the distance between the object and the axis of rotation as well as on the angular velocity. It increases linearly as a function of the distance with regard to the axis of rotation and quadratically as a function of the angular velocity. The bigger the radius in the rotor chamber is and the higher the speed is, the higher the centripetal acceleration is. However, the forces acting on the rotor also increase.

2.2.2 Area of application

Depending on the area of application of the centrifuge and also on the particle size, solids content, and volume throughput of the mixture of substances that is to be centrifuged, there are different types of centrifuges. The areas of application go from household use as a salad spinner or honey separator up to specialised technical applications in the clinical, biological, or biochemical context:

- For numerous clinical examinations, cellular material must be separated from the liquid to be analysed. The normal separation process can be sped up considerably by using laboratory centrifuges.
- In the metal-working industry, centrifuges are used for separating oil from metal cuttings. Dairies use centrifuges in order to separate cow's milk into cream and low-fat milk.
- Particularly big centrifuges are used in the sugar industry for separating the syrup from the crystalline sugar.
- Ultracentrifuges are predominantly used in biology and biochemistry in order to isolate particles, e.g. viruses. They are specifically designed for high speeds up to 500,000 rpm. The rotor moves in a vacuum in order to avoid air friction.

2 Layout and mode of operation

2.2.2.1 Speed, radius, and relative centrifugal force

The acceleration g , which the samples are subject to, can be increased by increasing the radius in the rotor chamber and by increasing the speed. These three parameters are interdependent and linked with each other via the following formula:

$$\text{Relative centrifugal force RCF} = 11.18 \times 10^{-6} \times r \times n^2$$

r = radius in cm

n = speed in rpm

RCF without any dimension

If two values are entered, the third value is determined by way of the stated formula. If, afterwards, the speed or the radius is changed, the resulting relative centrifugal force will be recalculated automatically by the control unit. If the RCF is changed, the speed will be adapted while the specified radius is maintained.

The speed-gravitational-field-diagram provides an overview of the relationship between speed, radius, and RCF (see chapter 11.2 - "Speed-gravitational-field-diagram").

2.2.2.2 Density

The laboratory centrifuge is suitable for the separation of constituents of different densities in mixtures with a maximum density of 1.2 g/cm^3 . All information concerning the speed of rotors and accessories refers to liquids with a density corresponding to this specification. If the density is above this value, the maximum permissible speed of the centrifuge must be reduced based on the following formula:












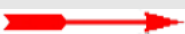




$$n = n_{\max} \times \sqrt{(1.2 / \text{Rho})}$$

Rho = density in g/cm^3

3 Safety

3.1 Marking of the unit

The following symbols are used for all types of centrifuges manufactured by Sigma:

	Dangerous voltage		On (Power)
	Hot surface		Off (Power)
	Caution! Risk of bruising		Name plate (see chapter 2.1.2 - "Name plate")
	Protective earth (ground)		CE mark in compliance with the directive 2006/42/EC
	Earth (ground)		Do not dispose as part of domestic waste
	Unplug the mains plug		GS mark (tested safety; only for Germany)
	Arrow indicating the direction of rotation		NRTL mark (only for the USA and Canada)
	4°C Guarantee		RCM mark (only for Australia)
	Centrifuge with heating system		

i
NOTE

Safety indications on the centrifuge must be kept readable at all times. If necessary, they must be replaced.

i
NOTE

Not all of the symbols/labels are used for this centrifuge type.

3 Safety

3.2 Explanation of the symbols and notes

In this operating manual, the following names and symbols to indicate hazards are used:



DANGER

This symbol stands for a **direct** hazard to the life and health of persons.

Non-observance of these symbols **causes** serious health problems up to life-endangering injuries.



DANGER

This symbol stands for a **direct** hazard to the life and health of persons due to electrical voltage.

Non-observance of these symbols **causes** serious health problems up to life-endangering injuries.



WARNING

This symbol stands for a **potential** hazard to the life and health of persons.

Non-observance of these symbols **can** cause serious health problems up to life-endangering injuries.



CAUTION

This symbol indicates a potentially hazardous situation

Non-observance of these notes can cause minor injuries or damage to property.



NOTE

This symbol indicates important information.

3.3 Responsibility of the operator

The operator is responsible for authorising only qualified personnel to work on the centrifuge (see chapter 3.4 - "Operating personnel").

The areas of responsibility of the personnel concerning the operation, maintenance, and care of the unit must be clearly defined.

The safety-conscious work of the personnel in compliance with the operating manual and the relevant EC and national health and safety regulations as well as with the accident prevention regulations must be checked at regular intervals (e.g. every month).

Under the international rules for health and safety at work, the operator is obliged to:

- take measures in order to prevent all danger to life or health during work.
- ensure that centrifuges are operated properly and entirely as intended (see chapter 1.2 - "Intended use").
- take protective measures against fire and explosion when working with hazardous substances.
- take measures for the safe opening of centrifuges.

3.4 Operating personnel

Persons operating the unit must

- be familiar with the fundamental regulations concerning workplace safety and accident prevention
- have read and understood this operating manual (and in particular the safety sections and warning notes) and confirmed this with their signature.

3.5 Informal safety instructions

- This operating manual is a part of the product.
- The operating manual must be kept at the location of use of the centrifuge. Ensure that it is accessible at all times.
- The operating manual must be handed over to any subsequent owner or operator of the centrifuge.
- Any changes made must be added to the operating manual.
- In addition to the operating manual, the general and local rules and regulations concerning the prevention of accidents and the protection of the environment must also be supplied.
- Safety and danger indications on the centrifuge must be kept readable at all times. If necessary, they must be replaced.

3 Safety

3.6 Safety instructions

3.6.1 Electrical safety

To reduce the risk of electrical shock, the centrifuge uses a three-wire electrical cord and plug to connect the equipment to earth-ground. To preserve this safety feature:



DANGER

- Ensure that the wall socket is properly wired and grounded.
- Check that the mains voltage agrees with the nominal voltage listed on the name plate.
- Do not place vessels containing liquid on the centrifuge lid or within the safety distance of 30 cm around the centrifuge. Spilled liquids may get into the centrifuge and damage electrical or mechanical components.
- Work on the power supply system must only be performed by certified electricians.
- Inspect the electrical equipment of the unit regularly. Defects such as loose or burnt cables must be eliminated immediately.

3.6.2 Mechanical safety

In order to ensure the safe operation of the centrifuge, observe the following:



WARNING

- Do not open the lid when the rotor is in motion!
- Do not reach into the rotor chamber when the rotor is in motion!
- Do not use the centrifuge if it was installed incorrectly.
- Do not use the centrifuge without panels.
- Do not use the centrifuge if the rotors and inserts show signs of corrosion or other defects.
- Only use the centrifuge with rotors and accessories that have been approved by the manufacturer. In case of doubt, contact the manufacturer (see chapter 7.3 - "Service contact").
- Do not hold your fingers between the lid and the housing when closing the lid. Risk of crushing!
- Defective lid relieving devices could cause the centrifuge lid to fall (contact the service department, if necessary). Risk of crushing!
- Do not hit or move the centrifuge during its operation.
- Do not lean against or rest on the centrifuge during its operation.
- Do not spin any substances that could damage the material of the rotors and buckets of the centrifuge in any way. Highly corrosive substances, for example, damage the material and affect the mechanical strength of the rotors and buckets.
- Stop the centrifuge immediately in the event of a malfunction. Eliminate the malfunction (see chapter 7 - "Malfunctions and error correction") or inform the service department of the manufacturer (see chapter 7.3 - "Service contact").
- Ensure that all repairs are performed only by authorised and specialised personnel.

**WARNING**

- Prior to any start-up, check the centrifuge, rotor, and accessories for signs of damage that can be discerned from the outside. Special attention must be paid to all of the rubber parts (e.g. motor cover, lid seal, and adapters) in terms of visible structural changes. Defective parts must be replaced immediately.
- Open the centrifuge when it is not in use so that moisture can evaporate.

3.6.3 Fire prevention

**DANGER**

- Do not spin explosive or inflammable substances.
- Do not use the centrifuge within hazardous locations.

3.6.4 Chemical and biological safety

If pathogenic, toxic, or radioactive samples are intended to be used in the centrifuge, it is in the responsibility of the user to ensure that all necessary safety regulations, guidelines, precautions, and practices are adhered to accordingly.

**DANGER**

- Infectious, toxic, pathogenic, and radioactive substances may only be used in special, certified containment systems with a bio-seal in order to prevent the material from being released.
- Take suitable precautions for your own safety if there is a risk of toxic, radioactive, or pathogenic contamination
- Materials that chemically react with each other with a high level of energy are prohibited.

**WARNING**

- Keep informed about local measures to avoid harmful emissions (depending on the substances to be centrifuged).
- Protective clothing is not required for the operation of the centrifuge. The materials to be centrifuged may, however, require special safety measures (e.g. centrifugation of infectious, toxic, radioactive, or pathogenic substances).

3 Safety

3.6.5 Safety instructions for centrifugation

For safe operation, observe the following before starting the centrifuge:



WARNING

- Ensure that the centrifuge was set up properly (see chapter 5 - "Set-up and connection").
- Maintain a safety distance of at least 30 cm (12 inches) around the centrifuge.
- Do not store any dangerous goods in the centrifuge area.
- Do not stay in the safety area longer than what is absolutely necessary for the operation of the centrifuge.
- Only use the centrifuge with rotors and accessories that have been approved by the manufacturer. We explicitly warn against the use of equipment of poor quality. Breaking glass or bursting vessels can cause dangerous imbalances at high speeds
- Ensure that rotor and buckets are correctly fitted (see chapter 6.2.2.1 - "Installation of the rotor").
- Observe the instructions on the installation of accessories (see chapter 6.2.2.3 - "Installation of accessories").
- The rotor must be loaded axial symmetrically at equal weights.
- If liquids with a density $> 1.2 \text{ g/cm}^3$ are used, reduce the speed (see chapter 2.2.2.2 - "Density").
- Do not use the centrifuge if the rotor is loaded asymmetrically.
- Do not use the centrifuge with tubes that are excessively long.



WARNING

3.6.6 Resistance of plastics

Chemical influences have a strong effect on the polymeric chains of plastics, and, therefore, on their physical properties. Plastic parts can be damaged if solvents, acids, or alkaline solutions are used.



NOTE

- Refer to the resistance data (see chapter 11.5 - "Resistance data")!

3.6.7 Rotors and accessories

3.6.7.1 Lifting and carrying rotors



All swing-out rotors applicable for this centrifuge and the angle rotor 12510 weigh more than 18 kg.

- Always lift the rotors with a lifting device or with a sufficient number of people helping you.

3.6.7.2 Service life

The rotors and accessories have a limited service life.



- Perform regular checks (at least once per month) for safety reasons!
- Pay special attention to changes, such as corrosion, cracks, material abrasion etc.
- After 10 years, they must be inspected by the manufacturer.
- After 50,000 cycles, the rotor must be scrapped for reasons of safety.
- If other data concerning the service life are engraved on the rotor or bucket, these data shall apply accordingly. For example, a bucket with the engraving "max cycles = 10,000" has a service life of 10,000 cycles, and a rotor with the engraving "Exp.Date 02/15" must be scrapped in February 2015 at the latest (see figure).



Fig. 8: Different service life – engraving on the bucket/rotor



- Refer to the table of rotors and accessories with a different service life (see chapter 11.4 - "Table of the service life of rotors and accessories")!

3 Safety

3.7 Safety devices

3.7.1 Lid lock device

The centrifuge can only be started when the lid is properly closed. The electrical lock must be locked. The lid can only be opened when the rotor has stopped. If the lid is opened by way of the emergency release system during operation, the centrifuge will immediately switch off and decelerate brakeless. If the lid is open, the drive is completely separated from the mains power supply, i.e. the centrifuge cannot be started (see chapter 7.1.1 - "Emergency lid release").

3.7.2 Standstill monitoring system

Opening of the centrifuge lid is only possible if the rotor is at a standstill. This standstill is checked by the microprocessor.

3.7.3 System check

An internal system check monitors the data transfer and sensor signals with regard to plausibility. Errors are detected with extreme sensitivity and displayed as error messages in a dialog box (see chapter 7.2 - "Table of error codes").

3.7.4 Earth conductor check

For the earth conductor check, there is an equipotential bonding screw on the rear panel of the centrifuge (see chapter 2.1.1 - "Functional and operating elements"). An earth conductor check can be carried out by authorized and specialized personnel using a suitable measuring instrument. Please contact the service department (see chapter 7.3 - "Service contact").

3.7.5 Imbalance monitoring system

A dialog box may pop up or emit a sound signal in order to indicate that the centrifuge is in the inadmissible imbalance range. If the rotor is loaded unevenly, the drive will be switched off in the acceleration phase or during the run.

3.7.6 Temperature monitoring system

If the temperature inside the rotor chamber rises above 50°C, the drive system will be switched off automatically. The centrifuge cannot be restarted until it has cooled.

3.7.7 Rotor monitoring system

When a rotor number and, if applicable, a bucket number are selected, the computer will automatically check whether the entered speed or the entered gravitational field are permissible for the selected rotor.

3.8 Measures in the event of hazards and accidents



DANGER

- If an emergency arises, switch off the centrifuge immediately!
- If in doubt, call the emergency doctor!

3.9 Remaining hazards

The unit was built state- of- the- art and according to the accepted safety rules. However, danger to life and limb of the operator, or of third parties, or impairments of the unit or other material assets cannot be completely excluded when the unit is being used.

- Use the unit only for the purpose that it was originally intended for (see chapter 1.2 - "Intended use").
- Use the unit only if it is in a perfect running state.
- Immediately eliminate any problems that can affect safety.

4 Storage and transport

4 Storage and transport

4.1 Dimensions and weight

	Sigma 8KS / 8KBS
Height (mm):	980
Height with open lid (mm):	1 690
Width (mm):	810
Depth (mm):	910
Weight (kg):	420

4.2 Storage conditions

The centrifuge can be stored in its original packaging for up to a year.

- Store the centrifuge only in dry rooms.
- The permissible storage temperature is between -20°C and +60°C.
- If you would like to store it for more than one year, or if you intend to ship it overseas, please contact the manufacturer.

4.3 Notes on transport



NOTE

The centrifuge weighs approx. 420 kg!

- Always lift the centrifuge with a lifting device (e.g. a fork lift).
- When lifting the centrifuge, always reach under the centrifuge from the rear side.
- For transport, use suitable packaging and, if at all possible, the original packaging(see chapter 4.4 - "Packaging").

Special equipment - water cooling system

Prior to transporting the centrifuge, the water circuit must be drained in order to avoid damage, e.g. due to freezing. Only specialised personnel are authorised to perform these tasks. Consultation with the manufacturer is obligatory!

4.4 Packaging



WARNING

When unloading the centrifuge, wear safety shoes in order to avoid injuries caused by the rolling centrifuge!

The centrifuge is packaged in a wooden crate.


- Remove the upper cover.
- Remove the packing material and all of the accessories.

Unloading the centrifuge with a forklift truck

One side panel of the crate has metal fittings (see the following illustration, detail A).

- Remove the two lower screws of each of the metal fittings and then the four outer screws in the upper half of the side panel (see the following illustration). Then, remove the side panel.

1 Metal fittings

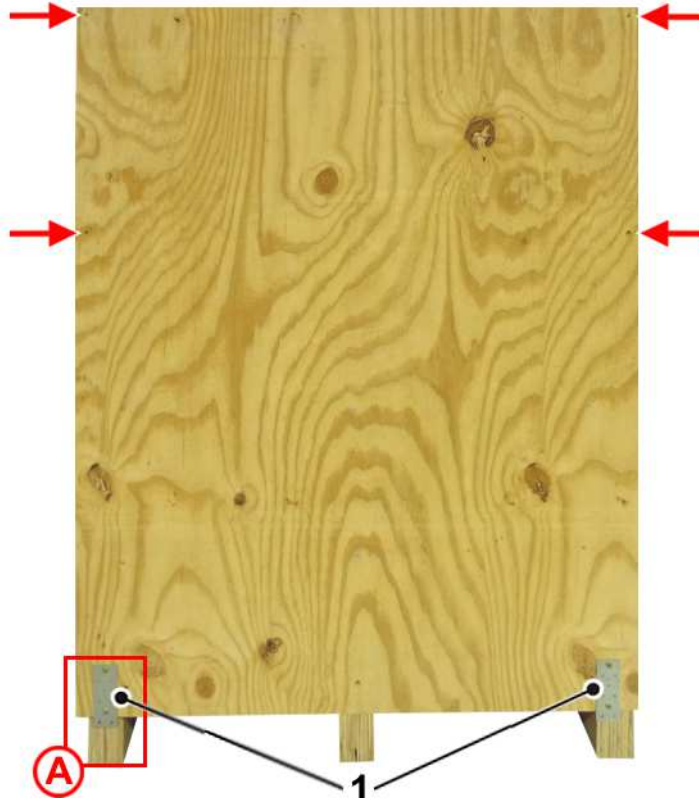

4 Storage and transport

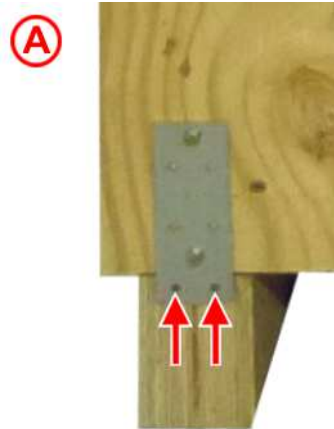


Fig. 9: Positions of the screws for removing the side panel

- Lift the centrifuge off the bottom of the crate by way of a forklift truck.



NOTE

The centrifuge weighs approximately 420 kg!

Unloading the centrifuge with a hand pallet jack

If no forklift truck is available, one side panel of the wooden crate can be used as a ramp so that the centrifuge can be unloaded from the crate by way of a hand pallet jack (ground clearance 85 mm max.).

One side panel of the crate has metal fittings (see the illustration above, detail A).

- Remove the two lower screws of the two metal fittings and then the four outer screws in the upper half of the side panel (see the illustration above). Then, remove the side panel. The metal fittings remain on the side panel.



NOTE

Do not remove the screws in the middle of the side panel. They hold the substructure of the ramp in place.

4 Storage and transport

- Position the side panel in front of the bottom of the crate as a ramp.



Fig. 10: Positioning of the ramp

- Secure the ramp on the bottom of the crate by way of the two metal fittings and the screws that have been removed beforehand (see the following illustration).


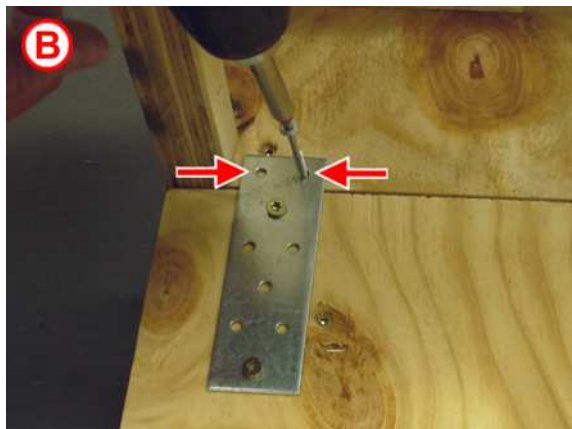


Fig. 11: Positions of the screws for securing the ramp (here: left side)

4 Storage and transport

- Move the pallet jack over the ramp and under the centrifuge and lift the centrifuge off the bottom of the crate.
- Unload the centrifuge carefully via the ramp. While doing so secure the centrifuge with a sufficient number of persons in order to prevent it from slipping off to the side.

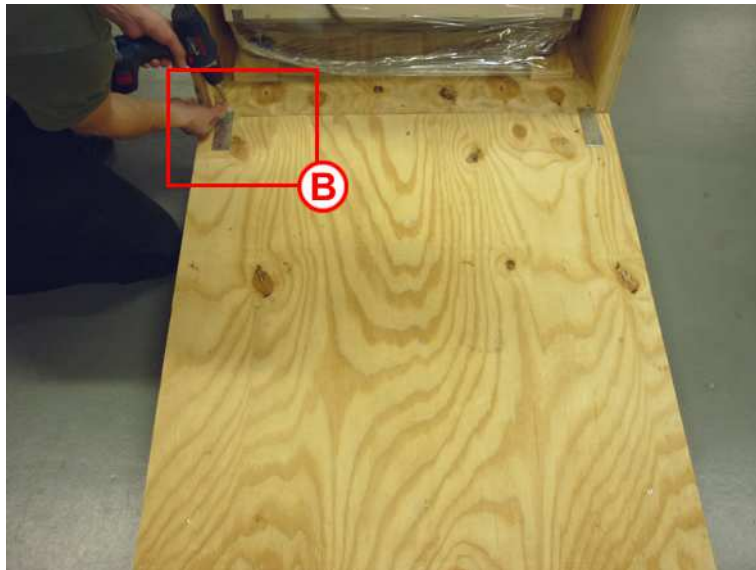

CAUTION

Do not unload the centrifuge via the ramp on its castors since, in this case, the centrifuge would be difficult to steer and could slip off the ramp.

Unloading the centrifuge with an electric pallet jack



NOTE

Due to its design and dimensions, an electric pallet jack is not suitable for unloading the centrifuge.

- Retain the packaging for any possible future transport of the centrifuge.

5 Set-up and connection

5.1 Installation site

Operate the centrifuge only in closed and dry rooms.

All the energy supplied to the centrifuge is converted into heat and emitted to ambient air.

- Air cooled centrifuges: Ensure sufficient ventilation.



NOTE

The centrifuge 8KS / SKBS with an air-cooled compressor should not be set up with its left side against a wall, since otherwise the hot air, which is emitted out the back, will be drawn in again as fresh air for cooling. As a result, the unit will switch off due to overheating (see chapter 7.2 - "Table of error codes").

- Water cooled centrifuges: Provide sufficient water throughput.
- Do not position the centrifuge near heat generators.
- Avoid direct sunlight (UV radiation).
- During transport from cold to warmer places, condensational water will collect inside the centrifuge. Allow sufficient time for drying (min. 24 h) before using the centrifuge again.

5.2 Set-up and alignment of the centrifuge



DANGER

Ensure that the centrifuge is disconnected from the power supply during the set-up and alignment.



NOTE

Pay attention to the layout plan (see chapter 11.8 - "Layout plan")!

5 Set-up and connection

5.2.1 Set-up

- Transport the centrifuge as closely as possible to the installation site with a lifting device (e.g. forklift).
- Set the centrifuge down.
- To place the centrifuge on the castors, open the front door with the supplied square spanner (part no. 930 114) by turning it clockwise by 90°. Open the front door to the left in order to access the two locking screws located at the front.
- Loosen the two hexagon lock nuts with the open spanner, size 24 (part no. 930 024), and screw in the locking screws with the open spanner, size 8/10 (part no. 930015), in the anticlockwise direction up to the stop until the adjustment feet are completely relieved of the load.
- Loosen the lock nuts of the two locking screws on the back from the side and from below with the open spanner, size 24 (part no. 930 024). Then, screw in the locking screws anticlockwise with the open spanner, size 13/14 (part no. 930 013), until the adjustment feet are completely relieved of the load.
- Transport the centrifuge on the castors to the installation site.



CAUTION

Do not use the lid handle to move the centrifuge because it could break off!



CAUTION

The castors of the centrifuge are made of steel without any plastic coating. Damage to the surface of the floor cannot be excluded.

5.2.2 Alignment



DANGER

For reasons of safety, the centrifuge must not be operated while it is set up on its castors!



NOTE

The centrifuge must be set-up stably and horizontally at the installation site.

- Set the centrifuge on its adjustable feet (see chapter 5.2.1 - "Set-up" in reverse order). To do so, turn all of the feet clockwise by hand until they touch the ground. Then, perform approximately two more clockwise turns with the open spanner until the castors are suspended in the air.

5 Set-up and connection

- Open the lid with the emergency lid release system. To do so, lift off the plugs on the left side (e.g. with a screwdriver) and unlock the lid locks by turning them clockwise with the supplied square spanner (see chapter 7.1.1 - "Emergency lid release").
- Align the centrifuge with a spirit level in two directions.
- Tighten the four lock nuts and close the front door.
- Plug in the mains power plug.

- 1 Locking screw
- 2 Lock nut



Fig. 12: Alignment at the front

- 3 Adjusting foot
- 4 Castor



Fig. 13: Alignment at the back

5 Set-up and connection

5.3 Installation of the cover holder

The cover of the swing-out rotor 11805 with a windshield can be placed in the lid of the centrifuge when loading or unloading the centrifuge. In addition, it is possible to install a cover holder on the right side of the centrifuge. The required accessories are included in the scope of supply:

- 1 Holder for the rotor cover
- 2 Hexagon socket head screw (M6x16)

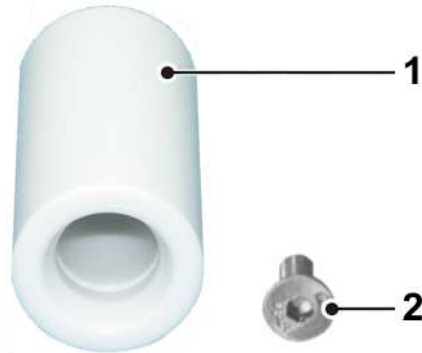


Fig. 14: Cover holder accessories

Installation

There is a hole for the cover holder in the middle of the upper edge on the right-hand side of the centrifuge.

- Insert the hexagon socket head screw through the hole of the cover holder and fasten the cover holder to the hole on the right-hand side of the centrifuge by way of the hexagon socket key (included in the scope of supply).

- 3 Holder in the centrifuge lid
- 4 External holder



Fig. 15: Storage options for the windshield cover

5.4 Power supply

5.4.1 Type of connection



DANGER

The operating voltage on the name plate must correspond to the local supply voltage!

SIGMA laboratory centrifuges are units of protection class I. The centrifuges of this type have a five-wire power cord (2.5 m) with a 5-channel Cekon-plug (16 A). Behind the front door there are three fuses with a rocker switch.

5.4.2 Customer-provided fuses

Typically, the centrifuge must be protected with 16 Amp L or B fuses that are to be provided by the customer.

5.4.3 Power isolating device

The electrical installation of the building must include a power isolating device in the form of a switch or circuit breaker. This device must be located near the centrifuge. Furthermore, it must be easily accessible for the operator and marked as the power isolating device for the centrifuge.

5.5 Condensate drain

The condensate drain is used to drain the condensate off that forms in the rotor chamber during the centrifugation. It consists of a hose with a plastic valve. This hose runs from the rotor chamber to a point on the left behind the front door of the centrifuge.



WARNING

- Do not open the condensate drain unless the rotor is at a complete standstill.

Draining the condensate off

- Unplug the mains power plug.
- Open the front door with the supplied plastic wrench (order no. 930 114).
- Connect the supplied hose connector and drain the condensate.
- Disconnect the hose connector by pressing the unlocking button.

5 Set-up and connection

i

NOTE

At the installation site, the condensate drain can be laid to the outside (see figure).

- 1 Drain outside of the centrifuge
- 2 Hose (part no. 80415)
- 3 Hose connector (unblocks the drain when connected)
- 4 Unlocking button
- 5 Plug (to be removed prior to the installation)
- 6 Centrifuge door (front, left side)
- 7 Clamp
- 8 Quick coupling
- 9 Lock nut (fasten by hand)
- 10 Condensate – rotor chamber

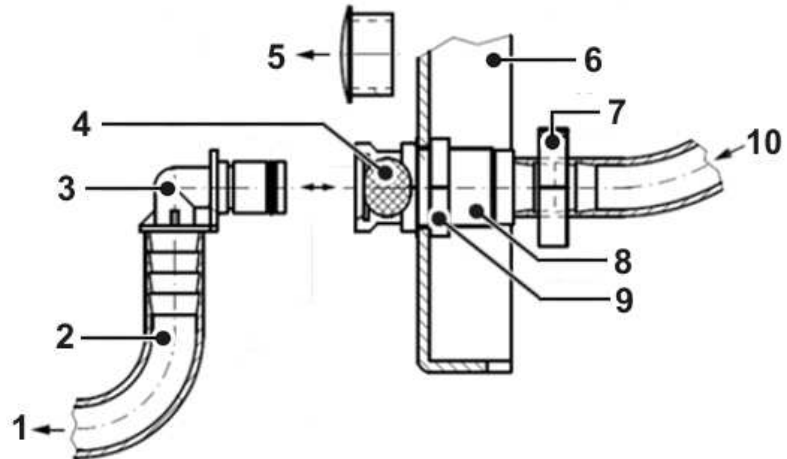


Fig. 16 Laying the condensate drain to the outside

5.6 Special equipment: water cooling system

Centrifuges with water cooling are equipped with a special refrigeration unit with a refrigerant (CFC-free), modified electronic system, and a special software version. The housing of the GMP version is completely closed. All of the other units have vent slots for the motor and electronic power system.

Operating conditions

- The centrifuge must only be operated with media-neutral water of drinking water quality (hardness <8) or with another suitable refrigerant. We recommend using a filter (to be provided by the customer).
- The operating pressure of the water must be between 1.5 and 5 bar. We recommend using a pressure reducer (to be provided by the customer).
- The maximum flow rate depends on the operating conditions (e.g. speed and rotor temperature). The minimum cooling water consumption at maximum power is approximately 5 l/min. At a standstill, the water supply is stopped by a solenoid valve. In the case of the GMP version (ref. no. 10857), the valve is opened in two steps:
 1. reduced flow, only for the internal air cooler
 2. unlimited flow, when the compressors are running.
- The maximum temperature at the water inlet is 20°C. The lower the water temperature, the more efficient the cooling.
- The centrifuge is equipped with a ¾-inch water connector .

- 1 Connector for the cooling water inlet
- 2 Connector for the cooling water outlet



Fig. 17: Connections for the water cooling system

6 Using the centrifuge

6 Using the centrifuge

6.1 Initial start-up



DANGER

- Before the initial start-up, please ensure that your centrifuge is properly set up and installed (see chapter 5 - "Set-up and connection").

6.2 Switching the centrifuge on

- Press the mains power switch.

The display then illuminates. The centrifuge is ready for operation.

6.2.1 Opening and closing the lid

The lid can be opened if the centrifuge is at a standstill and if the lid key is illuminated.

- Press the lid key in order to open the lid.

The centrifuge cannot be started if the lid is opened.

- To close, press with both hands slightly on the lid until both electrical locks are locked.



WARNING

Do not place your fingers between the lid and the housing when closing the lid. Risk of crushing!

6.2.2 Installation of rotors and accessories



WARNING

All swing-out rotors applicable for this centrifuge and the angle rotor 12510 weigh more than 18 kg.

- Always lift the rotors with a lifting device or with a sufficient number of people helping you.

6.2.2.1 Installation of the rotor

- Open the centrifuge lid by pressing the lid key.
- Unscrew the rotor tie-down screw from the motor shaft (counter-clockwise).
- Lower the rotor with its central bore straight down onto the motor shaft.
- Turn the adapter end of the motor shaft anti-clockwise with the supplied hexagon socket key (part no. 930 050) and tighten the internal thread of the rotor with 20 Nm with the rotor wrench (size 17/19, part no. 930 018).



WARNING

After frequent use, the rotor tie-down screw must be loosened by some turns, the rotor has to be lifted and fastened again. This must be done once a day or after 20 cycles. This ensures a proper connection between the rotor and the motor shaft.

- Follow the safety instructions and hazard warnings (see chapter 3 - "Safety")!

6.2.2.2 Installation of angle rotors with a hermetically sealed lid

- Lower the rotor onto the motor shaft.
- Turn the adapter end of the motor shaft anti-clockwise with the supplied hexagon socket key (part no. 930 050) and tighten the internal thread of the rotor with 20 Nm with the rotor wrench (size 17/19, part no. 930 018).



NOTE

- Slightly grease the rotor and lid seals after cleaning.

- Screw the rotor cover onto the rotor and tighten it with the supplied tool no. 17985 by hand.
- The rotors can be installed or removed with a closed lid after loosening the rotor tie-down screw.
- Follow the safety instructions and hazard warnings (see chapter 3 - "Safety")!



WARNING

The lid screw serves for the fastening of the lid onto the rotor only, not for the fastening of the rotor onto the motor shaft!

6 Using the centrifuge

6.2.2.3 Installation of accessories

- Only use inserts that are suitable for the rotor (see chapter 11.1 - "Range of accessories").
- All four buckets of the swing-out rotor need to be installed when spinning.
- Always load the axial symmetrical inserts/buckets of the rotors with the same accessories and fill to avoid imbalance.

Centrifugation with different tube sizes

Working with different tube sizes is possible. In this case, however, it is very important that the inserts are loaded axial symmetrically.

correct



incorrect

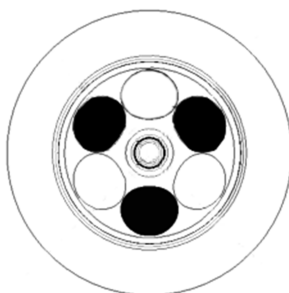


Fig. 18: Permissible and impermissible loading of a swing-out rotor with different tube sizes (example illustrations)

Centrifugation with low capacity

- Install the tubes axial symmetrically so that the buckets and their inserts are loaded evenly.
- It is not permissible to load angle rotors on only one axis.

correct



incorrect

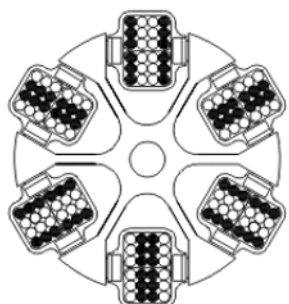
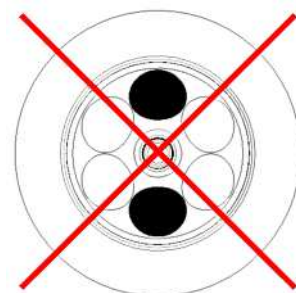


Fig. 19: Permissible and impermissible loading of an angle rotor and a swing-out rotor (example illustration)

6.2.2.4 Adapters

In order to ensure easy handling, even if vessels of various sizes are used, carrier systems were developed.

- Load the opposite adapters with the same number of vessels and with the same weights in order to avoid imbalance.
- If all of the compartments of a carrier are not used, the buckets must be loaded evenly. Loading the edges of a bucket only is not permissible.

6.2.2.5 Tubes

- Load the tubes outside of the centrifuge. Liquids in the buckets or multiple carriers cause corrosion .
- Fill the tubes carefully and arrange them according to their weight. Imbalances result in the excessive wear of the bearings.
- In high-speed angle rotors, the vessels must be filled up to their useful volume (= the volume stated for the vessel). If the vessels are only partially filled, they will deform. This may result in leaks at the seals that may become loose.
- When using glass tubes, the maximum value of 4,000 x g must not be exceeded (except special high-strength glass tubes; please refer to the information provided by the manufacturer).
- Follow the safety instructions and hazard warnings (see chapter 3 - "Safety")!

6 Using the centrifuge

6.2.2.6 Blood bag systems

- All six places on the rotor must be loaded with buckets.
- It is required to fill just two opposite buckets must be filled with one adapter for blood bags and two blood bag systems each.
- The opposite buckets, including the filled blood bag systems, must have an equal weight. If the number of blood bags is uneven, a balance adapter (part no. 17765 for 13865, part no. 17766 for 13866) must be used. For taring, several balance weights are available (see chapter 11.1 - "Range of accessories").
- The blood bags must be put into opposite buckets in a mirror-inverted way (see figure).
- In both adapter compartments, the main blood bag must be situated towards the middle. The opposite bucket must be loaded correspondingly (see figure, item 1).
- When using smaller bag systems or in the case of incompletely filled blood bags, adapters (e.g. part no. 17750) must be inserted together with the blood bag systems. This will help to avoid any slipping of the bags, which could result in an imbalance.

- 1 Proper loading
2 Possible loading
3 Improper loading

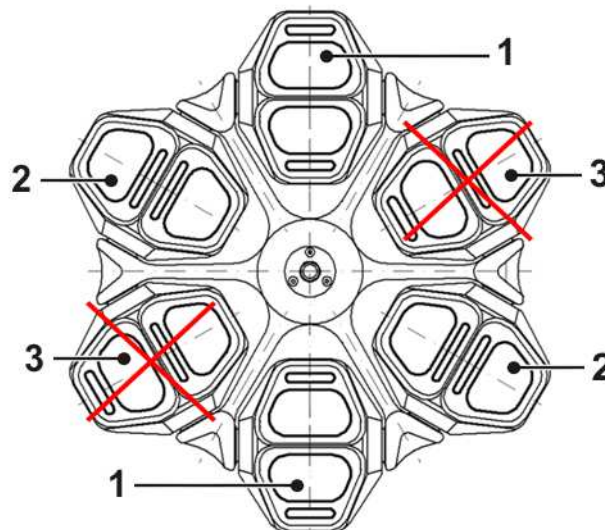


Fig. 20: Loading of blood bag systems



WARNING

The service life of the adapters for blood bags no. 13865 and 13866 is limited. Refer to the table of rotors and accessories with a different service life (see chapter 11.4 - "Table of the service life of rotors and accessories")!



NOTE

At maximum speed, the plastic adapters 13865 and 13866 must not be used at a temperature above 25°C. If they are used at a considerably lower speed (up to 2,500 rpm), the temperature can be higher than 25°C.

6.3 Control system "Spincontrol S"

6.3.1 User interface

The centrifuge is operated via three buttons with integrated light-emitting diodes and one function knob. The display is divided into several different fields. The various functions of the system can be called up by pressing and turning the function knob.

- 1 Start button
- 2 Display
- 3 Function knob
- 4 Stop button
- 5 Lid button

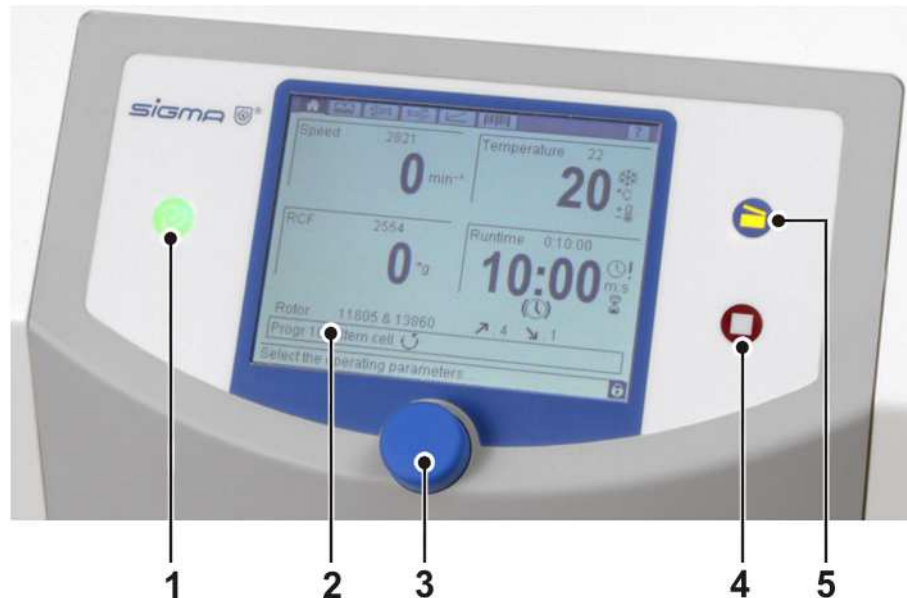


Fig. 21: User interface of the Spincontrol S control system

Display

The centrifuge display has the following fields:

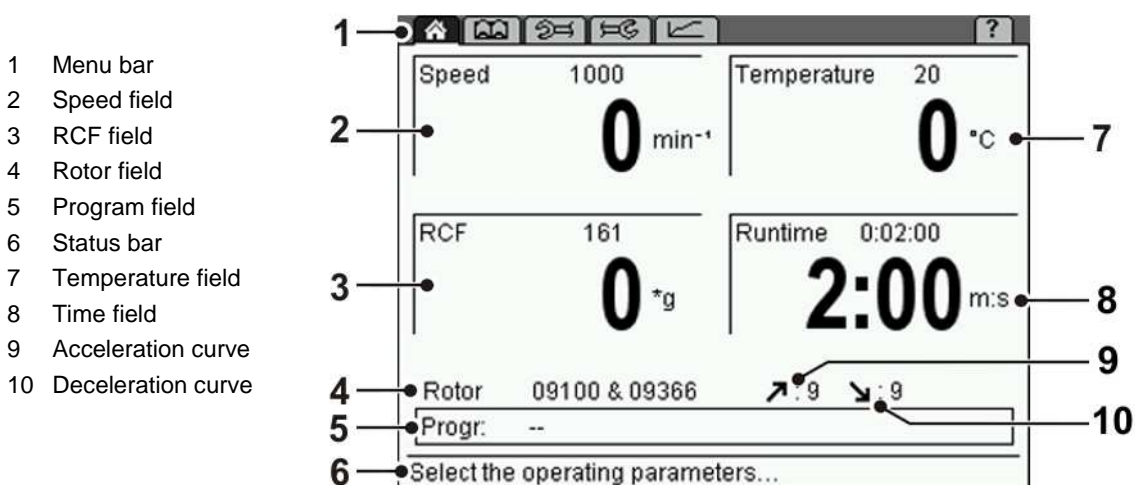


Fig. 22: Display of the Spincontrol S control system

6 Using the centrifuge

6.3.2 Manual mode

6.3.2.1 Starting a centrifugation run

The centrifuge is ready for operation when the start key is illuminated.

- Press the start key in order to start a centrifugation run.

6.3.2.2 Interrupting a centrifugation run

- Press the stop key in order to interrupt a centrifugation run. The centrifugation run will be terminated prematurely.

Quick stop

- Press the stop key for more than three seconds.

The centrifuge decelerates with the maximum deceleration curve.

After a quick stop, the centrifuge lid must be opened before a new centrifugation run can be started.

A quick stop can also be triggered during a normal deceleration, e.g. in order to speed up the deceleration.

When a quick stop is triggered, "Quick stop" will be displayed in the speed field.



NOTE

A quick stop can be performed even if the centrifuge is blocked against unauthorised use.

6.3.2.3 Interrupting a deceleration process

- Press the start key during a deceleration process in order to interrupt it and to restart the centrifuge.

6.3.2.4 Selection, display, and modification of data

The "Standard" menu is displayed.

- Turn the function knob in order to select a field. The selected field is inverted.
- Press the function knob. The display starts to flash and the modification mode is active.
- Turn the function knob in order to modify the set value of the selected field.
- Press the function knob again to confirm the entry and to quit the modification mode.

6.3.2.5 Standard menu

The "Standard" menu is symbolised by the icon "⚙️" on the menu bar and it is displayed a few seconds after the centrifuge has been switched on. In this menu, the parameters of a centrifugation run can be displayed and modified.

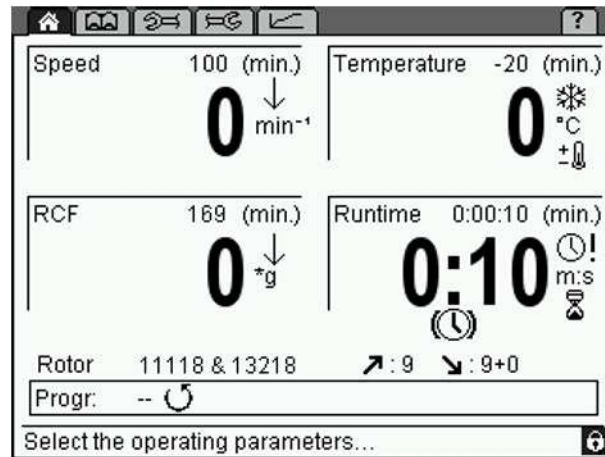


Fig. 23: Standard menu; here shown with all of the possible symbols

Speed

In the upper section of the field, the set speed of the centrifuge is displayed. The actual speed is displayed below this value. The values are stated in revolutions per minute ($\text{min}^{-1} = \text{rpm}$) and depend on the RCF values (see chapter 2.2.2.1 - "Speed, radius, and relative centrifugal force"). The maximum speed values depend on the rotor that is used.

Relative centrifugal force (RCF)

The relative centrifugal force is the acceleration that the sample is subjected to during the centrifugation run. The set value of this parameter is displayed in the upper section of this field, with the actual value shown below. The values are stated in g (gravitational acceleration) and they depend on the speed values (see chapter 2.2.2.1 - "Speed, radius, and relative centrifugal force"). The maximum RCF values depend on the rotor that is used.

Temperature

The set temperature is displayed in the upper section of the field, with the current rotor chamber temperature shown below. Temperatures between -20 °C and $+40\text{ °C}$ can be preselected.




NOTE

The centrifuge is not equipped with an active heater. As a result, temperatures above room temperature depend on the air friction of the running motor.

6 Using the centrifuge

Runtime

The set runtime is displayed in the upper section of this field, with the remaining runtime shown below. The runtime is defined as the period from the start of the centrifuge to the beginning of the deceleration phase. The maximum value is 99 h 59 min 59 sec.

In the "Setup" menu , it can be specified that the runtime is not to be started until the set speed is reached (see chapter 6.3.2.8 - "Setup menu"). In this case, the symbol "!" appears in the runtime field.

Continuous run

During the continuous run, the runtime of the centrifuge is unlimited and must be stopped manually. The centrifuge accelerates during the continuous run until the set speed is reached.

- Select the field "Runtime" and confirm the selection. The display flashes when it is activated.
- Turn the function knob from the time 0:00:10 anti-clockwise or from the time 99:59:59 clockwise. "Infinite" will be displayed. After the start of the centrifuge, the elapsed time will be displayed.
- Deactivate the continuous run by pressing the stop button or by entering a specific runtime.

Short run

A short run can be started if no run is active.

- Keep the start button pressed during the short run.

During the short run, the centrifuge accelerates with acceleration curve 9 (maximum acceleration) until the maximum speed of the rotor is reached. The runtime is counted and in the speed field the message "Short run" flashes.

When the start button is released, the centrifuge decelerates with the maximum deceleration curve to a standstill.



NOTE

The parameters speed, RCF, temperature, and runtime can be changed during the centrifugation.



NOTE

If the centrifuge is locked with level 02 or higher, it is not possible to start a short run.

Rotor: rotor selection list

This field shows the rotor that is currently being used.

- Select the field "Rotor" and confirm the selection. A list with all of the possible rotors without buckets is displayed.
- Select the desired rotor.
 - If an angle rotor is selected, additional information concerning this rotor will be displayed.
 - In the case of swing-out rotors, a list with all of the possible rotor/bucket combinations will be displayed. Select an item from the list so that the additional information concerning the combination will be displayed.
- Press the function knob in order to accept the data.

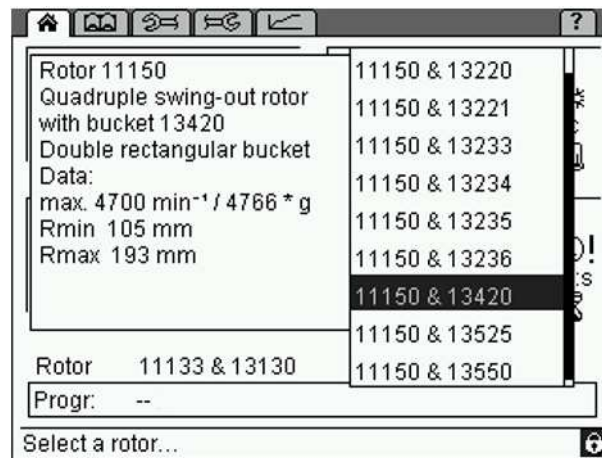


Fig. 24: Rotor selection list, here with potential rotor-bucket combinations and additional information

Automatic rotor identification

The centrifuge automatically identifies the rotor that is currently being used.

- If the system identifies a different rotor than the one that is set and if there are no different buckets for this rotor, the rotor input will be adapted automatically. The system will not display a message.
- If the system identifies a different rotor than the one that is set, and if there are different rotor/bucket combinations for this rotor, the system will automatically identify the correct rotor and select the rotor/bucket combination with the lowest speed. The system will display a corresponding message so that the combination can be adapted manually.
- If the system cannot identify the rotor, a message will be displayed. The rotor cannot be used in the centrifuge.

This prevents the maximum permissible speed from being exceeded.

6 Using the centrifuge

Acceleration ↗

This function is used to select an acceleration curve. One can select a linear rise (curves 0-9) or a quadratic rise (curves 10-19). The acceleration curves 20-29 can be programmed as desired (see chapter 11.3 - "Acceleration and deceleration curves").

Deceleration (brake) ↘

This function is used to select a curve that decelerates the centrifuge to a standstill. Deceleration curves are inverted images of the acceleration curves and are labelled with identical numbers. Deceleration curve no. 0 represents a brakeless deceleration.

Progr.: program list

This field in the "Standard" menu shows the program that is currently loaded. When the field is selected, the program list is displayed (for information on how to work with the programs, please see chapter 6.3.3 - "Program mode").

The program "RAPID_TEMP" (see below) cannot be deleted.

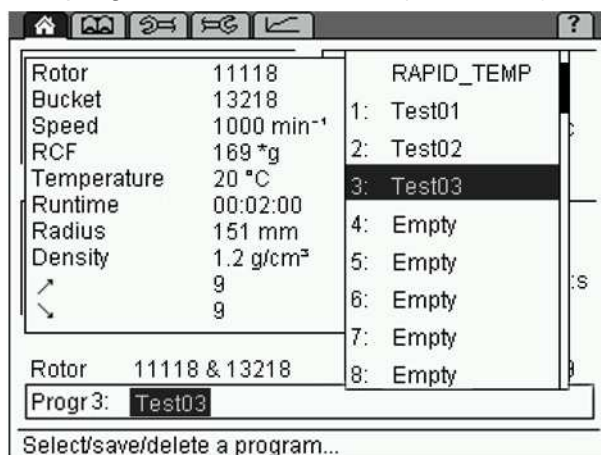



Fig. 25: Program list

Program "RAPID_TEMP"

Precooling or preheating at a standstill may distort the measurement results and cause the increased wear of the mechanical components. This is why the centrifuge has a special program that precools or preheats the rotor chamber rapidly to a preset value under defined conditions.

- Select the option "Progr" in the "Standard" menu  and confirm the selection. The program list will be displayed.
- Select the program "RAPID_TEMP" on the program list and confirm the selection. The display shows $\frac{1}{3}$ of the maximum rotor speed and the corresponding RCF value. The deceleration (brake) and acceleration curves correspond to curve 9 and the runtime field indicates "infinite" (continuous run).

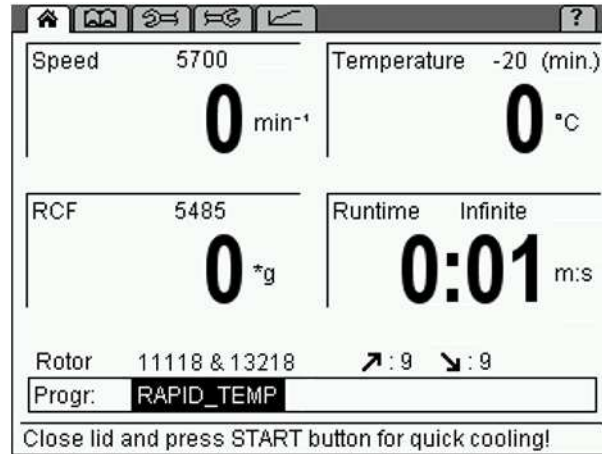


Fig. 26: Program "RAPID_TEMP"

i
NOTE

The program will only be loaded if the actual temperature is above the set temperature.

- Press the start button in order to start the rapid cooling run.
- During the rapid cooling run, the set temperature can be modified within the range below the actual temperature.

The current status of the program will be displayed on the status bar.

The "RAPID_TEMP" program will be stopped under the following conditions:

- The set value is reached. In this case, the "RAPID_TEMP" program stops with a sound signal and the standstill cooling system will be activated.
- The stop button is pressed. The "RAPID_TEMP" program will be stopped prematurely. No message will be issued when the set temperature is reached.
- A parameter is changed (except for the temperature) or another input is made. In this case, the "RAPID_TEMP" program will be aborted. No message will be issued when the set temperature is reached.

After the stop, the previous program will be reloaded or the changed parameters will be adopted as the new settings.

i
NOTE

The automatic lid opening function is suppressed after a rapid cooling run in order to prevent the system from reheating.

i
NOTE

The delta T temperature monitoring system (see chapter 6.3.2.7 - "Parameters menu", Process) remains inactive as long as the "RAPID_TEMP" program is active.

6 Using the centrifuge



If the "RAPID_TEMP" program is used, the temperature of the unloaded aluminium bucket will be displayed. If samples, which have not been precooled, are placed into the buckets, the displayed temperature will deviate from the actual sample temperature.

Progress indicator

The progress indicator provides a quick overview of the remaining runtime of the running centrifugation run. For this purpose, a green progress bar and percentage value are displayed in the program field.

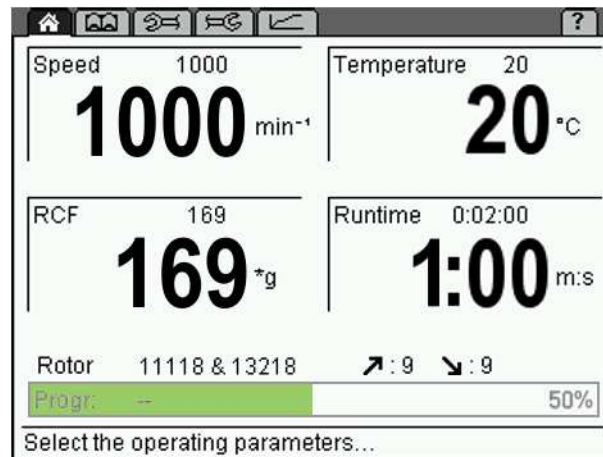


Fig. 27: Progress indicator during a centrifugation run

After the completion of the centrifugation run, the progress indicator remains at 100% until


- the lid is opened,
- a parameter is changed,
- a program is loaded, saved, or deleted, or
- a new centrifugation run is started.

6.3.2.6 Process library menu

The process library can be opened via the "📖" symbol on the menu bar. It provides the user with an overview of all of the stored programs as processes.

The processes are listed with their name, RCF, runtime, and temperature (exception: the "RAPID_TEMP" program will not be listed).

The order of the processes corresponds to their storage locations on the program list. Empty storage locations will not be displayed. If more than 11 programs have been stored, the user can scroll through the list.

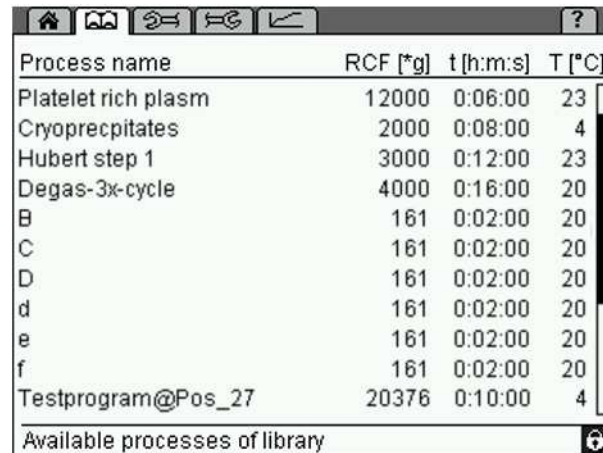

Process name	RCF [*g]	t [h:m:s]	T [°C]
Platelet rich plasm	12000	0:06:00	23
Cryoprecipitates	2000	0:08:00	4
Hubert step 1	3000	0:12:00	23
Degas-3x-cycle	4000	0:16:00	20
B	161	0:02:00	20
C	161	0:02:00	20
D	161	0:02:00	20
d	161	0:02:00	20
e	161	0:02:00	20
f	161	0:02:00	20
Testprogram@Pos_27	20376	0:10:00	4

Available processes of library

Fig. 28: Process library menu

Loading a process

- In the process library, select the desired process by turning the function knob. Press the function knob in order to confirm the selection.

The process will be loaded and the "Standard" menu 🏠 will be displayed.

Starting a process

- In the process library, select the desired process by turning the function knob. Then, press the start button.

The process will be loaded and started. The "Standard" menu 🏠 will be displayed.

6 Using the centrifuge

6.3.2.7 Parameters menu

The "Parameters" menu is symbolised by the "⚙️" symbol on the menu bar. It is used to specify various conditions for the centrifugation. These conditions are used to monitor the process and to control access to the centrifuge.

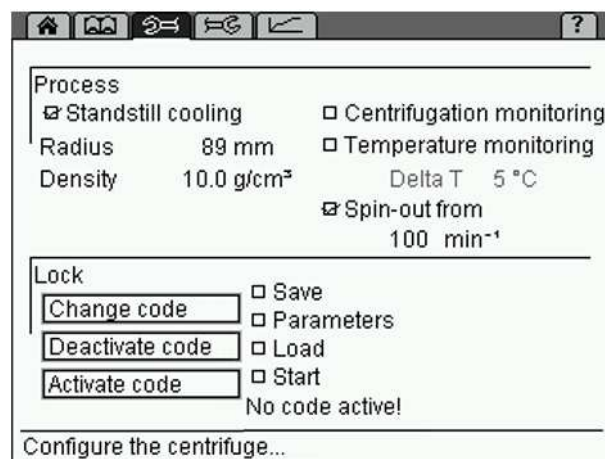


Fig. 29: Menu "Parameters"

Process

Standstill cooling

Depending on the substances to be centrifuged, it may make sense to precool the centrifuge. The precooling prevents the cooled samples in the uncooled centrifuge from heating up to an inadmissible temperature level.

If the standstill cooling function is activated, the centrifuge starts to precool after it is switched on. In the runtime field, the symbol "❄️" is displayed. The lid must be closed.



CAUTION

Unmoved air in the rotor chamber leads to an incorrect measuring and control behaviour and causes the compressor to freeze over. At temperatures below 0°C, aqueous liquids will freeze, thereby preventing sedimentation.

- Do not cool the rotor below 0°C at a standstill.



NOTE

The centrifuge is equipped with the "RAPID_TEMP" program. This program is used to precool the rotor chamber quickly under defined conditions (see chapter 6.3.2.5 - "Standard menu", Program list).

Radius

The radius determines the value of the relative centrifugal force (RCF) that the sample is subjected to. Normally, the maximum RCF value is displayed. If the value is reduced manually, a downward facing arrow "↓" will be displayed in the RCF field.

Density

This setting is useful for glass vessels. If the density of the liquid to be centrifuged is higher than 1.2 g/cm³, the value must be adapted manually in order to prevent the glass vessel from breaking. This will reduce the maximum possible final speed (see chapter 2.2.2.2 - "Density"). The reduction will be represented by a downward facing arrow "↓" in the speed field. Values between 1.2 and 10.0 g/cm³ are possible.

Centrifugation monitoring

The centrifugation monitoring function enables the continuous monitoring of the speed and runtime parameters during the centrifugation.

- Activate the centrifugation monitoring function by clicking.



NOTE

If the function is activated during a centrifugation run, the monitoring process will not be started until the start of the next centrifugation run.

The centrifugation monitoring function compares the speed values of the current run with the reference values that are stored in the control unit. After every run, it issues a corresponding message.

The runtime is considered faulty if the centrifugation run had to be stopped prematurely.

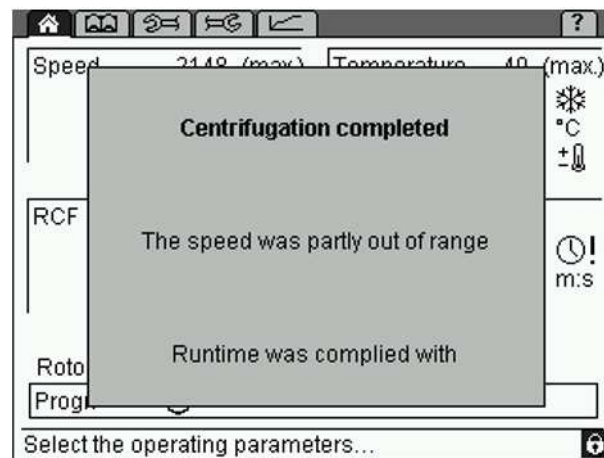


Fig. 30: Example of a centrifugation monitoring message



NOTE

If the centrifugation monitoring system is used in combination with free acceleration or deceleration curves, unjustified error messages may result in certain cases.

6 Using the centrifuge

Temperature monitoring

The control system includes a temperature monitoring function. If the actual temperature difference with regard to the set value deviates from the set temperature difference, the centrifuge will stop and an error message will be issued.

- Activate the temperature monitoring function by clicking.
- Adjust the desired "Delta T" limit value in steps of 1°C or 1°F and confirm the setting.

If the temperature monitoring function is active, the symbol "±" is displayed in the temperature field of the "Standard" menu.

In this case, the centrifuge can only be started if the current temperature is between the preset temperature in the "Standard" menu and the "Delta T" limit value.

If the temperature leaves this range during the centrifugation run, an error message will be displayed and the centrifuge will be decelerated to a complete standstill.



NOTE

The delta T temperature monitoring system remains inactive as long as the "RAPID_TEMP" program is active (see , "Program "RAPID_TEMP"").

Spin-out from...

If this function is active, the brake will be deactivated when the actual speed is below the set speed. As a result, the rotor will spin out in a brakeless manner.



NOTE

A spin-out, in particular with heavy rotors and at high speeds, can take a lot of time! (Depending on the rotor and load, the speed will be reduced by approximately 0.5 to 1 rpm per second.)

If the spin-out is active, "+0" is displayed next to the deceleration curve.

- The spin-out can be interrupted by a quick stop or by restarting the centrifuge.

Lock

In order to prevent any unauthorised use of the centrifuge, the following functions can be blocked:

- Saving of programs (level 01)
- Changing of parameters, short run (level 02)
- Loading of programs (level 03)
- Start button (level 04)

Blocking a function

- Select the function that is to be blocked. The lower levels will also be automatically selected (if, for example, the "Parameters" function is selected, the "Save" function will also be selected).
- Select the button "Activate code".
- Enter a four-digit code and confirm the entry.

The lock is now active. The symbol "🔒" will be displayed in the status line and the lock level will be indicated.

If changes are made after a function has been blocked, the system will ask for the code prior to executing the change.

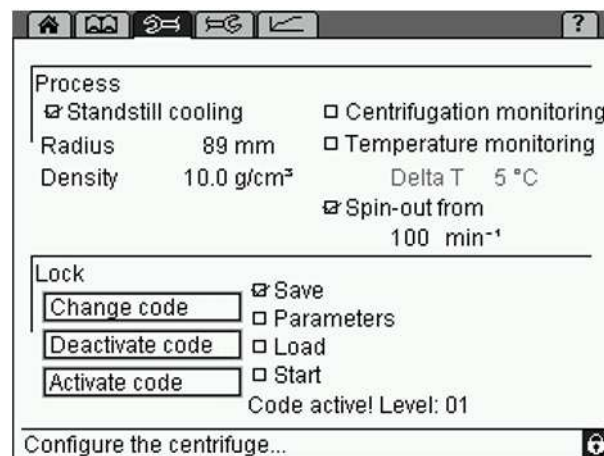


Fig. 31: "Parameters" menu with an active lock (level 01)

Unblocking a function

- Select the button "Deactivate code".
- Enter the code and confirm the entry.

The lock is now deactivated.

Changing the code

- Select the button "Change code".
- Enter the old code and confirm the entry.
- Enter the new code.
- For safety reasons, the code must be entered a second time.

The code is now changed.

6 Using the centrifuge

6.3.2.8 Setup menu

The "Setup" menu is symbolised by the "⚙️" symbol on the menu bar. It is used to perform basic settings concerning the control system of the centrifuge. It enables the optimum adaptation of the centrifuge to its specific area of application.

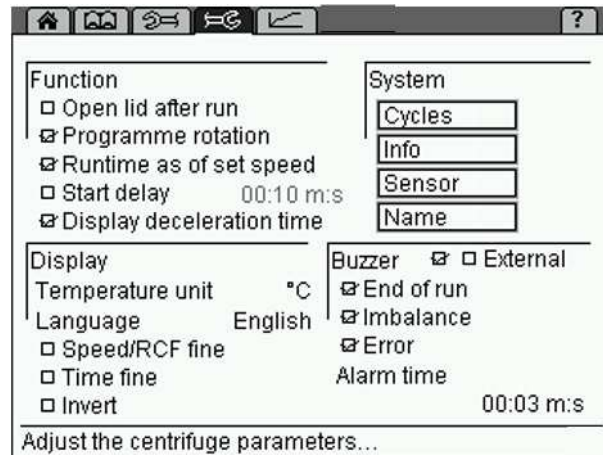


Fig. 32: Menu "Setup"

Function

Open lid after run

The automatic lid opening function ensures that the lid opens when the rotor has stopped.



CAUTION

When the lid is open, the cooling is not active. The sample may warm up!

Program rotation

See chapter 6.3.3.5 - "Automatic program rotation".

Runtime as of set speed

If this function is active, the runtime will not be measured until the set speed is reached. In the "Standard" menu, the symbol "⌚!" will be displayed in the runtime field.



NOTE

In the program mode, this function can be saved separately for every program. In this case, the symbol "⌚!" will be displayed under the runtime of the detailed program description.

Start delay

If the start delay function is active, the centrifuge will not start until the preset time has elapsed. The symbol "⌚" will be displayed in the runtime field.

Display deceleration time

In the activated mode, the deceleration time will be displayed instead of the runtime during and after the deceleration process. Below the time display, the symbol "⏸" will be displayed. During a deceleration process, the symbol flashes. Once the deceleration is complete, it is displayed in a permanent manner.

Display

Temperature unit

The temperature can be displayed in °C (Celsius) or °F (Fahrenheit).

Language

The control system can be used in various language versions.

If a language is selected by mistake, it can be changed on any screen as follows:

- Press and hold the stop button.
- Turn the function knob one notch to the left and then one notch to the right.
- Release the stop button. The "Language" window will be displayed.
- Select the desired language.

Speed/RCF fine and Time fine

This menu item can be used to preselect the set speed in steps of 1 rpm (instead of 100 rpm), the RCF value in steps of 1 x g (instead of 10 x g), and the set time in steps of 1 min or 1 sec (instead of 10 min or 10 sec).



NOTE

Regardless of the fine adjustment, the step size increases when the function knob is turned quickly.

Invert

If this function is activated, the display switches from the standard setting with a bright background and dark writing to a dark background with bright writing.

6 Using the centrifuge

System

Cycles

This field shows the number of cycles as well as the runtime of the rotor and buckets that are used.

Info

This item provides information on the software versions that are used in this centrifuge.

Sensor

The sensor mode is reserved for the service personnel.



NOTE

Values can neither be entered nor changed in the menus "Cycles", "Info", and "Sensor".

Name

In this field, an identification will be assigned to the centrifuge.

- The letters and characters can be entered when the cursor flashes in the text field. Turn the function knob in order to select a character and press it to confirm the selection. Then, press the knob again in order to enter the next character. Pressing the arrow button ← will delete the last character. The maximum number of characters is 19.
- When the name is complete, select the option "Accept" and confirm it.

Buzzer (signal)

With this function, a warning sound signal can be selected for

- the end of a centrifugation run,
- an imbalance message,
- an error message.

The duration of the warning signal can be specified.

External

This function is only available if the centrifuge is equipped with the option for the input and output of data (external signal, floating switch) (see chapter 6.3.4 - "Options for data input and output").

6.3.2.9 Curve menu

This menu is used to create and edit customised acceleration and deceleration curves (see chapter 11.3 - "Acceleration and deceleration curves"). It is symbolised by the "⏮" symbol on the menu bar.

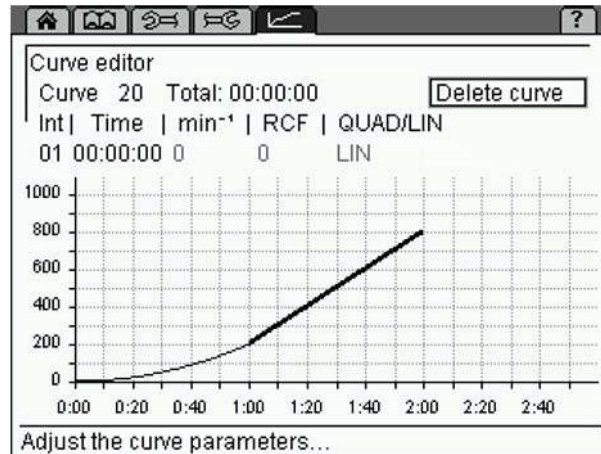


Fig. 33: Menu "Curve"

Creating or changing an acceleration curve



NOTE

During a centrifugation run, curves can only be displayed. They cannot be changed or deleted.

- Select the "Curve" menu. The curve editor will be displayed.
- Select a curve number between 20 and 29 in the input field "Curve". If the curve number is already used, the stored curve will be displayed.
- The input field "Int" is used to specify the interval number of the process. Up to ten intervals can be entered for a curve.
- Enter the interval time of the current interval into the input field "Time". While doing so, certain restrictions must be taken into consideration (see below).
- Enter the desired acceleration under "min⁻¹" (rpm) or "RCF". While doing so, certain restrictions must be taken into consideration (see below). The values are interdependent.
- In the first interval, "QUAD/LIN" can be used to select a linear or quadratic rise. All of the other intervals are linear.

The field "Total" shows the total runtime of the process. The maximum total runtime of a curve depends on the slope of the curve and on the final speed of the rotor.



NOTE

Only the last curve interval can be changed retroactively.

6 Using the centrifuge

Restrictions

- Acceleration and deceleration curves can include intervals with a positive slope as well as intervals with a negative slope and also intervals with 0 slope.
- The slope of the curve intervals can be 1 min⁻¹/sec (rpm/sec) minimum and 1000 min⁻¹/sec (rpm/sec) maximum.
- Quadratic curve intervals are only possible between 0 and 1000 min⁻¹ (rpm) maximum. If a final speed > 1000 min⁻¹ (rpm) is selected, this interval will automatically become linear above 1000 min⁻¹ (rpm).
- The possible runtime results from the maximum possible speed (depending on the rotor) and from the slope limitation.

Example 1: Start speed 0 min⁻¹ (rpm), final speed 100 min⁻¹ (rpm), runtime 1 hour not possible, since the necessary slope < 0.03 min⁻¹ (rpm), which means that it is beyond the defined range.

Example 2: Start speed 0 min⁻¹ (rpm), final speed 15000 min⁻¹ (rpm), runtime 10 sec. not possible, since the necessary slope is 1500 min⁻¹ (rpm), which means that it is beyond the defined range.

6.3.2.10 Option: Barcode menu



NOTE

In order to connect a barcode scanner, the centrifuge must be equipped with a second RS232 interface and a corresponding circuit board.

If the centrifuge is equipped with the barcode scanner function, the symbol "|||||" will be displayed on the menu bar. In order to use this function, a PC and a barcode scanner must be connected to the RS232 interface.

The "Barcode" menu shows all of the codes that are stored in the control system. Data structure 020 is displayed in the left-hand area of the screen, and data structure 001 is in the right-hand area. Free storage locations are symbolised by symbols without any code.



Fig. 34: Barcode menu

With the barcode scanner, the barcodes can be scanned directly at the centrifuge during the loading process. The centrifugation data can be enquired cyclically (automatically) by the PC. Once a centrifugation run is complete, the PC can detect this and save the result (centrifugation data with barcodes).

Barcode commands

Setting a barcode: "setbarcode"

Every code must be set separately. This is done automatically by way of the scanner. Barcodes will only be accepted under the following conditions:

- There is no active centrifugation run.
- There are not any barcodes stored in the system which have already been used for a centrifugation run.
- The barcode does not already exist.
- It is an ISBT 128 barcode with the data structure 001 (donation identification number) and there are fewer than 12 of these codes stored in the system.
- It is an ISBT 128 barcode with the data structure 020 (staff member identification number) and there are not any barcodes with this data structure stored in the system.

If the barcode is accepted, the data of the last centrifugation run will be automatically deleted and the "Barcode" menu will be displayed.



NOTE

In the event of a mains power failure, all the codes that have been scanned in so far will be lost and must be scanned again.

Reading out the stored barcodes: "getbarcodes"

All of the codes will be output in one response. This leads to the following possible response formats:

- No codes are stored
→ syntax "Barcodes none\r\n"
- Only codes with the data structure 001 are stored
→ syntax "Barcodes abc, def,...\r\n" (abc, def,... = barcodes)
- Only codes with the data structure 020 are stored
→ syntax "Barcodes abc (staff)\r\n" (abc = barcode)
- Codes with the data structures 001 and 020 are stored
→ syntax: "Barcodes abc (staff), def,...\r\n" (abc, def,... = barcodes)

6 Using the centrifuge

Deleting the stored barcodes: "deletebarcodes"

This command is used to delete all of the stored barcodes.
The "Barcode" menu will be updated.



NOTE

This command will not be accepted unless the centrifuge is at a standstill.

Reading out the data of the last centrifugation run: "getlastrun"

The response to this command includes information concerning the parameters and status of the last centrifugation run (for the exact format of the data see chapter 11.6 - "Serial Control Interface").

If barcodes have been stored, they will also be output as follows:

- Syntax data structure 001:
"Barcode;x;abc\r\n" (x = serial number starting at 1, abc = barcode)
- Syntax data structure 020:
"Barcode;Staff Member ID;abc\r\n" (abc = barcode)

The enquiry can be performed at any time (before, during, and after the centrifugation run). Depending on the time of the enquiry, some of the results may still be unknown.

6.3.2.11 Help menu

The help function is symbolised by the "?" symbol on the menu bar. It provides a short description of the control elements of the selected option.

Activating and deactivating the help function

- Select the question mark on the menu bar and press the function knob.
- Quit the help function by selecting the question mark and by pressing the function knob again.

Parameters can still be changed when the help function is activated.



Fig. 35: Menu "Help"

6.3.2.12 Changing the contrast

To change the contrast:

- Press and hold the stop button and turn the function knob one notch to the left. A dialog box will be displayed once the stop button is released.
- Adjust the contrast of the centrifuge display and confirm the change.

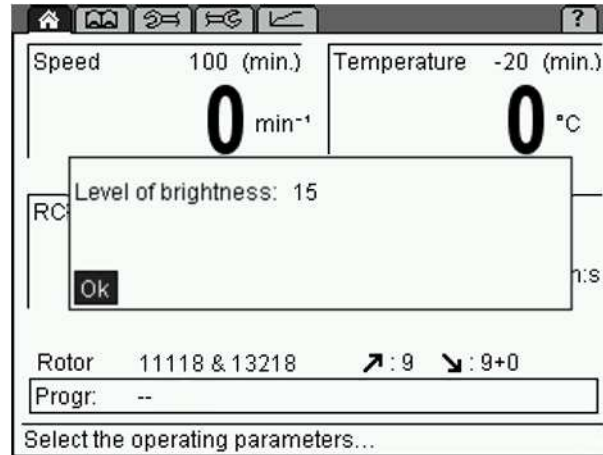


Fig. 36: Dialog box for changing the contrast


6.3.3 Program mode

A program contains all of the data that are required for a centrifuge run. Certain sedimentation results can be repeated under identical conditions. Programs can be loaded, executed, edited, and deleted when the centrifuge is at a standstill.

A maximum of 60 programs can be stored under the numbers 1-60. The program "RAPID_TEMP" does not occupy any storage location and cannot be deleted. It is used to bring the centrifuge to a specific temperature without any vessels.



"--" means that the values that are currently set are not a stored program.

The programs can be protected against unauthorised use, modification, or deletion with the aid of a code (see chapter 6.3.2.7 - "Parameters menu").

Stored programs are listed in the "Process library" menu .

6 Using the centrifuge

6.3.3.1 Saving a program

- Enter the parameters that are to be included in the program.
- Select the option "Progr" in the "Standard" menu  and confirm the selection. The program list will be displayed.
- Select a storage location from the program list.
- Save the program under the desired name. The letters and characters can be entered when the cursor flashes in the text field.
 - Turn the function knob in order to select a letter and press it to confirm the selection. Then, the next character can be selected.
 - Pressing the arrow button  will delete the last character.
- When the program name is complete, select "OK" and confirm it. The program will be saved and the "Standard" menu will be displayed.

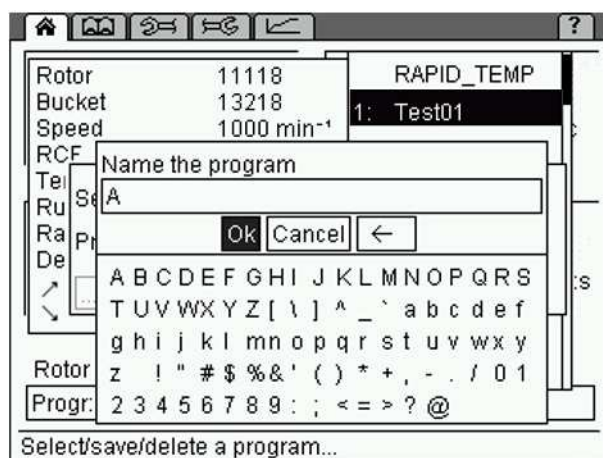




Fig. 37: Assignment of a program name prior to saving the program

6.3.3.2 Loading a program

- Select the option "Progr" from the "Standard" menu  and confirm the selection by pressing the function knob. The program list will be displayed.
 - Select the desired program from the list and confirm the selection by pressing the function knob.
- Or:
- Open the process library , select the desired program, and confirm the selection by pressing the function knob.
- The program will be loaded and the "Standard" menu will be displayed.

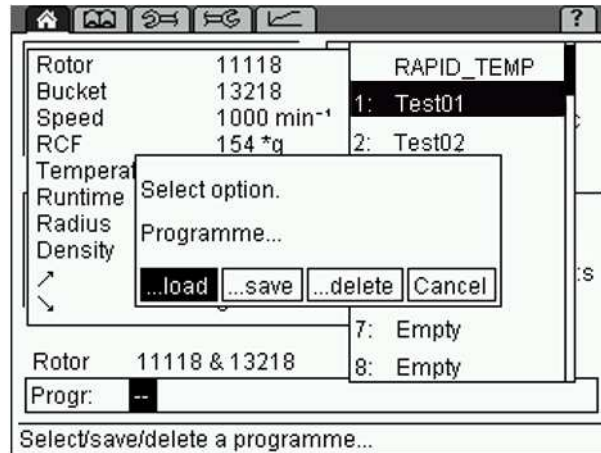




Fig. 38: Loading a program

6.3.3.3 Executing a program


- Select the option "Progr" from the "Standard" menu  and confirm the selection by pressing the function knob. The program list will be displayed.
- Select the desired program from the list and confirm the selection by pressing the function knob.
- Press the start button.

Or:

- Open the process library , select the desired program, and press the start button.

The program will be executed and the "Standard" menu will be displayed.

6.3.3.4 Deleting a program

- Select the option "Progr" in the "Standard" menu  and confirm the selection. The program list will be displayed.
- Select the program that is to be deleted.
- Select the option "Delete" and confirm it.

The program will be deleted and the "Standard" menu will be displayed.

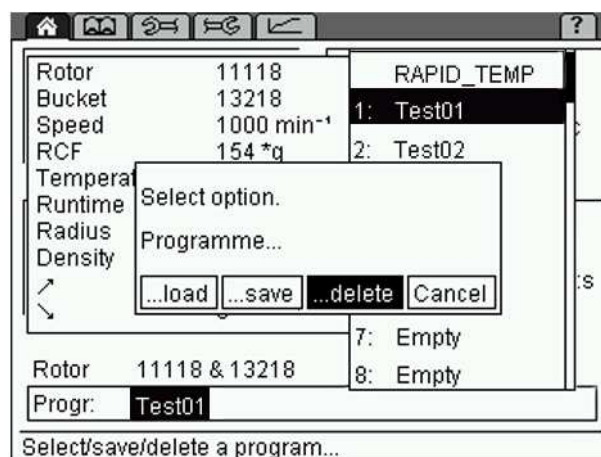



Fig. 39: Deleting a program

6 Using the centrifuge

6.3.3.5 Automatic program rotation

With the automatic program rotation, several programs can be executed directly one after the other.

- Activate the "Program rotation" function in the "Setup" menu .

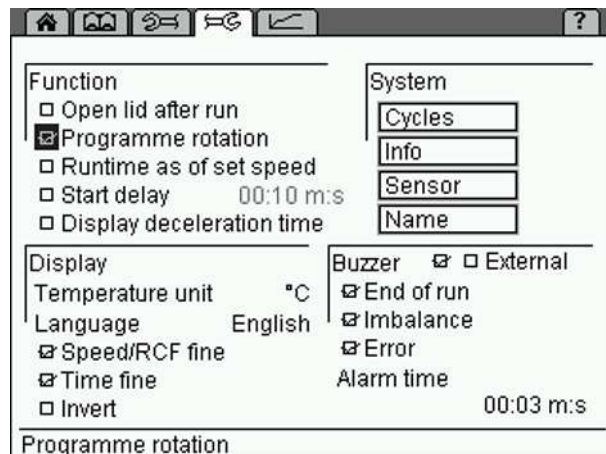


Fig. 40: Program rotation function

When a program is loaded while the program rotation function is active, this program will be used as the start program for the rotation. After the completion of the program, the next program on the program list will be loaded automatically. The rotation continues up to the next empty storage location and then restarts from the beginning (see the following illustration).

Example 1: Loading of Test04

Rotation: Test04, Test05, Test06, Test04,...

Example 2: Loading of Test05

Rotation: Test05, Test06, Test05,...

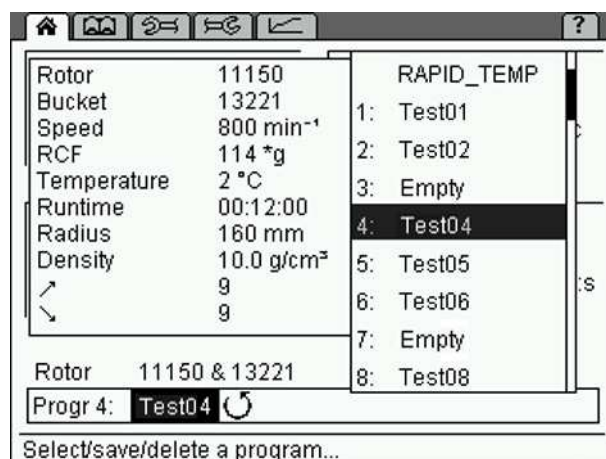


Fig. 41: Automatic program rotation

While the program rotation function is active, the arrow "↻" is displayed in the program line in the "Standard" menu.

6.3.4 Options for data input and output

- Connection for a serial interface (depending on the model, partly standard).
- External signal active DC 24 V, 0.5 A max. (part no. 17701)
- Floating switch AC 250 V max., 6 A (part no. 17702)
- Connection of a barcode scanner via a barcode extension board (see chapter 6.3.2.10 - "Option: Barcode menu")

6.4 Switching the centrifuge off

- Open the centrifuge when it is not in use so moisture can evaporate.
- Switch the centrifuge off by pressing the mains power switch.

7 Malfunctions and error correction

7 Malfunctions and error correction

7.1 General malfunctions

Malfunctions are indicated by a dialog box. If the acoustic signal is activated, it sounds when the error message is displayed.

- Eliminate the source of the problem (see table below).
- Acknowledge the error messages by pressing the lid key.



NOTE

Error messages can be eliminated by pressing the lid key. The error itself will not be eliminated, but the centrifuge can be operated again.

Type of error	Possible reason	Correction
No indication on the display	No power in the mains supply	Check fuse in the mains supply
	Power cord is not plugged in	Plug in power cord correctly
	Fuses have tripped	Reactivate temperature fuse (see chapter 5.4.1 - "Type of connection")
	Mains power switch off	Switch mains power switch on
Centrifuge cannot be started: start key LED is not illuminated	Several	Power off/on. If the error occurs again, contact service
Centrifuge cannot be started: lid key LED flashes	The lid lock is not closed correctly	Open and close the lid. If the error occurs again, contact service
Centrifuge decelerates during operation	Brief mains power failure	Press start key in order to restart the centrifuge
	System error	Power off/on. If the error occurs again, contact service
Centrifuge decelerates during operation, imbalance dialog box is displayed	<ul style="list-style-type: none"> – Improper loading – Centrifuge is inclined – Drive problem – Centrifuge was moved during run 	Balance load and restart the centrifuge. If the error occurs again, contact service(see chapter 7.1.1 - "Emergency lid release")
	– Ungreased load- bearing bolts	Clean and grease load- bearing bolts
Lid cannot be opened	Lid lock has not released	Unlock the lid manually and contact service
	Lid seal sticks	Clean the lid seal and apply talcum powder
Temperature value cannot be reached (only for refrigerated centrifuges)	Condenser dirty (only air-cooled units)	Clean the condenser. If the error occurs again, contact service

7.1.1 Emergency lid release

In the event of a power failure, it is possible to manually open the centrifuge lid.

- Switch off the mains power switch and disconnect the power cord from the socket.
- Remove the plugs (see figure, item 1) at the right side panel, e.g. with a screw driver.

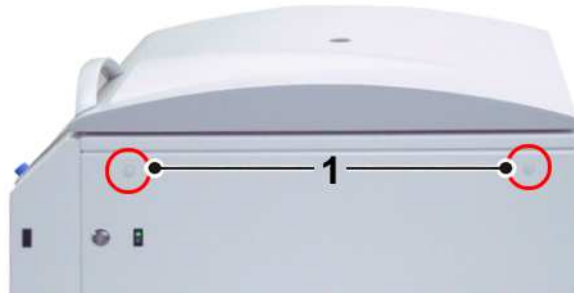


Fig. 42: Position of the openings for the emergency lid release

- Insert the supplied tube wrench (part no. 930 110) horizontally into the hole. The key will be guided through a funnel-shaped tube to the shaft of the lid lock motor.

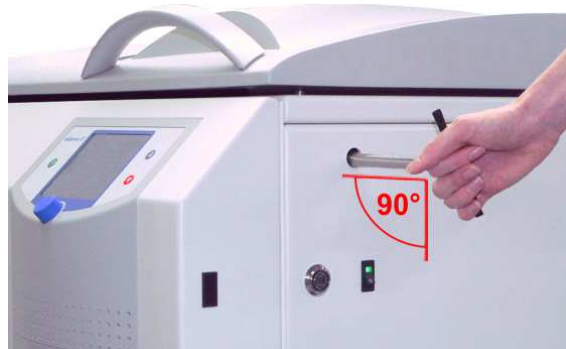


Fig. 43: The emergency lid release key must be inserted horizontally.

- Unlock the motorised lid locks as follows:
 - Turn the left lid lock anti-clockwise.
 - Turn the right lid lock clockwise.
- Then, reinsert the plugs.



WARNING

The lid may only be unlocked and opened when the rotor is at a standstill.

If the lid is opened via the emergency lid release system during a centrifuge run, the centrifuge will be switched off immediately and decelerate in an unbraked manner.

7 Malfunctions and error correction

7.2 Table of error codes

Error no.	Kind of error	Measures	Note
1-9	System error	<ul style="list-style-type: none"> Allow to slow down Power off/on 	All these errors stop the centrifuge or cause it to decelerate brakeless
10-19	Speedometer error	<ul style="list-style-type: none"> Allow to slow down Power off/on 	
20-29	Motor error	<ul style="list-style-type: none"> Power off Ensure ventilation 	
30-39	EEPROM error	<ul style="list-style-type: none"> Allow to slow down Power off/on 	With error 34, 35, and 36, the centrifuge will stop; with error 37 and 38 only an error message will be given
40-45	Temperature error (only for refrigerated centrifuges)	<ul style="list-style-type: none"> Allow to slow down Power off Allow to cool down Provide better ventilation (only air cooled centrifuges) Provide sufficient water throughput (only water cooled centrifuges) 	
46-49	Imbalance error	<ul style="list-style-type: none"> Allow to slow down Power off Eliminate the imbalance 	
50-59	Lid error	<ul style="list-style-type: none"> Press lid key Close lid Remove foreign matter from the opening of the lid lock device 	With error 50 and 51, the centrifuge will stop
60-69	Process error	<ul style="list-style-type: none"> Allow to slow down Power off/on 	With error 60 message "power failure during run", with error 61, the message "stop after power on"
70-79	Communication error	<ul style="list-style-type: none"> Allow to slow down Power off/on 	
80-89	Parameter error	<ul style="list-style-type: none"> Power off Allow to cool down Provide for better ventilation 	With error 83, error message only
90-99	Other errors	<ul style="list-style-type: none"> Check connections Provide sufficient water throughput (only water cooled centrifuges) 	



NOTE

If it is impossible to eliminate the errors, contact the service!

7.3 Service contact

In the event of queries, malfunctions, or spare part enquiries:

From Germany:

Contact

Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode (Germany)
Tel. +49 (0) 55 22 / 50 07-44 44
E-mail: support.lab@sigma-zentrifugen.de

Outside Germany:

Contact our agency in your country. All agencies are listed at
www.sigma-zentrifugen.de → [Sales Partners]



NOTE

- If you would like to utilise our service, please state the type of your centrifuge and its serial number.

8 Maintenance and service

8 Maintenance and service

The centrifuge, rotor, and accessories are subject to high mechanical stress. Thorough maintenance performed by the user extends the service life and prevents premature failure.



CAUTION

If corrosion or other damage occurs due to improper care, the manufacturer cannot be held liable or subject to any warranty claims.

- Use soap water or other water-soluble, mild cleaning agents with a pH value between 6 and 8 for cleaning the centrifuge and accessories (see also).
- Avoid corrosive and aggressive substances.
- Do not use solvents.
- Do not use agents with abrasive particles.
- Do not expose the centrifuge and rotors to intensive UV radiation or thermal stress (e.g. by heat generators).

8.1 Maintenance

8.1.1 Centrifuge

- Unplug the mains power plug before cleaning.
- Carefully remove all liquids, including water and particularly all the solvents, acids, and alkaline solutions from the rotor chamber using a cloth in order to avoid damage to the motor bearings.
- If the centrifuge has been contaminated with toxic, radioactive, or pathogenic substances, clean the rotor chamber immediately with a suitable decontamination agent (depending on the type of contamination).



WARNING

Take suitable precautions for your own safety if there is a risk of toxic, radioactive, or pathogenic contamination.

- Grease the motor shaft slightly after cleaning (grease for load-bearing bolts part no. 70284).

8.1.2 Condenser (only refrigerated centrifuges)

In order to cool the refrigerant that is compressed by the refrigeration unit, a lamellar condenser is used. It is cooled by air.

Dust and dirt obstruct the cooling flow of air. The dust on condenser pipes and lamellas reduces the heat exchange and thus the performance of the refrigeration unit.

This is why the installation site should be as clean as possible.

- Check the condenser at least once a month for dirt and clean it if necessary.
- If you have any queries, please contact service (see chapter 7.3 - "Service contact").

8.1.3 Accessories



CAUTION

For the care of the accessories, special safety measures must be considered as these are measures that will ensure operational safety at the same time!



WARNING

All swing-out rotors applicable for this centrifuge and the angle rotor 12510 weigh more than 18 kg.

- Always lift the rotors with a lifting device or with a sufficient number of people helping you.
- Immediately rinse off the rotor, buckets, or accessories under running water if they have come into contact with any liquids that may cause corrosion. Use a brush for test tubes to clean the bores of angle rotors. Turn the rotor upside down and allow it to dry completely.
- Clean the accessories outside the centrifuge once a week or preferably after each use. Adapters should be removed, cleaned and dried.



CAUTION

Do not clean the accessories in a dishwasher!

Cleaning in a dishwasher removes the anodised coating; the result is cracking in areas that are subject to stress.

- If the rotors or accessories have been contaminated with toxic, radioactive, or pathogenic substances, clean them immediately with a suitable decontamination agent (depending on the type of contamination). Take suitable precautions for your own safety if there is a risk of toxic, radioactive, or pathogenic contamination.
- Dry the accessories with a soft cloth or in a drying chamber at approx. 50°C.

8 Maintenance and service

8.1.3.1 Plastic accessories

The chemical resistance of plastic decreases with rising temperatures (see chapter 11.5 - "Resistance data").

- If solvents, acids, or alkaline solutions have been used, clean the plastic accessories thoroughly.



WARNING

Plastic accessories must not be greased!

8.1.4 Rotors, buckets and carriers

Rotors, buckets and carriers are produced with the highest precision, in order to withstand the permanent high stress from high gravitational fields. Chemical reactions as well as stress-corrosion (combination of oscillating pressure and chemical reaction) can affect or destroy the metals. Barely detectable cracks on the surface can expand and weaken the material without any visible signs.

- Check the material regularly (at least once a month) for
 - cracks
 - visible damage of the surface
 - pressure marks
 - signs of corrosion
 - other changes.
- Check the bores of the rotors and multiple carriers.
- Replace any damaged components immediately for your own safety.
- Grease the rotor tie-down screw after cleaning with grease for load bearing bolts (part no. 70284) and spread the grease with a cloth.

8.1.5 Load bearing bolts

Only greased load bearing bolts ensure the even swinging of the buckets and, therefore, the quiet run of the centrifuge. Non-greased bolts can lead to a system shut-down due to imbalances.

Apply a small quantity of grease (part no. 70284) to the load bearing bolts of the rotor and buckets after each cleaning.

8.1.6 Glass breakage



CAUTION

In the case of glass breakage, immediately remove all glass particles (e.g. with a vacuum cleaner). Replace the rubber cushions since even thorough cleaning will not remove all glass particles.

Glass particles will damage the surface coating (e.g. anodising) of the buckets, which will then lead to corrosion.

Glass particles in the rubber cushions of the buckets will cause glass breakage again.

Glass particles on the pivot bearing of the load-bearing bolts prevent the buckets and carriers from swinging evenly, which will cause an imbalance.

Glass particles in the rotor chamber will cause metal abrasion due to the strong air circulation. This metal dust will not only pollute the rotor chamber, rotor, and materials to be centrifuged but also damage the surfaces of the accessories, rotors, and rotor chamber.

In order to completely remove the glass particles and metal dust from the rotor chamber:

- Grease the upper third of the rotor chamber with e.g. Vaseline.
- Then, let the rotor rotate for a few minutes at a moderate speed (approx. 2000 rpm). The glass and metal particles will now collect at the greased part.
- Remove the grease with the glass and metal particles with a cloth.
- If necessary, repeat this procedure.

8.2 Sterilisation and disinfection of the rotor chamber and accessories

- Use commercially-available disinfectants such as, for example, Sagrotan[®], Buraton[®], or Terralin[®] (available at chemist's shops or drugstores).
- The centrifuge and the accessories consist of various materials. A possible incompatibility must be considered.
- Before using cleaning or decontamination agents that were not recommended by us, contact the manufacturer to ensure that such a procedure will not damage the centrifuge.
- For autoclaving, consider the continuous heat resistance of the individual materials (see chapter 8.2.1 - "Autoclaving").

Please contact us if you have any queries (see chapter 7.3 - "Service contact").



DANGER

If dangerous materials (e.g. infectious and pathogenic substances) are used, the centrifuge and accessories must be disinfected.

8 Maintenance and service

8.2.1 Autoclaving

The service life of the accessories essentially depends on the frequency of autoclaving and use.

- Replace the accessories immediately when the parts show changes in colour or structure or in the occurrence of leaks etc.
- During autoclaving, the caps of the tubes must not be screwed on in order to avoid the deformation of the tubes.



NOTE

It cannot be excluded that plastic parts, e.g. lids or carriers, may deform during autoclaving.

Accessories	Max. temp. (°C)	Min. time (min)	Max. time (min)	Max. cycles
Aluminium buckets	134-138	3	5	-
Aluminium rotors	134-138	3	5	-
Glass tubes	134-138	3	40	-
Polyallomer / polycarbonate rectangular carriers	115-118	30	40	-
Polyallomer / polycarbonate round carriers	115-118	30	40	-
Polyamide buckets	115-118	30	40	10
Polycarbonate / polyallomer lids for angle rotors	115-118	30	40	20
Polycarbonate tubes	115-118	30	40	20
Polyphenylsulfone (PPSU) caps for buckets	134-138	3	5	100
Polypropylene balance adapter for blood-bag systems	115-118	30	40	n.s.
Polypropylene copolymer tubes	115-121	30	40	20
Polypropylene rectangular carriers	115-118	30	40	-
Polypropylene rotors	115-118	30	40	20
Polypropylene round carriers	115-118	30	40	-
Polysulfone caps for buckets	134-138	3	5	100
Polysulfone lids for angle rotors	134-138	3	5	100
Rubber adapters	115-118	30	40	-
Stainless-steel balance weight for blood-bag systems	121	30	30	n.s.
Teflon tubes	134-138	3	5	100

8.3 Service

**DANGER**

In the event of service work that requires the removal of the panels, there is a risk of electric shock or mechanical injury. Only qualified specialist personnel is authorised to perform this service work.

The centrifuge is subject to high mechanical stress. In order to be able to withstand this high level of stress, high-quality components were used during the production of the centrifuge. Nevertheless, wear cannot be excluded and it may not be visible from the outside. Especially the rubber parts that are – among other things – part of the motor suspension, are subject to ageing.

This is why we recommend having the centrifuge checked by the manufacturer during an inspection once per year in the operating state and once every three years in the dismantled state. Rubber parts should be replaced after three years.

Information and appointments:

In Germany:

Contact

Sigma Laborzentrifugen GmbH

An der Unteren Söse 50

37520 Osterode (Germany)

Tel. +49 (0) 55 22 / 50 07-44 44

E-mail: support.lab@sigma-zentrifugen.de

Outside Germany:

Contact our agency in your country. All agencies are listed at www.sigma-zentrifugen.de → [Sales Partners]

**NOTE**

- If you would like to utilise our service, please state the type of your centrifuge and its serial number.

8 Maintenance and service

8.4 Return of defective parts

Although we exercise great care during the production of our products, it may be necessary to return a unit or accessory to the manufacturer.

In order to ensure the quick and economical processing of returns of centrifuges, spare parts, or accessories, we require complete and extensive information concerning the process. Please fill in the following forms completely, sign them, enclose them with the return package, and send them together with the product to:

Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode (Germany)

1. Declaration of decontamination

As a certified company and due to the legal regulations for the protection of our employees and of the environment, we are obliged to certify the harmlessness of all incoming goods. For this purpose, we require a declaration of decontamination.

- The form must be filled in completely and signed by authorised and specialised personnel only.
- Affix the original form in a clearly visible manner to the outside of the packaging.



NOTE

We will return the unit if no declaration of decontamination is provided!

2. Form for the return of defective parts

This form is for the product-related data. They facilitate the assignment, and they enable the quick processing of the return. If several parts are returned together in one packaging, please enclose a separate problem description for every defective part.

- A detailed problem description is necessary in order to perform the repair quickly and economically.



NOTE

If the form does not include a description of the malfunction, neither a refund nor a credit note can be issued. In this case, we reserve the right to return the parts to you at your expense.

- Upon request, we will prepare and submit to you a cost estimate prior to performing the repair. Please confirm such cost estimate within 14 days. If the cost estimate has still not been confirmed after 4 weeks, we will return the unit. Please note that you must bear the incurred costs.



The unit must be packaged in a transport-safe manner. Please use the original packaging, if at all possible.

If the product is dispatched to us in unsuitable packaging, you will be charged the cost for returning it to you in new packaging.

The forms can be downloaded online from
www.sigma-zentrifugen.de → [Service] → [Overhaul and repair].

9 Disposal

9 Disposal

9.1 Disposal of the centrifuge



In accordance with the directive 2002/96/EC, SIGMA centrifuges are marked with the symbol shown to the left. This symbol means that it is not permissible to dispose of the unit among household waste.

- You can return these centrifuges free of cost to Sigma Laborzentrifugen GmbH.
- Ensure that the unit is decontaminated. Fill in a declaration of decontamination (see chapter 8.4 - "Return of defective parts").
- Comply with any other applicable local rules and regulations.

9.2 Disposal of the packaging

- Use the packaging to return the centrifuge for disposal or
- dispose of the packaging, after having separated the individual materials.
- Comply with all local rules and regulations.

10 Technical data

Manufacturer	SIGMA Laborzentrifugen GmbH An der Unteren Söse 50 37520 Osterode (Germany)	
Type:	8KS	8KBS
<u>Connection requirements</u>		
Electr. connection:	see name plate	see name plate
Protection class:	I	I
IP code:	20	20
Connected load (kVA):	7.5	7.5
Power consumption (kW):	6.5	6.5
Max. current consumption (A):	11.5 (at 3x400V / 50 Hz)	11.5 (at 3x400V / 50 Hz)
Input fuse (AT):	16.0	16.0
<u>Performance data</u>		
Max. speed (rpm):	10,500	5,100
Max. capacity (ml):	12,000	12,000
Max. gravitational field (x g):	20,954	8,578
Max. kin. energy (Nm):	275,321	275,321
<u>Other parameters</u>		
Time range:	10 sec to 99 h 59 min 59 sec, short run, continuous run	10 sec to 99 h 59 min 59 sec short run, continuous run
Temperature range:	-20 to +40°C	-20 to +40°C
Storage locations:	60	60
<u>Physical data</u>		
Height (mm):	980	980
Height with open lid (mm):	1,690	1,690
Width (mm):	810	810
Depth (mm):	910	910
Weight (kg):	420	420
Noise level (dB(A)):	< 73 (at maximum speed)	< 73 (at maximum speed)
<u>Special equipment:</u> <u>Water cooling system</u>		
Tap connections (inch):	2 x ¾	2 x ¾
Inlet pressure (bar):	1.5 to 5.0	1.5 to 5.0
Max. flow rate (l/min):	5	5
Max. temperature at water inlet (°C):	20	20

10 Technical data

10.1 Ambient conditions

- The figures are valid for an ambient temperature of $+23^{\circ}\text{C} \pm 2^{\circ}\text{C}$ and a nominal voltage $\pm 10\%$. The minimum temperature is $< +4^{\circ}\text{C}$ and depends on the rotor type, speed, and ambient temperature.



NOTE

At a nominal voltage of 100V or 200V, a tolerance of +10% / -5% applies.

- For indoor use only.
- Allowable ambient temperature $+5^{\circ}\text{C}$ to $+35^{\circ}\text{C}$.
- Max. relative humidity of air 80% up to 31°C with a linear decrease to 67% relative humidity of air at 35°C .
- Maximum altitude 2,000 m above sea level.

10.2 Technical documentation

For environmental reasons, the comprehensive technical documentation of the centrifuge (e.g. circuit diagrams) and the safety data sheets of the manufacturers of refrigerant and lubricants are not attached to this documentation.

You can order these documents from our service department.

11 Appendix

11.1 Range of accessories



If the swing-out rotor 11805 with a windshield is used in combination with the blood bag bucket 13860, it can also be used without the windshield cover.

For all other rotor-bucket-combinations, the windshield must be closed.

Part no.	Description	Max. speed (rpm)	Max. gravitational field (x g)
11805	Swing-out rotor with wind shield, 6 place, for buckets no. 13845 ² , max. radius 25.2 cm, min. radius 14.1 cm	5 100	7 328
	13850, max. radius 29.5 cm, min. radius 17.0 cm	5 100	8 578
	13860, max. radius 28.7 cm, min. radius 14.5 cm	4 100	5 394
	Attention! Please refer to the note concerning the windshield (see above)!		
11806	Swing-out rotor, 6 place, for buckets no. 13845 ³ , max. radius 25.2 cm, min. radius 14.1 cm	3 700	3 857
	13850, max. radius 29.5 cm, min. radius 17.0 cm	3 700	4 515
	13855, max. radius 29.4 cm, min. radius 14.2 cm	3 300	3 579
	13860, max. radius 28.7 cm, min. radius 14.5 cm	3 700	4 393
12505-H ⁴	Angle rotor 6 x 500 ml, for bottles no. 13507, 15508 and adapters 14144 - 14155, sealable with hermetic aluminium lid no. 17891, Radius max. 16.6 cm, radius min. 4.1 cm, angle 30°	10 500	20 461
12510-H ⁵	Angle rotor 6 x 1000 ml, incl. steel buckets no. 13511, for bottles no. 13510, 15921 and adapters 14162 - 14182, sealable with hermetic aluminium lid no. 17892 Radius max. 20.6 cm, radius min. 5.5 cm, angle 25°	7 000	11 285
13845 ⁶	Bucket for microtiter plates, incl. plate holder no. 17981 and cap no. 17122, max. height of plates approx. 105 mm, 1 set = 2 pcs.		
13850	Round bucket, sealable with cap no. 17175, for adapters Ø 115 mm, no. 17775, 17777, 17780, adapters no. 13851, 13852 and bottles no. 13840, 15930, 15932 without 17175, 1 set = 2 pcs.		
13855	Double bucket 2 x 1000 ml, for tubes no. 13510, 15920 and 15921 with 13511, for adapters Ø 100 mm, no. 17650 - 17699, 1 set = 2 pcs.		

² Not for Sigma 8KBS

³ Not for Sigma 8KBS

⁴ Not for Sigma 8KBS incl. accessories

⁵ Not for Sigma 8KBS incl. accessories

⁶ Not for Sigma 8KBS

11 Appendix

Accessories for blood bags

Part no.	Description	Max. speed (rpm)	Max. gravitational field (x g)
13860	Double bucket 2 x 300 – 900 ml, incl. adapter no. 13865 for 2 blood bag systems, incl. 2 sets balance weights 5 g no. 17751 and 1 set balance weights 10 g no. 17752, 1 set = 2 pcs.		
13865	Adapter for 2 blood bags 500 - 800 ml incl. 2 sets balance weights 5 g no. 17751 and 1 set balance weights 10 g no. 17752, 1 set = 2 pcs., service life max. 1.000 cycles		
13866	Adapter for 2 blood bags 900 ml, 1 set = 2 pcs., service life max. 1.000 cycles		
17751	Balance weight 5 g, suitable in 13865, 13866, 1 set = 2 pcs.		
17752	Balance weight 10 g, suitable in 13865, 13866, 1 set = 2 pcs.		
17753	Balance weight 50 g, suitable in 17765, 17766, 1 set = 2 pcs.		
17754	Balance weight 50 g, suitable in 17765, 17766, 1 set = 2 pcs.		
17765	Balance adapter, replacement for one blood bag, 450 g, suitable in 13865, 1 set = 2 pcs.		
17766	Balance adapter, replacement for one blood bag, 600 g, suitable in 13866, 1 set = 2 pcs.		

Adapters and stainless steel tubes

Part no.	Description
13055	Stainless steel tube 50 ml, Ø 28.5 x 101.5 mm, sealable with cap no. 17054
13085	Stainless steel tube 85 ml, Ø 38 x 103 mm, sealable with cap no. 17185
13255	Stainless steel bottle 250 ml, Ø 61.4 x 125 mm, sealable with cap no. 17256
13507	Stainless steel tube 500 ml, Ø 69,5 x 152 mm, sealable with cap no. 17256
13510	Stainless steel bottle 1000 ml, Ø 100 x 167 mm, incl. cap no. 17520
13511	Stainless steel tube Ø 100 x 152 mm, for tubes no. 15920, 15921, use without cap in rotor 12510-H
13840	Stainless steel tube 1500 ml, Ø 115 x 165 mm, 1300 ml without cap
13851	Adapter for stainless steel bottle 1000 ml, no. 13510, adapters Ø 100 mm, no. 17650 - 17699, 1 set = 2 pcs.
13852	Adapter for 1000 ml bottles Ø 98 mm no. 15920, 15921, 1 set = 2 pcs.

The complete list of accessories can be downloaded from
www.sigma-zentrifugen.de.

11.1.1 Maximum speed for tubes

Some tubes, such as centrifuge glass tubes, microtubes, culture tubes, fluoropolymer tubes and especially high-volume tubes can be used in our rotors, buckets, and adapters at higher speeds than their breaking limit.

- Always fill the tubes up to their useful volume (= the volume that is stated for the tube).
- When installing the 500 ml bottles use the supplied supporting rings.



CAUTION

At speeds above 8,000 rpm, there is an increased risk of breakage, in particular for 250 and 500 ml bottles!

11.1.2 Rotor radii

The information in the accessories table concerning the radius refers to the values of the respective rotor as shown below. The radius calculation is described in (see chapter 2.2.2.1 - "Speed, radius, and relative centrifugal force").

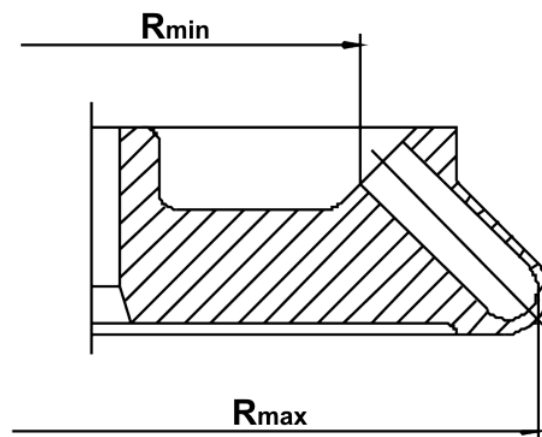


Fig. 44: Minimum and maximum radius of an angle rotor

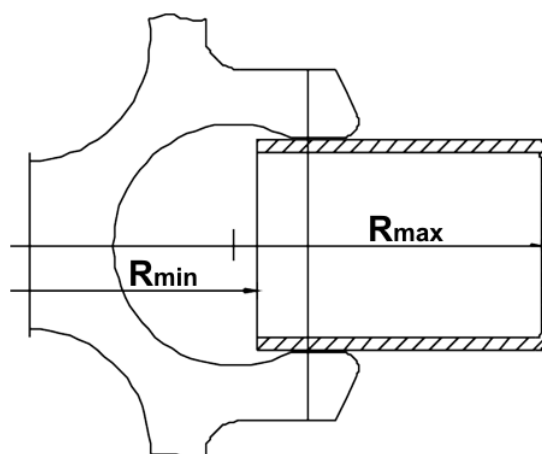


Fig. 45: Minimum and maximum radius of a swing-out rotor

11 Appendix

11.2 Speed-gravitational-field-diagram

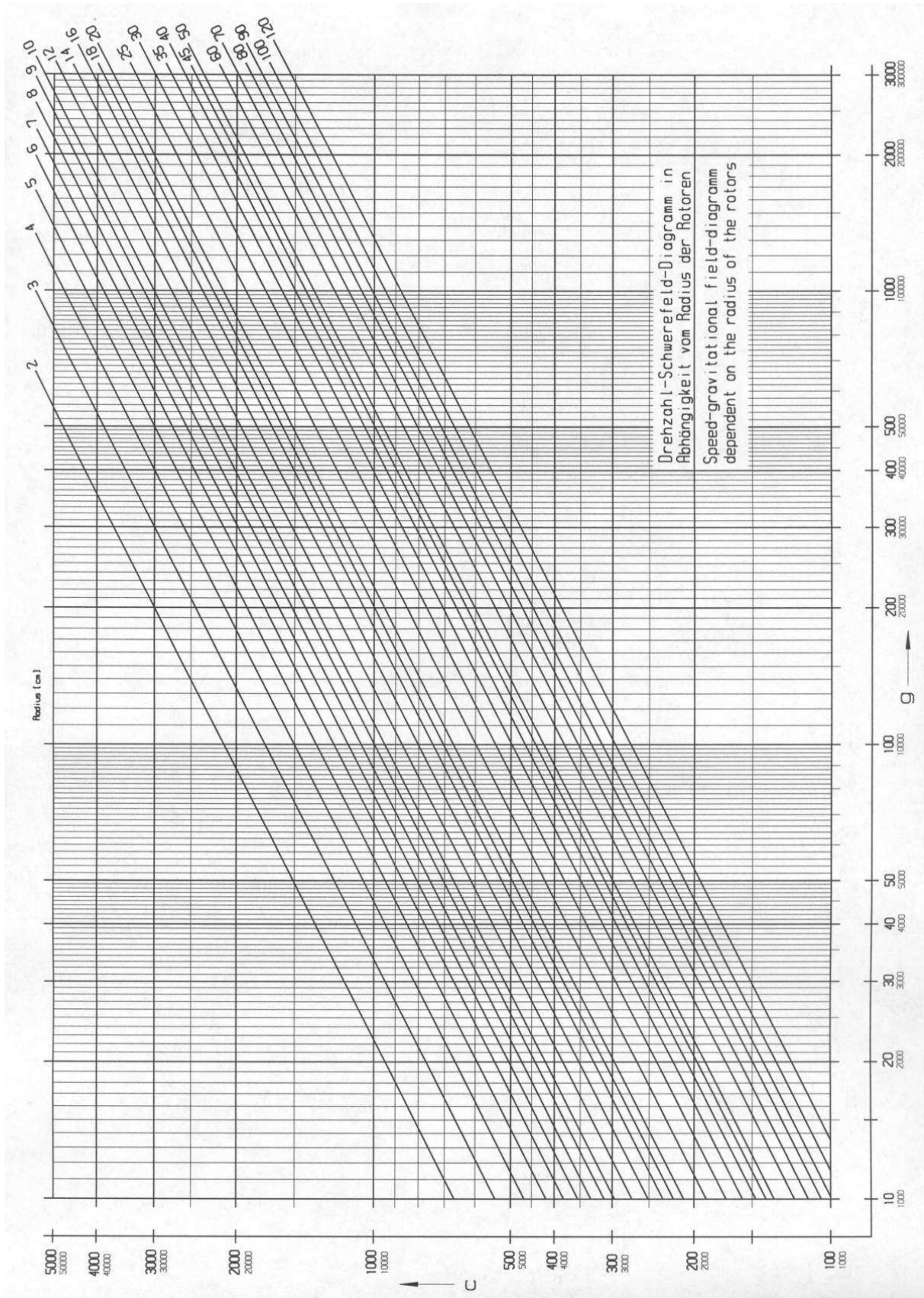


Fig. 46: Speed-gravitational-field-diagram

11.3 Acceleration and deceleration curves

Linear as well as quadratic curves are numbered in the direction of increasing acceleration (from right to left).

The deceleration curves are inverted images of the acceleration curves and are assigned the same numbers. An exception is curve 0. It decelerates brakeless (spin-out).

In general, the runtime, until the set speed is reached, depends on the moment of inertia of the rotor.

Linear curves

The slope of the fixed acceleration curves defines the time that is required to accelerate the rotor by 1,000 rpm.

Curve 9 is a special case compared to the other curves. The centrifuge accelerates with maximum power. The runtime, until the set speed is reached, depends solely on the moment of inertia of the rotor.

Linear curve no.	Slop
0	4 [rpm/sec]
1	6 [rpm/sec]
2	8 [rpm/sec]
3	17 [rpm/sec]
4	25 [rpm/sec]
5	33 [rpm/sec]
6	50 [rpm/sec]
7	100 [rpm/sec]
8	200 [rpm/sec]
9	1.000 [rpm/sec]

Fig. 47: Slope of linear curves

Quadratic curves

Curve 19 is a special case compared to the other curves. The centrifuge accelerates with maximum power. The runtime depends solely on the moment of inertia of the rotor.

Quadratic curve no.	Time until 1,000 rpm	Slope as of 1,000 rpm
10	500 sec	4 [rpm/sec]
11	333 sec	6 [rpm/sec]
12	250 sec	8 [rpm/sec]
13	118 sec	17 [rpm/sec]
14	80 sec	25 [rpm/sec]
15	60 sec	33 [rpm/sec]
16	40 sec	50 [rpm/sec]
17	20 sec	100 [rpm/sec]
18	10 sec	200 [rpm/sec]
19	2 sec	1.000 [rpm/sec]

Fig. 48: Slope of quadratic curves

11 Appendix

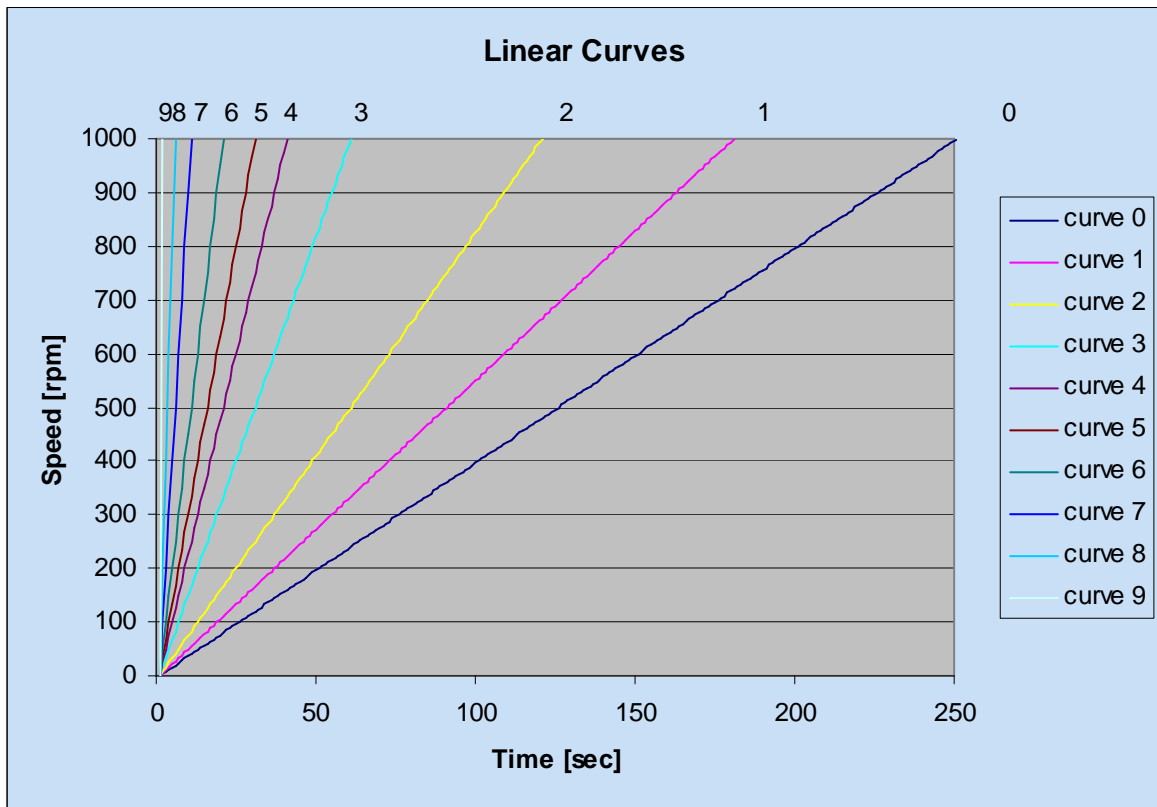


Fig. 49: Diagram of linear curves

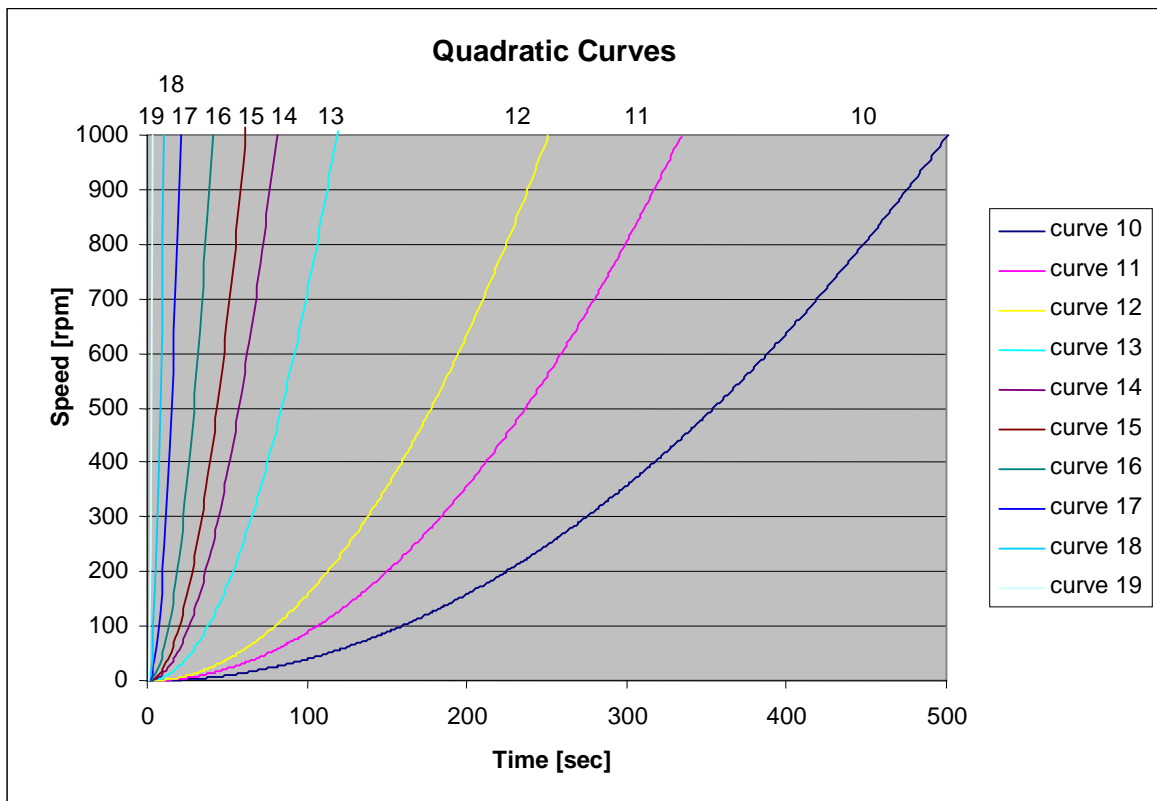



Fig. 50: Diagram of quadratic curves

11.4 Table of the service life of rotors and accessories

Rotors and accessories with a different service life					
If no other data concerning the service life are engraved on the rotor or accessory, rotors and buckets must be checked by the manufacturer after 10 years. After 50,000 cycles, rotors must be scrapped for safety reasons.					
Rotor / bucket	Cycles	Service life ("Exp.Date")	Autoclaving	Suitable for centrifuge	Remarks
9100	15,000			4-15C, 4K15C, 4-16, 4-16S, 4-16K, 4-16KS, 6-15, 6K15, 6-16, 6-16K	without engraving, only "spincontrol professional" and "spincontrol S"
11026		7 years		1-14, 1-14K	
12033		5 years		1-16 Edition, 1-16K Edition	
12082		7 years		1-14, 1-14K	
12083		7 years		1-14, 1-14K	
12084		7 years		1-14, 1-14K	
12085		7 years		1-14, 1-14K	
12092		5 years	20x	1-14, 1-14K	
12093		5 years	20x	1-14, 1-14K	
12094		5 years	20x	1-14, 1-14K	
12096		5 years	20x	1-14, 1-14K	
12101		5 years	20x	1-15, 1-15K, 1-15P, 1-15PK	
12124		5 years	20x	1-15, 1-15K, 1-15P, 1-15PK	
12126		5 years	20x	1-15, 1-15K, 1-15P, 1-15PK	
12134		5 years	20x	1-16, 1-16K	
12135		5 years	20x	1-16, 1-16K	
12137		5 years	20x	1-16, 1-16K	
12500		7 years		6-15, 6K15, 6-16, 6-16K	
12600		7 years		6-16S, 6-16KS	
13218	20,000			4-16, 4-16S, 4-16K, 4-16KS, 6-16, 6-16S, 6-16K, 6-16KS	
13296		5 years	10x	2-6, 2-6E, 2-7	
13299		5 years	10x	2-6, 2-6E, 2-7, 2-16P, 2-16KL, 2-16KHL	
13635	25,000			6-16, 6-16K, 6-16S, 6-16KS	
13845	20,000			8K, 8KS	
13850	10,000			8K, 8KS, 8KBS	
13860	35,000			8K, 8KS, 8KBS	
13864	1,000			8K, 8KS	without engraving
13865	1,000			8K, 8KS	without engraving
13866	1,000			8K, 8KS	without engraving

11 Appendix

11.5 Resistance data



NOTE

The data refer to resistance at 20°C.

- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant	Medium	Formula	Concentration [%]	High Density Polyethylene	Polyamide	Polycarbonate	Polyoxymethylene	Polypropylene	Polysulfone	Polyvinyl chloride, hard	Polyvinyl chloride, weak	Polytetrafluorethylene	Acrylonitrile-butadiene- caoutchouc	Aluminium
				HDPE	PA	PC	POM	PP	PSU	PVC	PVC	PTFE	NBR	AL
	Acetaldehyde	C ₂ H ₄ O	40	3	2	4	2	3	4	4	-	1	4	1
	Acetamide	C ₂ H ₅ NO	saturated	1	1	4	1	1	4	4	-	1	-	1
	Acetone	C ₃ H ₆ O	100	1	1	4	1	1	4	4	-	1	4	1
	Acrylonitrile	C ₃ H ₃ N	100	1	1	4	3	3	4	4	4	1	4	1
	Allyl alcohol	C ₃ H ₆ O	96	1	3	3	2	2	2	2	4	1	1	1
	Aluminium chloride	AlCl ₃	saturated	1	3	2	4	1	-	1	-	1	1	4
	Aluminium sulfate	Al ₂ (SO ₄) ₃	10	1	1	1	3	1	1	1	1	1	1	1
	Ammonium chloride	(NH ₄)Cl	aqueous	1	1	1	2	1	1	1	1	1	1	3
	Ammonium hydroxide	NH ₃ + H ₂ O	30	1	3	4	1	1	2	1	-	1	-	1
	Aniline	C ₆ H ₇ N	100	1	3	4	1	2	4	4	4	1	4	1
	Anisole	C ₇ H ₈ O	100	3	4	4	1	4	4	2	-	1	4	1
	Antimony trichloride	SbCl ₃	90	1	4	1	4	1	-	1	-	1	-	4
	Benzaldehyde	C ₇ H ₆ O	100	1	3	4	1	1	3	4	4	1	4	1
	Benzene	C ₆ H ₆	100	3	2	4	1	3	4	4	-	1	4	1
	Boric acid	H ₃ BO ₃	aqueous	1	3	1	2	1	-	-	-	1	1	1
	Butyl acrylate	C ₇ H ₁₂ O ₂	100	1	2	4	2	3	4	4	4	1	-	1
	Butyl alcohol, normal	C ₄ H ₁₀ O	100	1	1	2	1	1	2	2	4	1	1	1
	Calcium chloride	CaCl ₂	alcoholic	1	4	2	3	1	-	-	4	1	1	3
	Carbon disulfide	CS ₂	100	4	3	4	2	4	4	4	4	1	3	1
	Carbon tetrachloride (TETRA)	CCl ₄	100	4	4	4	2	4	4	4	4	1	3	1
	Chlorine	Cl ₂	100	4	4	4	4	4	4	4	4	1	-	3
	Chlorine water	Cl ₂ x H ₂ O		3	4	4	4	3	-	3	3	1	-	4
	Chlorobenzene	C ₆ H ₅ Cl	100	3	4	4	1	3	4	4	4	1	4	1
	Chloroform	CHCl ₃	100	3	3	4	4	3	4	4	4	1	4	3

11 Appendix

- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant	Medium	Formula	Concentration [%]	High Density Polyethylene	Polyamide	Polycarbonate	Polyoxymethylene	Polypropylene	Polysulfone	Polyvinyl chloride, hard	Polyvinyl chloride, weak	Polytetrafluorethylene	Acrylonitrile-butadiene- caoutchouc	Aluminium
				HDPE	PA	PC	POM	PP	PSU	PVC	PVC	PTFE	NBR	AL
	Chromic acid	CrO ₃	10	1	4	2	4	1	4	1	-	1	4	1
	Chromic potassium sulphate	KCr(SO ₄) ₂ x 12H ₂ O	saturated	1	2	1	3	1	-	1	-	1	-	3
	Citric acid	C ₆ H ₈ O ₇	10	1	1	1	2	1	1	1	1	1	1	1
	Citric acid	C ₆ H ₈ O ₇	50	1	3	1	2	1	-	-	-	1	1	1
	Copper sulphate	CuSO ₄ x 5H ₂ O	10	1	1	1	1	1	1	1	1	1	1	4
	Cyclohexanol	C ₆ H ₁₂ O	100	1	1	3	1	1	1	1	4	1	2	1
	Decane	C ₁₀ H ₂₂	100	-	1	2	1	3	-	-	-	1	2	1
	Diaminoethane	C ₂ H ₈ N ₂	100	1	1	3	1	1	-	3	4	1	1	1
	Diesel fuel	—	100	1	1	3	1	1	-	1	3	1	1	1
	Dimethyl formamide (DMF)	C ₃ D ₇ NO	100	1	1	4	1	1	4	3	-	1	3	1
	Dimethyl sulfoxide (DMSO)	C ₂ H ₆ SO	100	1	2	4	1	1	4	4	-	1	-	1
	Dimethylaniline	C ₈ H ₁₁ N	100	-	3	4	2	4	-	-	-	1	-	1
	Dioxane	C ₄ H ₈ O ₂	100	2	1	4	1	3	2	3	4	1	3	1
	Dipropylene glycol (mono)methyl ether	C ₄ H ₁₀ O	100	3	1	4	1	4	4	4	4	1	-	1
	Ethyl acetate	C ₄ H ₈ O ₂	100	1	1	4	1	1	4	4	4	1	4	1
	Ethylene chloride	C ₂ H ₄ Cl ₂	100	3	3	4	1	3	4	4	4	1	-	1
	Ferrous chloride	FeCl ₂	saturated	1	3	1	3	1	1	1	1	1	-	4
	Formaldehyde solution	CH ₂ O	30	1	3	1	1	1	-	-	-	1	2	1
	Formic acid	CH ₂ O ₂	100	1	4	3	4	1	3	3	1	1	2	1
	Furfural	C ₅ H ₄ O ₂	100	1	3	3	2	4	-	-	-	1	4	1
	Gasoline	C ₅ H ₁₂ - C ₁₂ H ₂₆	100	2	1	3	1	3	3	2	-	1	1	1
	Glycerol	C ₃ H ₈ O ₃	100	1	1	3	1	1	1	1	2	1	1	1
	Heptane, normal	C ₇ H ₁₆	100	2	1	1	1	2	1	2	4	1	1	1
	Hexane, n-	C ₆ H ₁₄	100	2	1	2	1	2	1	2	4	1	1	1
	Hydrogen chloride	HCl	5	1	4	1	4	1	1	1	-	1	2	4
	Hydrogen chloride	HCl	concentrated	1	4	4	4	1	1	2	3	1	4	4
	Hydrogen peroxide	H ₂ O ₂	3	1	3	1	1	1	1	1	-	1	3	3
	Hydrogen peroxide	H ₂ O ₂	30	1	4	1	4	1	1	1	-	1	3	3
	Hydrogen sulphide	H ₂ S	10	1	1	1	1	1	1	1	3	1	3	1
	Iodine, tincture of	I ₂		1	4	3	1	1	-	4	4	1	1	1

11 Appendix

- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant	Medium	Formula	Concentration [%]	High Density Polyethylene	Polyamide	Polycarbonate	Polyoxymethylene	Polypropylene	Polysulfone	Polyvinyl chloride, hard	Polyvinyl chloride, weak	Polytetrafluorethylene	Acrylonitrile-butadiene- caoutchouc	Aluminium
				HDPE	PA	PC	POM	PP	PSU	PVC	PVC	PTFE	NBR	AL
	Isopropyl alcohol	C ₃ H ₈ O	100	1	1	1	1	1	1	1	4	1	-	2
	Lactic acid	C ₃ H ₆ O ₃	3	1	3	1	2	1	1	2	-	1	1	1
	Magnesium chloride	MgCl ₂	10	1	1	1	1	1	1	1	1	1	1	1
	Mercuric chloride	HgCl ₂	10	1	4	1	3	1	1	1	1	1	1	4
	Mercury	Hg	100	1	1	1	1	1	1	1	3	1	1	3
	Methyl acetate	C ₃ H ₆ O ₂	100	1	1	4	2	1	-	4	4	1	-	1
	Methyl alcohol	CH ₄ O	100	1	2	4	1	1	3	1	3	1	2	1
	Methyl benzene	C ₇ H ₈	100	3	1	4	1	3	4	4	4	1	4	1
	Methyl ethyl ketone (MEK)	C ₄ H ₈ O	100	1	1	4	1	1	4	4	4	1	4	1
	Methylene chloride	CH ₂ Cl ₂	100	4	3	4	3	3	4	4	4	1	-	1
	Mineral oil	—	100	1	1	1	1	1	1	1	-	1	1	1
	Nitric acid	HNO ₃	10	1	4	1	4	1	1	1	-	1	4	3
	Nitric acid	HNO ₃	100	4	4	4	4	4	-	4	-	1	4	1
	Nitrobenzene	C ₆ H ₅ NO ₂	100	3	4	4	3	2	4	4	4	1	4	1
	Oleic acid	C ₁₈ H ₃₄ O ₂	100	1	1	1	2	1	-	1	-	1	3	1
	Oxalic acid	C ₂ H ₂ O ₄ x 2H ₂ O	100	1	3	1	4	1	1	1	1	1	2	1
	Ozone	O ₃	100	3	4	1	4	3	1	1	-	1	4	2
	Petroleum	—	100	1	1	3	1	1	1	1	3	1	1	1
	Phenol	C ₆ H ₆ O	10	1	4	4	4	1	4	1	3	1	3	1
	Phenol	C ₆ H ₆ O	100	2	4	4	4	1	3	4	3	1	3	1
	Phosphoric acid	H ₃ PO ₄	20	1	4	2	4	1	-	-	-	1	2	4
	Phosphorus pentachloride	PCl ₅	100	-	4	4	4	1	-	4	4	1	-	1
	Potassium hydrogen carbonate	CHKO ₃	saturated	1	1	2	1	1	-	-	-	1	-	4
	Potassium hydroxide	KOH	30	1	1	4	3	1	1	1	1	1	-	4
	Potassium hydroxide	KOH	50	1	1	4	3	1	1	1	1	1	-	4
	Potassium nitrate	KNO ₃	10	1	1	1	1	1	-	-	-	1	1	1
	Potassium permanganate	KMnO ₄	100	1	4	1	1	1	-	1	-	1	3	1
	Pyridine	C ₅ H ₅ N	100	1	1	4	1	3	4	4	4	1	4	1
	Resorcinol	C ₆ H ₆ O ₂	5	1	4	2	3	1	4	2	-	1	-	2
	Silver nitrate	AgNO ₃	100	1	1	1	1	1	1	1	1	1	2	4

11 Appendix

- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant		Concentration	High Density Polyethylene	Polyamide	Polycarbonate	Polyoxymethylene	Polypropylene	Polysulfone	Polyvinyl chloride, hard	Polyvinyl chloride, weak	Polytetrafluorethylene	Acrylonitrile-butadiene-caoutchouc	Aluminium	
			Medium	Formula	[%]	HDPE	PA	PC	POM	PP	PSU	PVC	PVC	PTFE
	Sodium bisulphite	NaHSO ₃	10	1	1	2	4	1	-	-	-	1	1	1
	Sodium carbonate	Na ₂ CO ₃	10	1	1	1	1	1	-	-	-	1	-	3
	Sodium chloride	NaCl	30	1	1	1	1	1	1	1	1	1	1	3
	Sodium hydroxide	NaOH	30	1	1	4	1	1	1	1	1	1	2	4
	Sodium hydroxide	NaOH	50	1	1	4	1	1	1	1	-	1	2	4
	Sodium sulfate	Na ₂ SO ₄	10	1	1	1	1	1	1	1	1	1	1	1
	Spirits	C ₂ H ₆ O	96	1	1	1	1	1	1	1	3	1	-	1
	Styrene	C ₈ H ₈	100	4	1	4	1	3	-	4	4	1	4	1
	Sulphuric acid	H ₂ SO ₄	6	1	4	1	4	1	1	1	-	1	2	3
	Sulphuric acid	H ₂ SO ₄	fuming	4	4	4	4	4	4	4	4	1	4	3
	Tallow	—	100	1	1	1	1	1	-	1	1	1	1	1
	Tetrahydrofuran (THF)	C ₄ H ₈ O	100	3	1	4	1	3	4	4	4	1	3	1
	Tetrahydronaphthalene	C ₁₀ H ₁₂	100	3	1	4	1	4	4	4	4	1	-	1
	Thionyl chloride	Cl ₂ SO	100	4	4	4	2	4	4	4	4	1	-	3
	Tin chloride	SnCl ₂	10	1	4	2	2	1	-	-	-	1	1	4
	Transformer oil	—	100	1	1	3	3	1	1	1	-	1	1	1
	Trichloroethane	C ₂ H ₃ Cl ₃	100	3	3	4	2	4	4	4	4	1	4	4
	Urea	CH ₄ N ₂ O	10	1	1	1	1	1	-	-	-	1	1	1
	Urine	—	100	1	1	1	1	1	-	1	1	1	-	2
	Vinegar	C ₂ H ₄ O ₂	10	1	4	1	1	1	1	1	1	1	2	1
	Vinegar	C ₂ H ₄ O ₂	90	1	4	4	4	1	3	1	4	1	-	1
	Wax	—	100	-	1	1		1	-	-	-	1	-	1
	Wines	—	100	1	1	1	2	1	1	1	1	1	-	4
	Xylene	C ₈ H ₁₀	100	3	1	4	1	4	4	4	4	1	4	1

11 Appendix

11.6 Serial Control Interface

V2.1 SERIAL CONTROL INTERFACE SPECIFICATION 

Serial Control Interface Specification

Spincontrol

Version: V2.1
Date: 14.12.2015
File: rs232_spincontrol.doc

PAGE 1

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



1 Contents

1	Contents	2
2	Introduction	3
3	General specifications	3
4	Communication protocol	3
4.1	Reset message	4
4.2	General user commands	4
4.2.1	Overview of commands	4
4.2.2	Control commands	4
4.2.3	Commands to change the setpoints	4
4.2.4	Commands to request process values	5
4.2.5	Commands to request setpoints	5
4.2.6	Commands to change parameters	5
4.2.7	Commands to read parameters	5
4.2.8	Commands to request the status	6
4.2.9	Other commands	6
4.3	Additional commands of Spincontrol S and Professional	9
4.3.1	Commands related to curves	9
4.3.2	Data of last run	10
4.3.3	Commands related to programs	14
4.3.4	Other commands	15
4.4	Additional commands of Spincontrol S	16
4.5	Commands of centrifuges for robot placement	16
4.5.1	Commands for control panel	16
4.5.2	Commands for motor driven lid or hatch	17
4.5.3	Commands for rotor positioning	17
4.5.4	Commands for bucket lifter unit	17
4.5.5	Commands for Rotor Cycle Counter	17
4.5.6	Commands for Servo Cycle Counters	17
4.6	Table of user commands	18
5	Examples	20
6	Hardware interface (optional accessory)	21
6.1	Pinning of the connector	21
6.2	Typical connection to a PC	21
7	History	22

2 Introduction

This document describes the hardware specification and software protocol to communicate with a serial RS232 connection to a sigma centrifuge with Spincontrol electronics.

The serial interface offers the possibility of firmware updates (by service technician), control and monitoring of centrifuge parameters and also the readout of service data like error list and cycles.

The communication data is ASCII coded for easy access with standard terminal software, e.g. "zoc" (<http://www.emtec.com/zoc/>) which offers an easy way to monitor and log the centrifuge process parameters.

The Spincontrol serial protocol is syntax compatible to the older Zent2 protocol used in Sigma Robot centrifuges. In contrast to the Zent2 protocol the character echo is not enabled by default. This protocol is also fully compatible to labworldsoft[®] (<http://www.labworldsoft.com/>), an innovative windows software application for laboratory automation which allows measuring, controlling and regulating of all centrifuge operations.

3 General specifications

Interface standard:	RS232
Baud rate:	9600
Parity:	No
Data bits:	8
Stopbits:	1
Data format:	ASCII

The serial communication works without hardware- or XON/XOFF software handshake.

4 Communication protocol

User commands consist of an ASCII-coded command string and - if needed - a parameter separated by a space. The command parser works non case sensitive.

The character received won't be echoed by the centrifuge processor normally, except if barcode menu is implemented in the centrifuge software. You can tell the centrifuge to echo each character by sending the "**echoon**" command. The user command and the return string of the centrifuge will always be terminated with the characters '0x0A' and '0x0D' (CR and LF).

The command "**cmderror**" can be used to ensure the correct execution of the last command.

The centrifuge outputs a prompt to indicate that it's ready to receive commands. The default prompt is "SIGMA>", but if a name is given to the centrifuge it will be expanded (to give a pc the possibility to distinguish several centrifuges) to "SIGMA xyz>" where "xyz" is the given name.

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



4.1 Reset message

Centrifuges output a message after reset. Detailed output differs by model, but all models output the reset reason first and output is done when the prompt appears. Reset reasons are:

- ~hwreset
loss of power
- ~wdreset
the watch dog timer forced a reset
- ~exreset
reset by external reset pin
- ~swreset
reset initiated by software

4.2 General user commands

The following categories of user commands are available for all models.

4.2.1 Overview of commands

An overview about available commands is output by sending "?" or "??". Both commands are equal and output of available commands depends on model.

? outputs the command list
?? outputs the command list

4.2.2 Control commands

These commands cause an immediate action.

start starts the centrifuge with the set values
stop stops the centrifuge with the pre-adjusted deceleration
fstop stops the centrifuge with the maximal deceleration
door opens the door (only possible when the rotor is stationary and centrifuge is not equipped with a motor driven hatch/lid, see chapter 4.5.2 Commands for motor driven lid or hatch)
reset resets the centrifuge. This command has the same effect as power-on
reseterr resets an error message of type "Log" and "Warning"

4.2.3 Commands to change the setpoints

Commands to change setpoints (*OUT_SP_n y*)

setspeed or **OUT_SP_1** sets the speed
settemp or **OUT_SP_2** sets the temperature (only centrifuges with cooling/heating)
settime or **OUT_SP_3** sets the runtime

4.2.4 Commands to request process values

Commands to request process values (*IN_PV_n*)

<i>speed</i>	or	<i>IN_PV_1</i>	requests the actual rotor speed
<i>temp</i>	or	<i>IN_PV_2</i>	requests the actual temperature (only centrifuges with cooling/heating)
<i>time</i>	or	<i>IN_PV_3</i>	requests the remaining time

4.2.5 Commands to request setpoints

Commands to request setpoints (*IN_SP_n*)

<i>getsetspeed</i>	or	<i>IN_SP_1</i>	requests the set rotor speed
<i>getsettemp</i>	or	<i>IN_SP_2</i>	requests the set temperature (only centrifuges with cooling/heating)
<i>getsettime</i>	or	<i>IN_SP_3</i>	requests the set time

4.2.6 Commands to change parameters

Commands to change parameters (*OUT_PAR_n y*)

<i>setaccel</i>	or	<i>OUT_PAR_1</i>	sets the acceleration
<i>setdecel</i>	or	<i>OUT_PAR_2</i>	sets the deceleration

For Spincontrol Comfort, Spincontrol Professional, Spincontrol L and Spincontrol S the parameter of these commands is the curve nr to be used for acceleration or deceleration. For Spincontrol universal, Spincontrol easy and Spincontrol basic a "0" sets the soft mode and a "1" sets the normal mode. For setdecel there is also the parameter "-1" which sets the free spinout mode.

4.2.7 Commands to read parameters

Commands to request parameters (*OUT_PAR_n*)

<i>getaccel</i>	or	<i>IN_PAR_1</i>	requests the acceleration
<i>getdecel</i>	or	<i>IN_PAR_2</i>	requests the deceleration

11 Appendix

V 2. 1

SERIAL CONTROL INTERFACE SPECIFICATION



4.2.8 Commands to request the status

status requests the status of the centrifuge. The value is displayed decimal.

Value	Normal centrifuge or with motor driven lid	Centrifuge with hatch in the lid
0	Rotor is spinning or door is opening / closing.	Rotor is spinning and the centrifuge is not in positioning mode
1	Rotor is stationary: the door can be opened	- Rotor is stationary or - during positioning (not locked) and/or hatch is not open. The hatch can be opened and the rotor is ready for positioning
2	The door is opened	The hatch is open and the rotor is locked. Ready for loading or unloading.
3	An error has occurred	An error has occurred

status1 advanced status of the centrifuge. The value is displayed hexadecimal.

Bit	Status	Normal centrifuge or with motor driven lid	Centrifuge with hatch in the lid
1..0	00	Door is opening/closing	Hatch is opening/closing or undefined or lid is open
	01	Door is open	Hatch is open
	10	Door is close	Hatch is close
	11	Not used	Not used
3..2	00	Wait	Wait
	01	Door can be opened	Hatch can be opened
	10	Door can be closed	Hatch can be closed
	11	Not used	Hatch can be opened or closed
4	0	No imbalance	No imbalance
	1	Centrifuge shut down with imbalance (only set while centrifuge breaks)	Centrifuge shut down with imbalance (only set while centrifuge breaks)
5	0	Rotor is stopped	Rotor is stopped
	1	Rotor is spinning	Rotor is spinning
6	0	No error	No error
	1	Centrifuge shut down with an error	Centrifuge shut down with an error

status2 advanced status of the centrifuge, only centrifuges for roboter placement. The value is displayed hexadecimal.

Bit	Status	Centrifuge with motor driven lid	Centrifuge with hatch in the lid
0	1	Not implemented	Lid is closed
Bit	Status	Centrifuge without bucket lifter unit	Centrifuge with bucket lifter unit
1	0	Always	Bucket is not at its lower end position
	1	Not implemented	Bucket is at its lower end position
2	0	Always	Bucket is not at its upper end position
	1	Not implemented	Bucket is at its upper end position

4.2.9 Other commands

curr Displays all current parameters tabularly: speed, temp (only centrifuges with cooling/heating), status, status1
 The optional parameter "/tn" outputs the data continuously where 'n' defines the repeat rate in seconds. Entering a '.' stops monitoring. The parameter are separated by '\t'

Example with 5 seconds repeat rate:

```
curr /t5
speed temp status status1
3017 22 0 0020
3009 22 0 0020
3005 22 0 0020
3003 22 0 0020
3002 22 0 0020
.
```

V 2 . 1

SERIAL CONTROL INTERFACE SPECIFICATION



cmderror Displays the error status of the last command .
 The centrifuge returns '1' if no error occurred.

syserror Displays the error status (current error number) of the centrifuge
 The centrifuge returns '0', if no error occurred

In case of error numbers 90, 93 and 95, additional 3 parameters are returned as decimal value with information about states of the servo units (currently only available in models with bucket lifter unit). In detail:

	Parameter 2	Parameter 3	Parameter 4
b15/14 = status of	01 = rotor lock unit	10 = slider unit	11 = bucket lifter unit
b13/12	-	-	-
b11	unknown state	unknown state	unknown state
b10	no catch	-	-
b9	time out slow	time out slow	time out slow
b8	time out fast	time out fast	time out fast
b7	switch error	switch error	switch error
b6	-	-	-
b5	-	over current while closing	-
b4	-	-	-
bit 3 (1 = S4 active)	locked switch	closed switch	up switch
bit 2 (1 = S3 active)	catched switch	nearly closed switch	nearly up switch
bit 1 (1 = S2 active)	-	nearly open switch	nearly down switch
bit 0 (1 = S1 active)	unlocked switch	open switch	down switch

geterr the same as "syserror" (for compatibility with Zent2)

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



info Displays software version and other service information like this (8K):

```

info
Centrifuge Name: 8K
Part No.: 10855, Version: 001
Device Name: Zent5 Controlboard
Part No.: 70926, Version: 001
Software Part No.: 26490 - Software Version: 009
CompilationDate: Mar 27 2007 (14:16:22)
TotalCycles: 70
TempOffset: -8
ImbalOffset: -1
SIGMA Laborzentrifugen GmbH Osterode
www.sigma-zentrifugen.de
Err      Para      Code      Timestamp
02       10       125       0
15       10       124       0
12       10       100       0
12       8        55        0
02       8        40        0
    
```

or this (2-6):

```

info
Cent:    2-6
PN:      10220
Dev:     2-6 Controlboard
PN:      70925
Ver:     001
SW PN:   26487
SW Ver:  017
comp:    Nov 7 2008 (08:58:05)
    
```

The exact output format of this command may vary between different centrifuge types.

geterrpara Because the error list with all parameters is output by "info" command on models with Spincontrol S and Spincontrol Professional only, for some other models the geterrpara command is implemented which outputs error list with all parameters (implementation depends on software version).

echoon This command activates the character echo. Every character will be echoed and the following messages are sent as acknowledge for every single command:

Return Message	Description
OK	Command successful
CNF	Command not found
NEA	Not enough arguments (e.g. set speed value missing)
ERR	Command not possible
CYCLES	"start" command received but max. cycles of rotor or bucket reached -> start command must be sent again as confirmation to ignore cycles

echooff This command de-activates the character echo.

getcurvelist Optional command. Returns curve list with Curve number, Acceleration in rpm/s, Deceleration in rpm/s, if implemented.

Output format:

```
Curve, Accel, Decel  
0,100,100  
1,1600,1600
```

getrotor Requests the selected rotor by rotor list index.

getrotorlist Optional command. Returns rotor list with Rotor, Bucket, minimum Radius, maximum Radius, maximum Speed and maximum Temperature, if implemented.

Output format:

```
Rotor, Bucket, Rmin, Rmax, Nmax, Tmax  
11037,13035,49,133,4000,40  
11171,13299,38,142,4000,40  
11171,13296,65,133,4000,40  
12072,0,80,139,4000,40  
12073,0,58,139,4000,40
```

setrotor Selects a new rotor by rotor list index.

4.3 Additional commands of Spincontrol S and Professional

Models with Spincontrol Professional or Spincontrol S have additional commands.

4.3.1 Commands related to curves

getcurve This command returns the data of a free programmable curve. With the parameter n you can choose the curve between 20 and 29:

```
getcurve 22  
CurveNr: 20  
Interval 1: Time: 130 Speed: 100 LIN  
Interval 2: Time: 60 Speed: 148  
Interval 3: Time: 60 Speed: 194  
Interval 4: Time: 60 Speed: 257  
Interval 5: Time: 60 Speed: 327  
Interval 6: Time: 60 Speed: 526  
Interval 7: Time: 60 Speed: 800  
Interval 8: Time: 30 Speed: 1000  
Interval 9: Time: 600 Speed: 100  
Interval 10: Time: 40 Speed: 100  
TotalTime: 1160
```

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



setcurve This command sets new data for the free programmable curves. The command is followed by the parameter:

setcurve [curveNr],[Lin/Quad],[Int1Time],[Int1Speed] ,[Int2Time],[Int2Speed],etc.

Notice this command is only possible, if no free curve is running!

```
SIGMA> setcurve 22,0,130,112,60,148,60,193,60,256,60,326,60,524,60,1010
OK
SIGMA> getcurve 22
CurveNr: 22
Interval 1: Time: 130 Speed: 112 LIN
Interval 2: Time: 60 Speed: 148
Interval 3: Time: 60 Speed: 193
Interval 4: Time: 60 Speed: 256
Interval 5: Time: 60 Speed: 326
Interval 6: Time: 60 Speed: 524
Interval 7: Time: 60 Speed: 1010
Interval 8: Time: 0 Speed: 0
Interval 9: Time: 0 Speed: 0
Interval 10: Time: 0 Speed: 0
TotalTime: 490
SIGMA>
```

4.3.2 Data of last run

getlastrun This command triggers output of parameters and results of last spin in csv-Format
 This command is only available for Spincontrol S.

While there was no spin since last reset, only centrifuge ID, stored barcodes and string "No data available" will be output. Else Data will be output as follows:

Item	1 st column	2 nd column	3 rd column	Condition
Centrifuge ID	Centrifuge name	not assigned yet		Always
		xy...z		No name assigned
				Name assigned
				3 rd column is empty
Barcode of data structure "Staff Member Identification Number"	Barcode	Staff Member ID	Barcode content	Barcode exists in memory
Barcodes of data structure "Donation Identification Number"	Barcode	Barcode number (1-12)	Barcode content	One row for each barcode set (0 to 12 rows)
Used program (only Spincontrol S from Version number > 050)	Program			Always
		Program number	Program name	Existing program used
			Program name Changed during run	Existing program used, but it was changed during run
	- empty column -		RAPID_TEMP	RAPID_TEMP used

V 2 . 1

SERIAL CONTROL INTERFACE SPECIFICATION



Item	1 st column	2 nd column	3 rd column	Condition
			RAPID_TEMP Changed during run	RAPID_TEMP used, but it was changed during run
			--	No program used
Status	Status of run			Always
		Completed		Run finished already
		Not started		Spin did not start
		Still running		Still running
			Interrupted by error xy	Error during run
			Speed was partly out of setting	Speed error detected by run observation
			Stopped by user	Stop button pressed or shortrun
			Not started	Spin did not start
Blank line			OK	Run OK
Start Time	Start time of last run			Always
		abcd hours, ef minutes, gh seconds ago		Output depends on time since start
			OK Not started	Spin did start Spin did not start
Kind	Kind of last run			Always
		Short run		Short run
		Normal run		Normal run
			Not started OK	Spin did not start Always
Total Time	Total time			Only if started and finished already
		abcd hours, ef minutes, gh seconds		Output depends on total time
			Interrupted OK	Stop button pressed or shortrun Run OK
				Only if normal run started
Run Time	Run time	Infinite		Run time set to infinite
		abcd hours, ef minutes, gh seconds		Output depends on set run time
			Interrupted	Stop button pressed or error happened
			Changed during run	Parameter was changed during run
			Still running OK	Still spinning Run time OK
				Only if normal run started and Item was used
Runtime as of Set Speed	Runtime as of set speed	Active		Only if normal run started and Item was used
			Changed during run	Item was changed during run
			OK	Item OK

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



Item	1 st column	2 nd column	3 rd column	Condition
Deceleration Time	Deceleration time			Only if deceleration time was displayed on screen
		abcd hours, ef minutes, gh seconds		Output depends on deceleration time
			OK	Always
Speed	Speed			If started
		abcde 1/min		Depends on set speed
			Speed was partly out of setting	Speed error detected by run observation
			Not Reached	Set speed was not reached
			Not reached yet	Set speed still not reached
			Changed during run	Set speed was changed during run
			OK	Speed OK
RFC	RCF			If started
		abcde *g		Depends on set RCF
			Speed was partly out of setting	Speed error detected by run observation
			Not Reached	Set RCF was not reached
			Not reached yet	Set RCF still not reached
			Changed during run	Set RCF was changed during run
			OK	RCF OK
Temp	Temperature			Only models with Cooling/Heating if started
		-ab +/- 2 degree Celsius		Output depends on set temperature and set temperature unit
			Not Reached	Set temperature was not reached
			Not reached yet	Set temperature still not reached
			Changed during run	Set temperature was changed during run
			OK	Temperature OK
Rotor	Rotor			If started
		abcde		Output depends on set rotor
			OK	Always
Bucket	Bucket			Only if started and a rotor with bucket is set
		abcde		Output depends on set bucket
			OK	Always
Acceleration	Acceleration			If started
		Curve 9 (Short run)		Short run

V 2 . 1

SERIAL CONTROL INTERFACE SPECIFICATION



Item	1 st column	2 nd column	3 rd column	Condition
		Curve x		Normal run, output depends on set acceleration curve
			Changed during run	Set acceleration curve was changed during run
			OK	Set acceleration curve unchanged
Deceleration	Brake			If started
		Curve 9 (Short run)		Short run
		Curve x (Quick stop)		Quick stop
		Curve x		Normal run, output depends on set deceleration curve
			Changed during run	Set deceleration curve was changed during run
			Still running	Still spinning
			OK	Set deceleration curve unchanged
Spinout	Spinout			Only if normal run started and Item was used (not quick stop)
		From abcd 1/min		Output depends on set Spin out speed
			Changed during run	Item was changed during run
			Still running	Still spinning
		OK	Item OK	

The columns are separated by semicolon.

11 Appendix

V 2. 1

S E R I A L C O N T R O L I N T E R F A C E S P E C I F I C A T I O N



4.3.3 Commands related to programs

setpara This command sets all necessary parameters for a centrifugation at once. Because it's implemented to enable scanning all the centrifugation parameters using a 1D barcode scanner, its command parameters are NOT separated by colons as usual (Code128 has max. data length of 48 characters). Therefore setting up the command parameters has to follow this specification strictly, to guaranty setting the centrifugation parameters correctly.

Parameter number	Meaning	Accepted values
1 to 5	Rotor	Only rotors listed in the centrifuges rotor menu are accepted. 5 characters are mandatory, so fill up rotor number with leading zeros if necessary!
6 to 10	Bucket	Only buckets listed in the centrifuges rotor menu are accepted, but only if they fit to the rotor sent in 1 to 5. 5 characters are mandatory, so fill up bucket number with leading zeros if necessary! If a rotor without buckets is used, set to '00000'.
11 to 13	Radius in mm	'000' (=Rmax) and all values from Rmin to Rmax 3 characters are mandatory, so fill up radius with leading zeros if necessary!
14 to 16	Density in g/cm ³ * 10	'012' to '100' (=1.2g/cm ³ to 10.0g/cm ³) 3 characters are mandatory, so fill up density with leading zeros if necessary!
17	's' for speed, 'r' for RCF	's', 'S', 'r', 'R'
18 to 22	Speed or RCF	Speed: '00100' to maximum speed of rotor and density RCF: Minimum to maximum RCF of rotor and radius 5 characters are mandatory, so fill up speed/RCF with leading zeros if necessary!
23	Sign for temperature value	+', '-' (only centrifuges with cooling/heating, ignored else)
24 to 25	Temperature value	Minimum temperature of centrifuge to maximum temperature of rotor (only centrifuges with cooling/heating, ignored else). 2 characters are mandatory, so fill up temperature with leading zeros if necessary!
26 to 31	Run time in seconds	'000000' (infinite run) and '000010' to '359999' 6 characters are mandatory, so fill up time with leading zeros if necessary!
32 to 33	Acceleration curve	'00' to '19' (always) plus '20' to '29' (but only if the corresponding curve is stored in centrifuge). 2 characters are mandatory, so fill up curve number with leading zero if necessary!
34 to 35	Deceleration curve	'00' to '19' (always) plus '20' to '29' (but only if the corresponding curve is stored in centrifuge). 2 characters are mandatory, so fill up curve number with leading zero if necessary!
36 to 37	Spin-out speed * 100rpm	'00' (no spin-out), '01' to '10' (=spin-out speed from 100rpm to 1000rpm). 2 characters are mandatory, so fill up with leading zero if necessary!
38	Flag "runtime as of set speed"	'0', '1'

Due to this, the parameter length is fixed to 38, so the whole command takes 46 characters. If the length is not exactly 46 characters, the command will be ignored.
 This command is only available for Spincontrol S.

- getpara** This command returns all necessary parameters for a centrifugation at once. It's implemented to enable copying a parameter setting into another centrifuge (in combination with "setpara"). Therefore output is in the same format as expected by command "setpara".
This command is only available for Spincontrol S.
- setprog** This command is to store the actual centrifugation parameters to a program with the given number and name. Therefore, two parameters are mandatory, separated by comma. First parameter specifies the program number, valid from 1 to 60.
Attention: already stored program on this position will be overwritten!
The second parameter is a string with at least one, but up to 19 ASCII characters and specifies the program name.
This command is only available for Spincontrol S.
- getprog** This command returns set program number (1 to 60) and program name. It's implemented to copy programs at the same position with the same name on another centrifuge. Therefore output is in the same format as expected by command "setprog". If no program is set, output is "0,--", if RapidTemp is set, output is "0,RapidTemp".
This command is only available for Spincontrol S.
- getlibr** This command returns all stored user programs, one program per line in format program number (1 to 60), comma, program name, comma and program parameters as returned by command getpara. It's implemented to copy all programs at the same position with the same name and same parameters on another centrifuge. But a corresponding "setlibr" function is not implemented yet.
This command is only available for Spincontrol S.
- loadprog** This command loads a program of the centrifuge. It's only accepted if no centrifugation is in progress.
One parameter is mandatory and specifies the program to load, where valid programs are:
- 0 (only for centrifuges with refrigerator/heater) = RapidTemp program.
Command is only accepted if
→ set temperature is below actual temperature (centrifuge with refrigerator only)
→ set temperature is above actual temperature (centrifuge with heater only)
→ set temperature is different to actual temperature (centrifuge with refrigerator and heater)
→ never (centrifuge without refrigerator or heater)
- 1 to 60 = corresponding program stored in centrifuge.
Command is not accepted if the program doesn't exist.
This command is only available for Spincontrol S.

4.3.4 Other commands

- getname** Displays name of the centrifuge (given by centrifuge menu Setup System Name).
- getprocess** This command gives an overview about the currently set process data (rotor number, bucket number, spd in rpm, time in seconds – 0 is endless, temperature in °C [only centrifuges with cooling/heating], acceleration curve number, deceleration curve number) as well as information about rotor spinning (run = 1) or not (run = 0) and if an error appeared (err = 1) or not (err = 0). It also contains a crc (xor all data) to enable check of correct transmission.

```
getprocess
rotor,bucket,spd,time,temp,acc,dec, run, err,crc
11805, 13850, 200, 0, 20, 9, 29, 0, 0, 207
```

11 Appendix

V 2. 1

S E R I A L C O N T R O L I N T E R F A C E S P E C I F I C A T I O N



4.4 Additional commands of Spincontrol S

- setbarcode** This command adds one barcode (Code128 = ISBT128) to the next centrifugation run. It's possible to add up to 13 barcodes to one run (12 codes of data structure "Donation Identification Number" from blood bags plus 1 code of data structure "Staff Member Identification Number"). The syntax is "setbarcode abc" where abc = content of the barcode. The content of the barcodes is not verified in any matter, only the kind of data structure is checked. Each barcode has to be set separately.
This command is accepted only if
- no centrifugation is in progress and
 - there are no barcodes from an older run in memory (use "deletebarcodes" command to delete them) and
 - Barcode is of data structure "Donation Identification Number" or "Staff Member Identification Number" and
 - same barcode is not already stored and
 - less than 12 barcodes of data structure "Donation Identification Number" are set already if it's a barcode of data structure "Donation Identification Number" and
 - No barcode of data structure "Staff Member Identification Number" is set already if it's a barcode of data structure "Staff Member Identification Number".
- Using this command will also delete all memorized data of the last run (see getlastrun command) to inhibit invalid combination of barcodes with old run data.
If barcode menu is implemented in the centrifuge software, it will come up and display the stored barcodes.
This command is only available for Spincontrol S.
- getbarcodes** This command is always accepted and triggers output of existing barcodes in memory as follows:
"Barcodes abc, def, ghi, ..."
where abc = content of first barcode, def = content of second barcode, ghi = content of third barcode and so on for one up to 12 barcodes of data structure "Donation Identification Number".
If a barcode of data structure "Staff Member Identification Number" is stored, it's output as first barcode with the extension (staff), so output is
"Barcodes abc (staff), def, ghi, ..."
If no barcode is stored, output is
"Barcodes none"
This command is only available for Spincontrol S.
- deletebarcodes** This command deletes all existing barcodes from memory. It's only accepted (even if no barcodes are stored) if no centrifugation is in progress.
If barcode menu is just on display, it will be updated.
This command is only available for Spincontrol S.

4.5 Commands of centrifuges for robot placement

run n Starts the centrifuge with speed n [rpm].

Note: If the centrifuge is equipped with a hatch in the lid, this command closes the hatch and the centrifuge begins to start the run when the hatch is closed.

If the centrifuge is equipped with a bucket lifter unit, the bucket is moved to its lowest position before the run starts.

4.5.1 Commands for control panel

lock Lock buttons and navigation on control unit (control possible via RS232 only)

unlock Unlock buttons and navigation on control unit

4.5.2 Commands for motor driven lid or hatch

close closes the lid / hatch

door opens the lid / hatch

4.5.3 Commands for rotor positioning

setpos n n=0: unlock the rotor
 n>0: go to position n
 the lid must be close for positioning

Note: If the centrifuge is running, this command stops the run automatically and the rotor goes to position n. If the centrifuge is equipped with a hatch in the lid, the hatch opens automatically during positioning. If the centrifuge is equipped with a bucket lifter unit, the bucket is moved to its lowest position before positioning starts.

pos Outputs the position of the rotor in positioning mode

4.5.4 Commands for bucket lifter unit

lift move the bucket to its upper end position

Note: The command is not accepted while the rotor is spinning during run or positioning.

release move the bucket to its lower end position

4.5.5 Commands for Rotor Cycle Counter

An additional rotor cycle counter is implemented for free use by the user. This counter can only be read out by serial interface. The maximum count value is 4294967295. The value will be set to 0 in case of overflow.

rcycle Displays the current rotor cycle counter.

bcycle Displays the current bucket cycle counter. ONLY Spincontrol L

erasercycle Resets the rotor cycle counter to "0". ONLY Spincontrol Universal

4.5.6 Commands for Servo Cycle Counters

liftercycles get cycles of bucket lifter unit

lockcycles get cycles of rotor lock unit

slidercycles get cycles of slider unit

11 Appendix

V 2 . 1

S E R I A L C O N T R O L I N T E R F A C E S P E C I F I C A T I O N



4.6 Table of user commands

The following table contains the available user commands.

Command name	2. name	Parameters	Return values	Unit	Format ¹	Only models with
?	??		list of commands		ASCII	
bcycle			1	cycles	UINT	Robot placement
close						Robot placement
cmderror			1	0..1	UINT	
curr			4	rpm, °C, status, status1	UINT, INT, UINT, HEX	
deletebarcodes						Spincontrol S
door						
echooff				ch. 4.2.9		
echoon				ch. 4.2.9		
erasercycle						Robot placement and Spincontrol Universal
fstop						
getaccel	IN_PAR_1		1	acc. curve nr	UINT	
getbarcodes			ch. 4.4	Barcodes	ASCII	Spincontrol S
getcurve		1	ch. 4.3.1	ch. 4.3.1	ch. 4.3.1	Spincontrol Professional, S
getcurvelist			list	see ch. 4.2.9	see ch. 4.2.9	optional
getdecel	IN_PAR_2		1	dec. curve nr	UINT	
geterr	syserror		1 or 4	error	UINT	
geterrpara			list		ASCII	Except Spincontrol Professional, S
getlastrun			ch. 4.3.2	ch. 4.3.2	ch. 4.3.2	Spincontrol Professional, S
getlibr			see ch. 4.3.3	see ch. 4.3.3	see ch. 4.3.3	Spincontrol S V051
getname			1	Name	ASCII	Spincontrol Professional, S
getpara			1	see ch. 4.3.3	see ch. 4.3.3	Spincontrol S V051
getprocess			10	see ch. 4.3.4	see ch. 4.3.4	Spincontrol Professional, S
getprog			2	see ch 4.3.3	see ch 4.3.3	Spincontrol S V051
getrotor			1	Rotor List Index	UINT	
getrotorlist			list	see ch. 4.2.9	see ch. 4.2.9	optional
getsetspeed	IN_SP_1		1	rpm	UINT	

V 2 . 1

SERIAL CONTROL INTERFACE SPECIFICATION



Command name	2. name	Parameters	Return values	Unit	Format ¹	Only models with
getsettemp	IN_SP_2		1	°C	INT	Temperature control
getsettime	IN_SP_3		1	sec.	UINT	
info			list	ch. 4.2.9		
lift						Bucket lifter unit
liftercycles			1	cycles	UINT	Bucket lifter unit
loadprog		1		Program number	UINT	Spincontrol Professional, S
lock						Robot placement
lockcycles			1	cycles	UINT	Bucket lifter unit
pos			1	1.4	UINT	Robot placement
rcycle			1	cycles	UINT	Robot placement
release						Bucket lifter unit
reset						
reseterr						
run		1		rpm	UINT	Robot placement
setaccel	OUT_PAR_1	1		acc. curve nr	UINT	
setbarcode		1		Barcode	ASCII	Spincontrol S
setcurve		ch. 4.3.1	ch. 4.3.1	ch. 4.3.1	ch. 4.3.1	Spincontrol Professional, S
setdecel	OUT_PAR_2	1		dec. curve nr	UINT	
setpara		1		see ch. 4.3.3	see ch. 4.3.3	Spincontrol Professional, S
setpos		1		1.4	UINT	Robot placement
setprog		2		see ch. 4.3.3	see ch. 4.3.3	Spincontrol Professional, S
setrotor		1		Rotor List Index	UINT	
setspeed	OUT_SP_1	1		rpm	UINT	
settemp	OUT_SP_2	1		°C	INT	Temperature control
settime	OUT_SP_3	1		sec.	UINT	
slidercycles			1	cycles	UINT	Bucket lifter unit
speed	IN_PV_1		1	rpm	UINT	
slart						
status			1	ch. 4.2.8	UINT	
status1			1	ch. 4.2.8	HEX	
status2			1	ch. 4.2.8	HEX	Robot placement
stop						
temp	IN_PV_2		1	°C	INT	Temperature control
time	IN_PV_3		1	sec.	UINT	
unlock						Robot placement

¹ UINT = decimal unsigned integer value; INT = decimal signed integer value; HEX = hexadecimal value

PAGE 19

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



5 Examples

Note: All commands have to be send without quotation marks and brackets!

[CR] and [LF] are ASCII coded control characters (Carriage Return and Linefeed)

setting the setspeed to 1000 rpm:

```
`setspeed 1000[CR][LF]`
```

starting the centrifuge:

```
`start[CR][LF]`
```

requesting the actual rotorspeed:

```
`speed[CR][LF]` answerstring: `1000[CR][LF]`
```

running the centrifuge at 2000 RPM for 2 minutes. Temperature: 5°C:

```
`setspeed 2000[CR][LF]`
```

```
`settemp 5[CR][LF]`
```

```
`settime 120[CR][LF]`
```

```
`start[CR][LF]`
```

requesting the actual status of the centrifuge:

```
`status[CR][LF]` answer string: `0[CR][LF]` (rotor is spinning)
```

or: `1[CR][LF]` (rotor is stationary)

requesting all actual values of the centrifuge:

```
`curr[CR][LF]` answerstring: `speed temp status status1[CR][LF]  
2000 5 1 0004[CR][LF]`
```

to request the actual parameters periodically you have to put a '/tn' behind the command. The 'n' stands for the repeat rate in seconds. Entering a '.' stops monitoring.

requesting the actual values of the centrifuge periodically every 5 seconds:

command:

```
`curr /t5[CR][LF]`
```

answer of the centrifuge:

```
`speed temp status status1[CR][LF]  
2000 5 1 0004[CR][LF]`
```

5 seconds later:

```
2001 5 1 0004[CR][LF]`
```

Stop requesting the actual values of the centrifuge periodically:

command:

```
`curr /t.[CR][LF]`
```

6 Hardware interface (optional accessory)

6.1 Pinning of the connector


Pinning of the 9 pin SUB-D (male) connector on the cover of the centrifuge

PIN 2: TxD (transmit)

PIN 3: RxD (receive)

PIN 5: GND (ground)

PIN 1,4,6,7,8,9: Not Connected

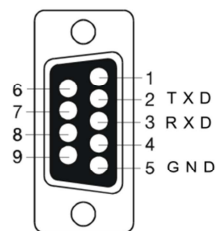

6.2 Typical connection to a PC

Typical connection to a personal computer or a terminal with serial RS232 interface:

9-pin male socket at PC / terminal:

Standard 1:1 serial cable (9 pin female ⇔ 9 pin female)

Required: PIN 2 ⇔ PIN 2 TxD (Cent.) ⇔ RxD (PC)
 PIN 3 ⇔ PIN 3 RxD (Cent.) ⇔ TxD (PC)
 PIN 5 ⇔ PIN 5 GND (Cent.) ⇔ GND (PC)

25-pin male socket at PC / terminal:

Standard 1:1 serial cable (9 pin female ⇔ 9 pin female) + 9-25 way Adaptor, 9 pin male ⇔ 25 pin female

or

Serial cable (9 pin female ⇔ 25 pin female):

Pin-Pin Configuration (only boldface printed required):

D Sub 9 : 1 **2 3** 4 **5** 6 7 8 9
 D Sub 25: 8 **3 2** 20 **7** 6 4 5 22

Note: Do not use a "Null-Modem" cable/adaptor with crossed RxD/TxD signals

11 Appendix

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



7 History

Autor	Datum	Version	Änderungen
D. Hanke	07.11.2008	0.1	Initial Version
D. Hanke	11.11.2008	0.2	status2 and geterr added
D. Hanke	22.02.2010	0.3	Add chapters 5.1 – 5.3
Ch. Seydel	13.12.2010	0.4	Add getname
Ch. Seydel	25.03.2011	0.5	Add bicycle, remove erasercycle
Ch. Seydel	05.03.2012	0.6	Add getlastrun, Spincontrol L and spincontrol S
Ch. Seydel	02.05.2012	0.7	Add Centrifuge ID to getlastrun command
Ch. Seydel	27.05.2012	0.8	Change answer to getlastrun command: - Change order: 1) Centrifuge name 2) Status of run 3) Blank line - Delete Radius and Density
Ch. Seydel	15.10.2012	0.9	Add Barcode commands and add Barcodes to getlastrun command
Ch. Seydel	28.03.2013	1.0	reseterr command limited to reset only errors of type "Log" and "Warning"
Ch. Seydel	06.12.2013	1.1	- update "setbarcode"-, "getbarcodes"-, "deletebarcodes"- and "getlastrun"-commands - add "loadprog" - add the prompt to chapter 4 Communication protocol
Ch. Seydel	20.01.2014	1.2	- add commands "setpara" and "setprog" to chapter 4.9 - getlastrun output strings for "Runtime as of set speed" and "Quick stop" changed to same string used in GUI
Ch. Seydel	03.02.2014	1.3	changed "Break time" to "Deceleration time" in "getlastrun" command
Ch. Seydel	27.08.2014	1.4	- add chapter 4.5.4 Commands for bucket lifter unit - add bucket position to status2 command

PAGE 22

V 2.1

SERIAL CONTROL INTERFACE SPECIFICATION



Autor	Datum	Version	Änderungen
Ch. Seydel	27.08.2014	1.5	<ul style="list-style-type: none"> - add lock and unlock commands to chapter 5.2 Commands for centrifuges with rotor positioning - add door command to chapter 5.1 Commands for centrifuges with motor driven lid or hatch and exception to door command at chapter 4.3 Control commands - add getprocess command for Spincontrol S and Professional - add chapter 4.1 Reset message
Ch. Seydel	29.09.2014	1.6	OK-Return message to "echoon" changed into capital letters (was Ok before)
Ch. Seydel	30.04.2015	1.7	Add hint "(only centrifuges with cooling/heating)" to commands settemp, temp, getsettemp and to corresponding parameters of commands curr, setpara, getprocess
Ch. Seydel	04.05.2015	1.8	Add comment "only set while centrifuge breaks" to bit5 of status1
Ch. Seydel	28.09.2015	1.9	<ul style="list-style-type: none"> - update "getlastrun" - add "geterrpara" - add "liftercycles", "slidercycles", "lockcycles" - update "geterr" 90/93/95 - add CYCLES to return messages - re-arrange chapters
D. Hanke	03.12.2015	2.0	<ul style="list-style-type: none"> - add "getrotor", "setrotor" - add "getcurvelist", "getrotorlist"
Ch. Seydel	14.12.2015	2.1	Add getpara, getprog and getlibr

11 Appendix

11.7 EC declaration of conformity



EC – DECLARATION OF CONFORMITY

The product named hereinafter was developed, designed, and manufactured in compliance with the relevant, fundamental safety and health requirements of the listed EC directives and norms.

In the event of modifications that were not authorised by us or if the product is used in a manner that is not in line with the intended purpose, this declaration will be rendered void.

<i>Product name:</i>	Laboratory centrifuge
<i>Product type:</i>	Sigma 8KS Sigma 8KBS, Sigma 8KHBS
<i>Order number:</i>	10620, 10621, 10622, 10623, 90916, 91032 10625, 10626, 91080, 91143, 91155
<i>Directives:</i>	2006/42/EC Machinery Directive 2014/35/EU Low Voltage Directive 2014/30/EU EMC Directive 2011/65/EU RoHS Directive
<i>Normes:</i>	EN 61010-2-020:2006 EN 61000-3-2:2006, A1:2009, A2:2009 EN 61000-3-3:2013 EN 61326-1:2013

Sigma Laborzentrifugen GmbH

An der Unteren Söse 50
37520 Osterode
Germany

Authorised representative
for CE matters:
Eckhard Tödteberg

Osterode, 01/03/2016




General Manager

CE_8KS_2016-03-01_en

11 Appendix

11.8 Layout plan

here: water cooled centrifuge, GMP version, part no. 10623

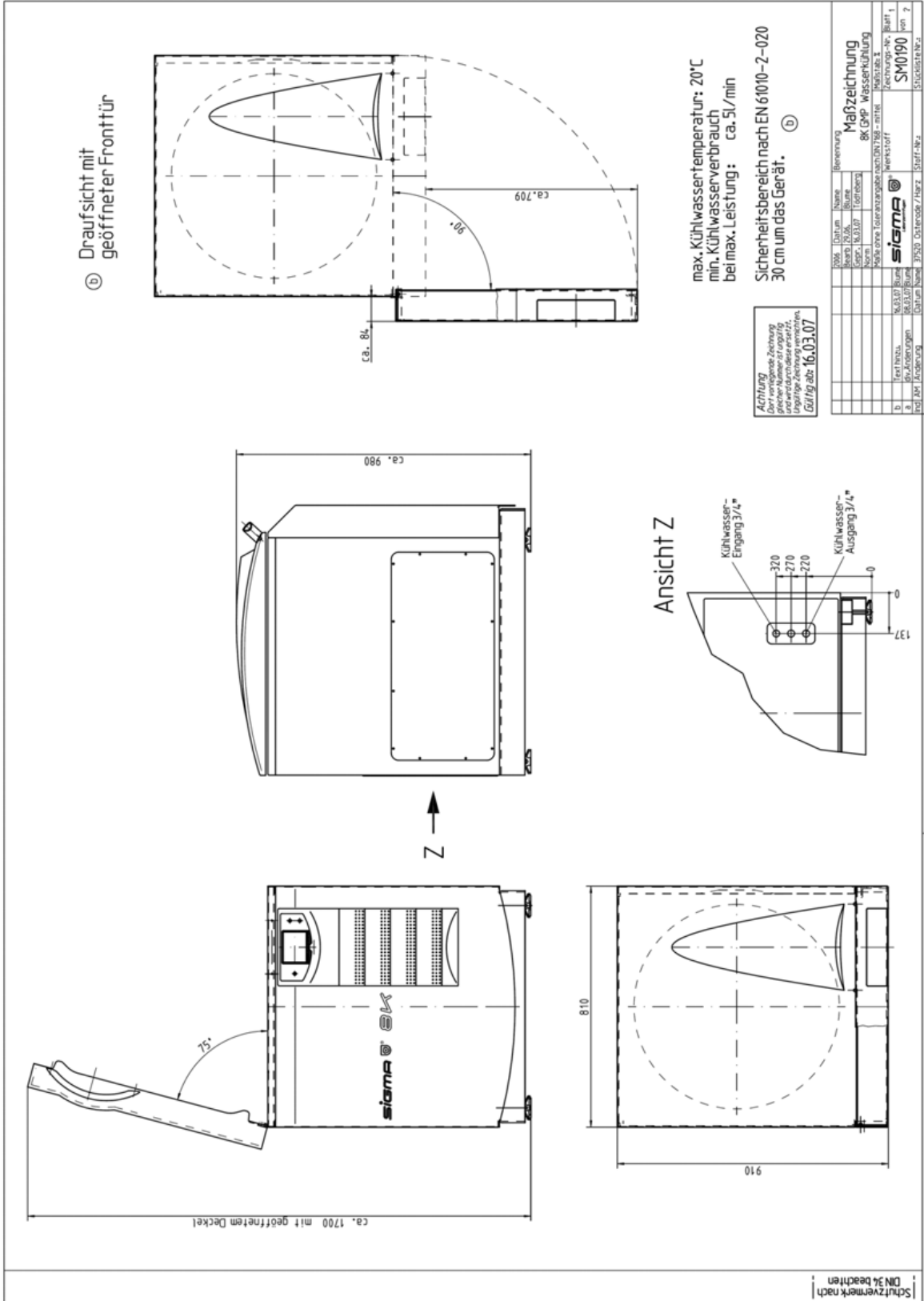

11 Appendix

12 Index

A	
Acceleration	48
Acceleration curve	48, 59, 87
Acceleration curve (changing).....	59
Acceleration curve (creating)	59
Accessories for blood bags.....	84
Accessories, cleaning and care.....	73
Acid.....	22, 72, 74
Acoustic signal.....	68
Activating the help function	62
Adapter	73
Adapters	41
Adapters and stainless steel tubes.....	84
Alignment.....	32
Alignment of the centrifuge	31
Alkaline solutions	22, 72, 74
Ambient conditions.....	82
Ambient temperature	82
Anodised coating	73
Application examples	15
Autoclaving	76
Automatic lid opening function	56
Automatic program rotation.....	66
Automatic rotor identification.....	47
B	
Barcode commands	61
Barcode menu	60
Barcodes (deleting).....	62
Barcodes (reading out)	61
Barcodes (setting).....	61
Blocking a function.....	55
Blood bag systems.....	42
Brake	48
Brakeless deceleration.....	48
Brakeless deceleration from set speed	54
Brief mains power failure.....	68
Bucket.....	41
Buckets, cleaning and care	74
Buzzer/signal	58
C	
Capacity.....	81
Carrier.....	41
Carrier systems.....	41
Carriers, cleaning and care	74
CE mark in compliance with the directive 2006/42/EC	17
Centrifugation monitoring	53
Centrifugation of blood bag systems	42
Centrifugation of infectious, toxic, radioactive, or pathogenic substances.....	21
Centrifugation principle	15
Centrifugation with different tube sizes.....	40
Centrifugation with low capacity	40
Centrifuge cannot be started.....	68
Centrifuge decelerates during operation	68
Centrifuge is inclined.....	68
Centrifuge was moved during run	68
Centrifuge, cleaning and care	72
Centrifuges	
- definition.....	9
Changes in colour	76
Changes in structure.....	76
Changing an acceleration curve.....	59
Changing the code.....	55
Chemical and biological safety.....	21
Chemical resistance of plastic.....	74
Cleaning agents.....	75
Cleaning the bores of angle rotors	73
Cleaning the centrifuge	72
Closing the lid	20, 38
Communication error	70
Condensate	35
Condensate drain.....	35
Condensation.....	31
Condenser dirty	68
Condenser, cleaning and care	73
Connected load.....	81
Contamination.....	21, 73
Continuous heat resistance.....	75
Continuous run	46
Contrast (changing)	63
Control system.....	43
Cooling water connection.....	13
Copyright	10

Index

Corrosion	20, 23, 41, 72, 73, 74, 75	EC declaration of conformity	10
Cost estimate	78	EC-Declaration of Conformity	119
Cracking	73	EEPROM error.....	70
Cracks	74	Electrical connection	81
Creating an acceleration curve.....	59	Electrical safety.....	20
Current consumption.....	81	Emergency lid release	69
Curve menu	59	Equipotential bonding screw	24
Customer-provided fuses	35	Error correction	68
Cycles.....	39, 58	Error message	68
D		Executing a program.....	65
Damage of the surface.....	74	Explanation of the symbols and notes.....	18
Dangerous goods.....	22	Explosive substances	21
Dangerous materials.....	75	External signal	58
Data of the last centrifugation run (reading out).....	62	F	
Deactivating the help function	62	Fire preventions	21
Deceleration (brake)	48	Flow rate.....	37
Deceleration curve.....	48, 59, 87	Form for the return of defective parts	78
Declaration of conformity	10	Function.....	56
Declaration of Conformity.....	119	Functional and operating elements	11
Declaration of decontamination.....	78	Fuses have tripped	68
Declaration of Decontamination	80	G	
Decontamination agent	73, 75	General conditions	9
Deformation of tubes.....	76	Glass breakage.....	75
Deleting a program	65	Glass particles	75
Deleting the stored barcodes	62	GMP version	37
Density.....	22, 53	Gravitational field	81
Dialog box.....	24	Grease for load bearing bolts.....	72, 74
Different service life of rotors and accessories	23, 42, 89	H	
Dimensions	26	Hazard warnings.....	9, 41
Direct hazard to the life and health.....	18	Heater.....	45
Directive 2002/96/EC	80	Help function.....	62
Disinfectants	75	Highly corrosive substances	20
Disinfection of the rotor chamber and accessories	75	I	
Display.....	43, 57	Imbalance	40, 41, 42, 74
Display deceleration time	57	Imbalance error.....	70
Disposal of the centrifuge.....	80	Imbalance monitoring system	24
Disposal of the packaging	80	Importance of the operating manual.....	9
Draining the condensate off	35	Important information	18
Draining the water circuit.....	26	Improper loading.....	68
Drive problem	68	Infectious substances	75
E		Inflammable substances	21
Earth conductor check	24	Informal safety instructions	19
		Initial start-up	38

Input fuse	81	Menu Barcode	60
Inspection by the manufacturer	77	Menu Curve	59
Installation of accessories	40	Menu Help	62
Installation of angle rotors with a hermetically sealed lid	39	Menu Parameters	52
Installation of rotors and accessories	39	Menu Processlibrary	51
Installation of the rotor	39	Menu Setup	56
Intended use	9	Menu Standard	45
Interrupting a centrifugation run	44	Mode of operation	15
Interrupting a deceleration process	44	Modification mode active	44
Invert	57	Motor error	70
IP code	81	Motor shaft	72
K		Multiple carrier	41
Kinetic energy	81	N	
L		Name plate	35
Language	57	No indication on the display	68
Layout of the centrifuge	11	No power in the mains supply	68
Layout plan	121	Noise level	81
Leaks	41, 76	Non-greased bolts	74
Lid cannot be opened	68	Notes on safety and hazards	9
Lid error	70	O	
Lid lock device	24	Online download of forms	79
Lid lock has not released	68	Open lid after run	56
Lid seal sticks	68	Opening the lid	38
Lifting and carrying rotors	23	Operating conditions for water cooled centrifuges	37
Limit "Delta T"	54	Operating personnel	19
Linear curves	87	Operating voltage	35
Load bearing bolts, cleaning and care	74	Operational safety	73
Loading a process	51	Option for the input and output of data	58
Loading a program	64	Option: Barcode menu	60
Lock	55	Options for data input and output	67
M		Overheating	31
Mains plug	14	Overseas shipping	26
Mains power switch off	68	P	
Mains voltage	20	Packaging	27
Maintenance	72	Parameter error	70
Malfunctions and error correction	68	Parameters menu	52
Manual mode	44	Pathogenic substances	21, 72, 75
Manufacturer	81	Plastic accessories, cleaning and care	74
Marking of the unit	17	Potential hazard to the life and health	18
Maximum speed for tubes	85	Potentially hazardous situation	18
Measures in the event of hazards and accidents	25	Power consumption	81
Mechanical safety	20	Power cord	35
		Power cord is not plugged in	68

Index

Power isolating device	35	Rotors, cleaning and care	74
Power supply	35	Runtime	46
Pressure marks.....	74	Runtime as of set speed	56
Prevention of accidents.....	9	S	
Problem description	78	Safety area	22
Process.....	52	Safety devices	24
Process error	70	Safety distance	20, 22
Process library	51	Safety instructions.....	9, 41
Program "RAPID_TEMP".....	48	Safety instructions for centrifugation	22
Program list.....	48, 64, 65	Safety, chemical and biological.....	21
Program mode	63	Safety, electrical	20
Program rotation	56	Safety, mechanical.....	20
Progress indicator	50	Safety-conscious work	19
Protection class	35, 81	Saving a program.....	64
Q		Scope of supply	10
Quadratic curves.....	87	Selection, display, and modification of data	44
Quick stop.....	44	Serial number	71, 77
R		Service.....	77
Radioactive substances	21, 72	Service contact	71
Radius	16, 52	Service life	72
Range of accessories	83	Service life of rotors and accessories.....	23, 42, 89
RAPID_TEMP program.....	48	Service life of the accessories.....	76
RCF	45	Service life of the adapters for blood bags no.	13865 and 13866
Reading out the data of the last centrifugation	run.....	42	
62		Service work	77
Reading out the stored barcodes	61	Setting a barcode.....	61
Refrigerant.....	37	Set-up	32
Refrigeration unit.....	37	Set-up and connection	31
Relative centrifugal force (RCF).....	16, 45	Setup menu	56
Remove glass particles and metal dust from	the rotor chamber	Set-up of the centrifuge.....	31
75		Short run.....	46
Responsibility of the operator.....	19	Solvents.....	22, 72, 74
Return of centrifuges, spare parts, or	accessories	Sound signal	24
78		Sound signal (warning)	58
Return of defective parts	78	Spare part enquiries.....	71
Risk of electrical shock.....	20	Special equipment: water cooling system	37, 81
Rocker switch	35	Speed	16, 45, 81
Rotor chamber	72	Speed/RCF fine	57
Rotor identification, automatic.....	47	Speed-gravitational-field-diagram	86
Rotor monitoring system	24	Speedometer error.....	70
Rotor radii	85	Spin-out from set speed	54
Rotor selection list.....	47	Standard menu	45
Rotor tie-down screw	39, 74	Standards and regulations	10
Rotors and accessories with a different	service life		
23, 42, 89			

Standstill cooling.....	52	Temperature unit.....	57
Standstill monitoring system	24	Temperature value not reached	68
Start delay.....	57	Thermal stress	72
Starting a centrifugation run	44	Time fine	57
Starting a process.....	51	Time range.....	81
Sterilisation of the rotor chamber and accessories	75	Toxic substances	21, 72
Storage and transport	26	Tubes.....	41
Storage conditions	26	Type.....	81
Storage locations	81	Type of connection.....	35
Stress-corrosion see		Type of the centrifuge	71, 77
corrosion	74	U	
Structural changes.....	21	Unblocking a function.....	55
Supply voltage	35	Ungreased load- bearing bolts	68
Supporting rings.....	85	Useful volume	
Switching the centrifuge off	67	- volume that is stated for the tube	41, 85
Switching the centrifuge on	38	User interface	43
System.....	58	Using smaller blood bag systems.....	42
System check	24	UV radiation	31, 72
System error	68, 70	V	
T		Vessels	41
Table of error codes.....	70	W	
Table of rotors and accessories with a different service life	23, 42, 89	Warning signal	58
Technical data	81	Warranty and liability.....	9
Technical documentation	82	Water circuit (draining).....	26
Temperature	45	Water connector	37
Temperature at the water inlet	37	Water cooling system.....	37, 81
Temperature error.....	70	Draining the water circuit	26
Temperature inside the rotor chamber	24	Water operating pressure.....	37
Temperature monitoring.....	54	Wear	77
Temperature monitoring system.....	24	Weight	26
		Work on the power supply system	20