


CERTIFICATE OF CALIBRATION

ISSUED BY: LAMBDA CALIBRATION LTD

DATE OF ISSUE: 15 October 2019 CERTIFICATE No: 576962

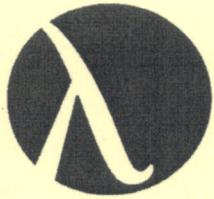


Lambda
CALIBRATION LTD

Units 11-13
Chorley Central Business Park
Stump Lane, Chorley
Lancashire PR6 0BL
Tel: 01257 244670

Page 1 of 3

APPROVED SIGNATORY

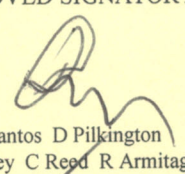
E Santos D Pilkington
D Whalley C Reed R Armitage

Customer: DJB Labcare Ltd
Address: 20 Howard Way, Interchange Park
Newport Pagnell
MK16 9QS

Item Number: 17070127 (4046)
Description: Digital Thermometer
Model/Range: TMD-56
Manufacturer: Amprobe
Date of Cal: 15 Oct 2019
Calibrated by: Mohammed Abid
Procedure Name: Amprobe, Digital Thermometer, TMD-56 (DJB Labcare)
Rev/Basis: 03:E-650, Based on BS EN 60584.1
Temp/Humidity: 20.0°C ± 2°C <80%rh

The Results on the following pages are: As Found

All Measurements are Traceable to National Standards.

Note 1: The unit under test was calibrated using a multifunction calibrator.
Note 2: Where the reported value lies within the stated tolerances then this will be indicated by 'A', if outside then by 'B'.
Note 3: Values quoted in the "UUT Indicated Value" column are not necessarily quoted to the same resolution as the actual displayed value on the UUT.
Note 4: Any supplied test leads have been checked as part of the Visual/Operational test but have not been used during calibration.
Note 5: Temperature indicating instruments that contain an internal reference junction for use with thermocouples are calibrated with the reference junction enabled.

Engineers' Notes: Batteries replaced.

Equipment Used: Multi-function Calibrator: LMMC-02 / LMMC-04 / LMMC-10 ✓ / LMMC-14

The reported expanded uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a coverage probability of approximately 95%. The uncertainty evaluation has been carried out in accordance with UKAS requirements. Unless otherwise stated any reported summary of the results does not take the measurement uncertainty into consideration. This certificate is issued in accordance with the laboratory accreditation requirements of the United Kingdom Accreditation Service. It provides traceability of measurement to the SI system of units and / or to units of measurement realised at the National Physical Laboratory or other recognised national metrology institutes. This certificate may not be reproduced other than in full, except with the prior written approval of the issuing laboratory. Unless otherwise stated the reported activities were carried out at the address detailed in the header; and the results relate only to the items calibrated.

CERTIFICATE OF CALIBRATION

ISSUED BY: LAMBDA CALIBRATION LTD

UKAS ACCREDITED CALIBRATION LABORATORY No: 0495

CERTIFICATE No:

576962

Page 2 of 3

Parameter	UUT Range	UUT Indicated Value	Applied Value	Acceptance Limits Low	Acceptance Limits High	Summary
Visual/Operational Test						
Result of Operator Evaluation						OK
Measurement of Thermocouples (Electrical Simulation)						
Channel T1						
Type T		-190.0°C	-190.0	-190.8	-189.2	A
		-80.0°C	-80.0	-80.7	-79.3	A
		-50.0°C	-50.0	-50.7	-49.3	A
		-30.0°C	-29.8	-30.3	-29.7	A
		-10.0°C	-10.0	-10.3	-9.7	A
		0.0°C	0.3	-0.3	0.3	A
		4.0°C	4.0	3.7	4.3	A
		37.0°C	37.0	36.7	37.3	A
		50.0°C	50.0	49.7	50.3	A
		100.0°C	100.2	99.7	100.3	A
		150.0°C	150.0	149.6	150.4	A
		200.0°C	200.0	199.6	200.4	A
		250.0°C	250.0	249.6	250.4	A
		300.0°C	300.0	299.6	300.4	A
		390.0°C	390.1	389.5	390.5	A
		100.0°F	100.3	99.3	100.7	A
Type K		0.0°C	0.0	-0.3	0.3	A
		500.0°C	500.0	499.4	500.6	A
		1000.0°C	999.7	999.2	1000.8	A
Type J		20.0°C	20.0	19.7	20.3	A
		1100.0°C	1100.0	1099.2	1100.8	A
Type E		20.0°C	20.0	19.7	20.3	A
		900.0°C	900.0	899.3	900.8	A
Type N		20.0°C	20.0	19.6	20.4	A
		1100.0°C	1100.1	1099.1	1101.0	A
Type R		500.0°C	500.0	497.8	502.3	A
		1100.0°C	1100.0	1097.5	1102.6	A
Type S		500.0°C	500.0	497.8	502.3	A
		1100.0°C	1100.0	1097.5	1102.6	A

CERTIFICATE OF CALIBRATION

ISSUED BY: LAMBDA CALIBRATION LTD

UKAS ACCREDITED CALIBRATION LABORATORY No: 0495

CERTIFICATE No:

576962

Page 3 of 3

Parameter	UUT Range	UUT Indicated Value	Applied Value	Acceptance Limits Low	Acceptance Limits High	Summary
Channel T2						
Type T		-190.0°C	-190.0	-190.8	-189.2	A
		-80.0°C	-80.0	-80.7	-79.3	A
		-50.0°C	-50.0	-50.7	-49.3	A
		-30.0°C	-29.9	-30.3	-29.7	A
		-10.0°C	-9.9	-10.3	-9.7	A
		0.0°C	0.1	-0.3	0.3	A
		4.0°C	4.0	3.7	4.3	A
		37.0°C	37.0	36.7	37.3	A
		50.0°C	50.0	49.7	50.3	A
		100.0°C	100.1	99.7	100.3	A
		150.0°C	150.0	149.6	150.4	A
		200.0°C	200.0	199.6	200.4	A
		250.0°C	250.0	249.6	250.4	A
		300.0°C	300.0	299.6	300.4	A
		390.0°C	390.0	389.5	390.5	A
		100.0°F	100.1	99.3	100.7	A
Type K		0.0°C	0.0	-0.3	0.3	A
		500.0°C	500.0	499.4	500.6	A
		1000.0°C	999.8	999.2	1000.8	A
Type J		20.0°C	19.9	19.7	20.3	A
		1100.0°C	1100.0	1099.2	1100.8	A
Type E		20.0°C	20.0	19.7	20.3	A
		900.0°C	900.0	899.3	900.8	A
Type N		20.0°C	20.1	19.6	20.4	A
		1100.0°C	1100.2	1099.1	1101.0	A
Type R		500.0°C	500.0	497.8	502.3	A
		1100.0°C	1100.0	1097.5	1102.6	A
Type S		500.0°C	500.0	497.8	502.3	A
		1100.0°C	1100.0	1097.5	1102.6	A

End of Calibration Data

Estimated Uncertainty of Measurement:

Electrical Measurement of Thermocouples

Type: B	+500°C to +1820°C	±(0.64°C)
Type: C	+0°C to +2320°C	±(0.48°C)
Type: E	-250°C to +1000°C	±(0.53°C)
Type: J	-210°C to +1200°C	±(0.30°C)
Type: K	-200°C to -250°C	±(0.66°C)
Type: K	-200°C to +1300°C	±(0.32°C)
Type: L	-200°C to +900°C	±(0.31°C)
Type: N	-200°C to +1300°C	±(0.40°C)
Type: R	+0°C to +1767°C	±(0.61°C)
Type: S	+0°C to +1767°C	±(0.57°C)
Type: T	-250°C to -200°C	±(0.69°C)
Type: T	-200°C to +400°C	±(0.32°C)